



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## TECHNICAL DATA SHEET



bioauxilium

BETTER TOOLS. REAL DISCOVERIES.

# THUNDER™ Phospho-SLP-76 (S376) + Total SLP-76 TR-FRET Cell Signaling Assay Kit

**CATALOG NUMBERS** KIT-SLP76PT-500  
400 points for phospho-SLP-76  
and 100 points for total SLP-76

Store at -80°C  
For research use only.  
Not for use in diagnostic procedures.

## PRODUCT DESCRIPTION

This assay kit measures intracellular levels of **phospho-SLP-76 (S376)** and **total SLP-76** protein in cell lysates using a simple, rapid and sensitive immunoassay based on the homogeneous (no-wash) THUNDER™ TR-FRET technology. The kit is compatible with both adherent and suspension cells.

## SPECIFICITY

This assay kit contains two specific and selective antibody pairs, one that recognizes **SLP-76** phosphorylated at **Ser376** and another that recognizes **total** (both phosphorylated and unphosphorylated) **SLP-76**.

## SPECIES REACTIVITY

Human (Swiss-Prot Acc.: Q13094; Entrez-Gene Id: 3937).

Other species should be tested on a case-by-case basis.

## TR-FRET ASSAY PRINCIPLE

The **Phospho-SLP-76 (S376) + Total SLP-76** assay kit is a homogeneous time-resolved Förster resonance energy transfer (TR-FRET) sandwich immunoassay (Figure 1). The THUNDER™ Cell Signaling assay workflow consists of 3 steps (Figure 2). Following cell treatment, cells are first lysed with the specific Lysis Buffer provided in the kit. Then **Phospho-SLP-76 (S376)** and **Total SLP-76** in the cell lysates are detected in separate wells with two pairs of fluorophore-labeled antibodies in a simple "add-incubate-measure" format (single-step reagent addition; no wash steps). For detection of the phosphorylated protein, one antibody is labeled with a donor fluorophore (Europium chelate; Eu-Ab1) and the second with a far-red acceptor fluorophore (FR-Ab2). The same approach is used for the second antibody pair detecting the total protein (Eu-Ab3 and FR-Ab4). The binding of the two matched labeled antibodies to distinct epitopes on the target protein (either **phospho-SLP-76 (S376)** or **total SLP-76**) takes place in solution and brings the two dyes into close proximity. Excitation of the donor Europium chelate molecules with a flash lamp (320 or 340 nm) or a laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which in turn emit a TR-FRET signal at 665 nm. Residual energy from the Eu chelate generates light at 615 nm. The signal at 665 nm is proportional to the concentration of **Phospho-SLP-76 (S376)** and **Total SLP-76** in the cell lysate. Data can be expressed as either the signal at 665 nm or the 665 nm/615 nm ratio.

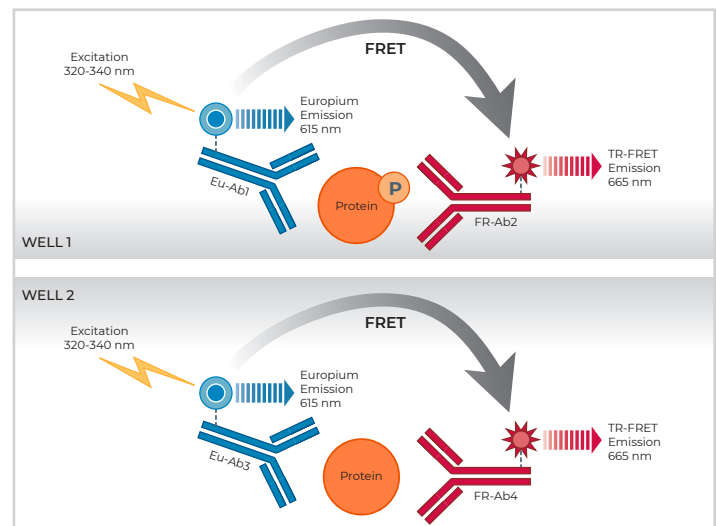


Fig. 1 Schematic representation of the TR-FRET cell signaling assay principle.

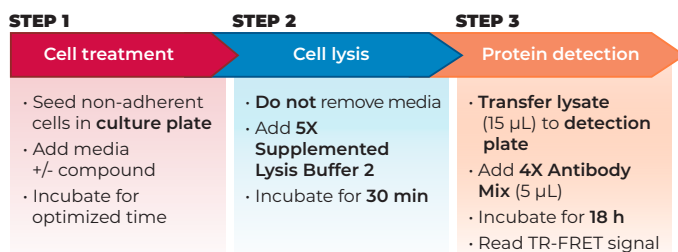


Figure 2 Assay workflow using the 2-plate (transfer) proto-

## KIT COMPONENTS

	500 points*
Eu-labeled phospho-SLP-76 (S376) antibody (Eu-Ab1)	20 µL
Acceptor-labeled phospho-SLP-76 (S376) antibody (FR-Ab2)	80 µL
Eu-labeled total-SLP-76 antibody (Eu-Ab3)	5 µL
Acceptor-labeled total-SLP-76 antibody (FR-Ab4)	20 µL
Lysis Buffer 2 (5X)	5 mL
Detection Buffer (10X)	250 µL
Positive control cell lysate	200 µL
Phosphatase Inhibitor Cocktail (100X)	250 µL

\* The number of assay points is based on an assay volume of 20 µL in half-area 96-well or low-volume 384-well assay plates using the kit components at the recommended concentrations (refer to the User Manual).

# TECHNICAL DATA SHEET


Phospho-SLP-76 (S376) + Total SLP-76

## VALIDATION DATA

This assay kit has been validated for the relative quantification of phospho-SLP-76 (S376) and total SLP-76 in Jurkat cell lysates using the 2 plate assay protocol.


- Non-adherent cells were cultured in RPMI+10% FBS before being centrifuged and resuspended at the desired density in RPMI without serum.
- Following cell treatment, cells were lysed with the 5X **Lysis Buffer 2** (to a final concentration of 1X) supplemented with the phosphatase inhibitors sodium fluoride (5 mM) and sodium orthovanadate (10 mM).
- Following a **30-min** incubation at room temperature (RT) on an orbital shaker (400 rpm), lysates (15  $\mu$ L) were then transferred to a 384-well assay plate followed by addition to separate wells of either the labeled antibodies Eu-Ab1 and FR-Ab2 (5  $\mu$ L) for detection of phospho-SLP-76 (S376) or Eu-Ab3 and FR-Ab4 (5  $\mu$ L) for detection of total SLP-76.
- The plate was incubated at RT for **18 hours** and the TR-FRET signal was recorded at 665 and 615 nm (EnVision<sup>®</sup>; lamp excitation).

## STIMULATION OF PHOSPHO-SLP-76 (S376) IN JURKAT

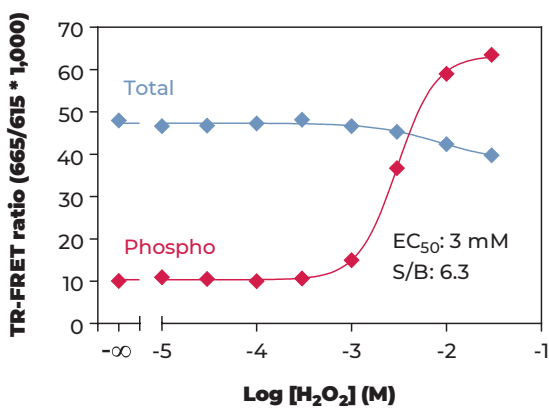


Jurkat cells (400,000 cells/well; in triplicate) were incubated with serial dilutions of H<sub>2</sub>O<sub>2</sub> for 15 min at RT. Data show that treatment of Jurkat cells with H<sub>2</sub>O<sub>2</sub> stimulates phosphorylation of SLP-76 at S376 but does not affect the levels of Total SLP-76.

## JURKAT CONTROL LYSATE TITRATION (QC TEST) PHOSPHO-SLP-76 (S376)

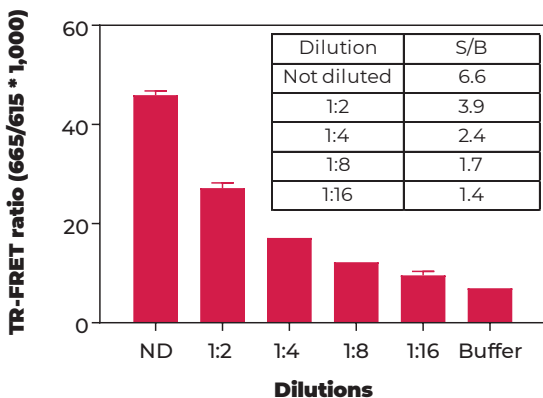

## JURKAT CONTROL LYSATE TITRATION (QC TEST) TOTAL SLP-76

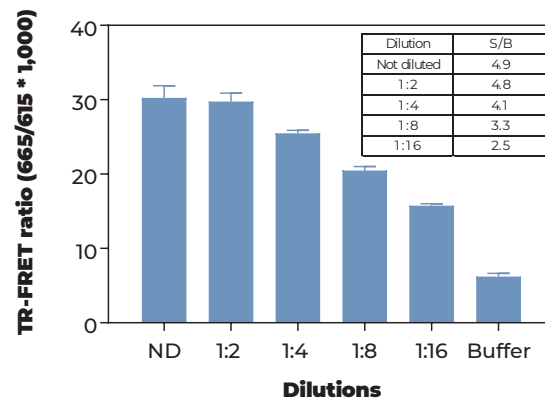

Quality Control: the Phospho-SLP-76 (S376) + Total SLP-76 assay kit is routinely tested against H<sub>2</sub>O<sub>2</sub> treated Jurkat lysates. Jurkat cells were cultured in a T175 flask, centrifuged and resuspended at 20 million cells/mL, and stimulated with 30 mM of H<sub>2</sub>O<sub>2</sub> for 15 min at RT. Following cell lysis using 5X Lysis Buffer 2 (final concentration 1X), lysates were serially diluted with 1X Lysis Buffer 2 and tested in triplicate and in separate wells for phospho-SLP-76 (S376) and total SLP-76. Data show a linear relationship between lysate dilutions and TR-FRET ratio values. Note that due to the very high sensitivity of the Total SLP-76 kit, lysates from the T175 flask required at least a 1:2 pre-dilution in order to be within the dynamic assay range.



bioauxilium

FOR MORE INFORMATION ON DEVELOPING AND OPTIMIZING TR-FRET CELL SIGNALING ASSAYS, CONSULT THE USER MANUAL.

7171 Frederick-Banting, Suite 2115, Montréal (Québec) H4S 1Z9 CANADA [bioauxilium.com](http://bioauxilium.com) 1-866-610-6933 x 103