



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.202.1401
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet
JMJD3 Homogeneous Assay Kit
Catalog #50416-1
Size: 96 reactions

DESCRIPTION: The *JMJD3 Homogeneous Assay Kit* is designed to measure JMJD3 activity for screening and profiling applications. JMJD3 is a JmjC-domain protein that exhibits demethylation activity toward H3-K₂₇Me² and H3-K₂₇Me³. The *JMJD3 Homogeneous Assay Kit* comes in a convenient AlphaLISA[®] format (Scheme 1), with biotinylated histone H3 peptide substrate, primary antibody, demethylase assay buffer, and purified JMJD3 for 96 enzyme reactions. The key to the *JMJD3 Homogeneous Assay Kit* is a highly specific antibody that recognizes demethylated substrate. With this kit, only three simple steps on a microtiter plate are required for methyltransferase detection. First, a sample containing JMJD3 enzyme is incubated with the biotinylated substrate. Next, acceptor beads and primary antibody are added, then donor beads, followed by reading the Alpha-counts.

COMPONENTS:

Catalog #	Component	Amount	Storage	
50115	JMJD3 (KDM6B)	4 µg	-80°C	(Avoid freeze/thaw cycles!)
52140F	Primary antibody 6	5 µl	-80°C	
	Biotinylated histone H3 peptide substrate	300 rxns	-80°C	
52408	4x HDM Assay Buffer 3	3 ml	-80°C	
52301	4x Detection buffer	2 ml	-20°C	

MATERIALS REQUIRED BUT NOT SUPPLIED:

AlphaLISA[®] anti-rIgG acceptor beads, 5 mg/ml (PerkinElmer #AL104C)
AlphaScreen[®] Streptavidin-conjugated donor beads, 5 mg/ml (PerkinElmer #6760002S)
Optiplate -96 (PerkinElmer #6005290)
AlphaScreen[®] microplate reader
Adjustable micropipettor and sterile tips

APPLICATIONS: Great for studying enzyme kinetics and HTS applications.

CONTRAINDICATIONS: Green and blue dyes that absorb light in the AlphaScreen signal emission range (520-620 nm), such as Trypan Blue. Avoid the use of the potent singlet oxygen quenchers such as sodium azide (NaN₃) or metal ions (Fe²⁺, Fe³⁺, Cu²⁺, Zn²⁺ and Ni²⁺). The presence of >1% RPMI 1640 culture medium leads to a signal reduction due to the presence of excess biotin and iron in this medium. MEM, which lacks these components, does not affect AlphaScreen assays.


OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com

SAFETY: This product is for research purposes only and not for human or therapeutic use. This product should be considered hazardous. Do not ingest, inhale, get in eyes, on skin, or on clothing. If so, wash thoroughly.

STABILITY: At least 6 months from date of receipt when stored as directed.

Scheme 1: Our histone demethylase assays utilize highly specific antibodies that recognize demethylated histone substrate. First, a sample containing the enzyme is incubated with a biotinylated substrate. Next, acceptor beads and primary antibody are added, then donor beads, followed by reading the Alpha-counts, as shown below.

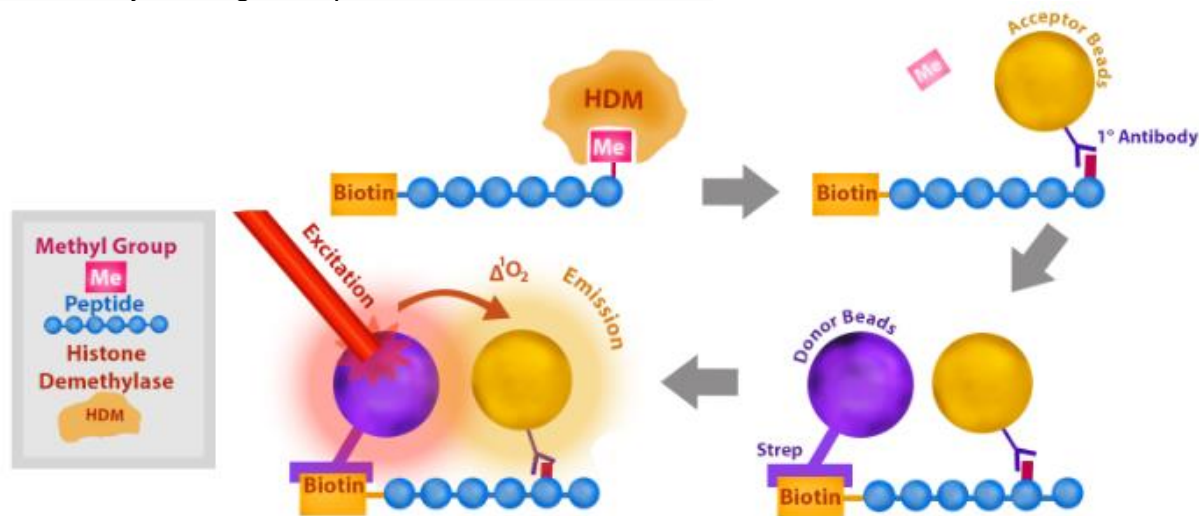

REFERENCE: Swigut, T., Wysocka, J. *Cell* 2007; **131**(1): 29-32.

ASSAY PROTOCOL:

All samples and controls should be tested in duplicate. We recommend preincubating the enzyme with inhibitor, however, it is acceptable to add the substrate mixture and inhibitor followed by diluted JMJD3 without the preincubation step.

Step 1:

- 1) Re-suspend lyophilized **Biotinylated histone H3 peptide substrate** in 950 µl of distilled water.
- 2) Prepare serial dilutions of the test inhibitors in **1x HDM Assay Buffer 3** (Scheme 2). Add 6 µl of inhibitor solution to each well designated "Test Sample". For the wells designated

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com

"Blank" and "Positive Control" add 6 μ l of the same solution without inhibitor (typically **1x HDM Assay Buffer 3** with respective concentration of DMSO).

- 3) Thaw **JMJD3** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full content of the tube. Aliquot **JMJD3** enzyme into single use aliquots. Store remaining undiluted enzyme in aliquots at -80°C immediately. *Note: **JMJD3** is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme.*
- 4) Dilute **JMJD3** in **1x HDM Assay Buffer 3** at 1.50 ng/ μ l (9 ng/6 μ l). Keep diluted enzyme on ice until use. Discard any unused diluted enzyme after use.
- 5) Preincubate 6 μ l of diluted **JMJD3** with 6 μ l of diluted inhibitor(s) for up to 30 minutes at room temperature, with slow shaking. For the wells designated as "Blank", add 6 μ l **1x HDM Assay Buffer 3**.
- 6) Prepare master mix: N wells \times (5 μ l **4x HDM Assay Buffer 3** + 3 μ l **Biotinylated substrate**).
- 7) Initiate reaction by adding 8 μ l of master mix prepared as described above. Incubate at room temperature for one hour. *Note: All incubations are done with slow shaking on a rotator platform.*

Scheme 2: The test compound is first serially diluted in 100% DMSO with the highest concentration at (X) mM. Each intermediate compound dilution (in 100% DMSO) will then get directly diluted 30x fold into **1x HDM Assay Buffer 3** for 3.3x concentration (DMSO). From this intermediate step, 6 μ l of compound is added to 6 μ l of demethylase enzyme dilution is incubated for 30 minutes at room temperature. After this incubation, 8 μ l of master mix (peptide substrate + 4x HDM Assay Buffer) is added. The final DMSO concentration is 1% for all wells.

Reagent	Blank	Positive Control	Test Inhibitor
JMJD3 (1.5 ng/ μ l)	-	6 μ l	6 μ l
1x HDM Assay Buffer 3	6 μ l	-	-
Test Inhibitor/Activator	-	-	6 μ l
1x HDM Assay Buffer 2 (3.3% DMSO)	6 μ l	6 μ l	-
4x HDM Assay Buffer 3	5 μ l	5 μ l	5 μ l
Biotinylated Substrate	3 μ l	3 μ l	3 μ l
Total	20 μl	20 μl	20 μl

Step 2:

Note: Protect your samples from direct exposure to light!

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**


Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com

- 1) Dilute anti-rabbit Acceptor beads (PerkinElmer #AL104C) (1:500) and Primary antibody 6 (1:250) with 1x Detection buffer in one step. Add 10 μ l of acceptor beads/antibody mixture per well. Incubate 30 min at room temperature.

Step 3:

- 1) Dilute **Streptavidin-conjugated donor beads** (PE #6760002S) 125-fold with **1x Detection buffer**. Add 10 μ l of donor beads per well. Shake on a rotator platform for 30 minutes at room temperature.
- 2) Read Alpha-counts.

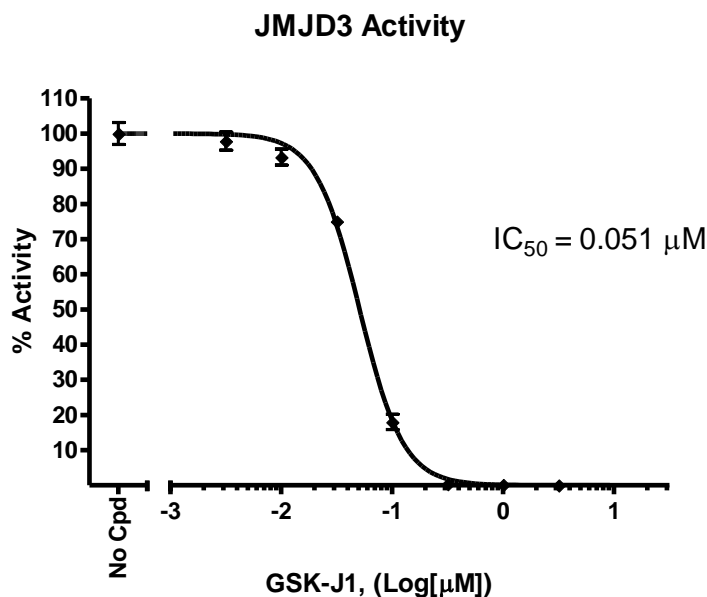Example of Assay Results:

JMJD3 enzyme activity, measured using the *JMJD3 Homogeneous Assay Kit*, BPS Bioscience Cat. #50416-1. *Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at info@bpsbioscience.com*

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.202.1401
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

RELATED PRODUCTS:

<u>Product Name</u>	<u>Catalog #</u>	<u>Size</u>
JMJD3 Assay Kit, Homogeneous	50416-2	384 reactions
JMJD2A Assay Kit, Homogeneous	50413	384 reactions
JMJD2B Assay Kit, Homogeneous	50414-2	384 reactions
JMJD2C Assay Kit, Homogeneous	50415	384 reactions
JMJD2D Assay Kit, Homogeneous	79838	384 reactions
JMJD2E Assay Kit, Homogeneous	50417	384 reactions
JMJD2C Assay Kit, Chemiluminescence	50405	96 reactions
JMJD2D Assay Kit, Chemiluminescence	50418	96 reactions
JMJD3 recombinant protein	50115	100 µg
JMJD2A recombinant protein	50123	100 µg
JMJD2B recombinant protein	50111	100 µg
JMJD2C recombinant protein	50105	100 µg
JMJD2D recombinant protein	50117	100 µg
JMJD2E recombinant protein	50118	100 µg

AlphaScreen® and AlphaLISA® are registered trademarks of PerkinElmer, Inc.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**
Or you can Email us at: [**info@bpsbioscience.com**](mailto:info@bpsbioscience.com)
Please visit our website at: [**www.bpsbioscience.com**](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B
 San Diego, CA 92121
Tel: 1.858.202.1401
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Alpha-counts signal of positive control reaction is same as "blank" value.	JMJD3 has lost activity	Enzyme loses activity upon repeated freeze/thaw cycles. Use fresh JMJD3, BPS Bioscience #50115. Store enzyme in single-use aliquots. Increase time of enzyme incubation. Increase enzyme concentration.
	Streptavidin Donor beads or anti-rIgG acceptor beads fail to show significant signal.	Reorder Streptavidin Donor beads or anti-rIgG acceptor beads from Perkin Elmer.
	Incorrect settings on instruments	Refer to instrument instructions for correct settings to increase sensitivity of light detection.
Alpha-counts signal is erratic or varies widely among wells	Inaccurate pipetting/technique	Run duplicates of all reactions. Use a multichannel pipettor. Use master mixes to minimize errors.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401**, Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com