



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court West, Ste. B  
San Diego CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**Data Sheet**  
**PRMT5 Chemiluminescent Assay Kit**  
Catalog #52002L  
Size: 96 reactions

**DESCRIPTION:** The *PRMT5 Chemiluminescent Assay Kit* is designed to measure PRMT5 activity for screening and profiling applications. The *PRMT5 Chemiluminescent Assay Kit* comes in a convenient format, with a 96-well plate precoated with histone H4 peptide substrate, the antibody against methylated arginine3 residue of Histone H4, the secondary HRP-labeled antibody, S-adenosylmethionine, methyltransferase assay buffer, and purified PRMT5 enzyme for 96 enzyme reactions. The key to the *PRMT5 Chemiluminescent Assay Kit* is a highly specific antibody that recognizes methylated R3 residue of Histone H4. With this kit, only three simple steps on a microtiter plate are required for methyltransferase detection. First, S-adenosylmethionine is incubated with a sample containing assay buffer and methyltransferase enzyme. Next, primary antibody is added. Finally, the plate is treated with an HRP-labeled secondary antibody followed by addition of the HRP substrate to produce chemiluminescence that can then be measured using a chemiluminescence reader.

**COMPONENTS:**

Catalog #	Components	Amount	Storage	
51045	PRMT5/MEP50 human enzyme	20 µg	-80°C	<b>Avoid freeze/thaw cycles!</b>
52120	100 µM S-adenosylmethionine	250 µl	-80°C	
52150-3	Primary antibody 4-3	100 µl	-80°C	
52131H	Secondary HRP-labeled antibody 2	10 µl	-80°C	
	4x PRMT5 assay buffer (add DTT before experiment)	3 ml	-20°C	
79556	Blocking buffer	50 ml	+4°C	
	HRP chemiluminescent substrate A (transparent bottle)	6 ml	+4°C	
	HRP chemiluminescent substrate B (brown bottle)	6 ml	+4°C	
	96-well plate precoated with histone substrate	1	+4°C	

**MATERIALS REQUIRED BUT NOT SUPPLIED:**

TBST buffer (1 x TBS, pH 8.0, containing 0.05% Tween-20)  
DTT (Dithiothreitol), 0.5M (Sigma, Cat. # D0632)  
Luminometer or fluorescent microplate reader capable of reading chemiluminescence  
Adjustable micropipettor and sterile tips  
Rotating or rocker platform

*\*Since May of 2017, this kit has been improved with a more specific antibody and reformulated buffers. The previous version of PRMT kit #52002L can still be purchased upon special request.*

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court West, Ste. B  
San Diego CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**APPLICATIONS:** Great for studying enzyme kinetics and HTS applications.

**CONTRAINDICATIONS:** DMSO >1%, strong acids or bases, ionic detergents, high salt

**STABILITY:** One year from date of receipt when stored as directed.

**REFERENCE:** Stopa, N., *et al.* 2015. *Cell Mol Life Sci.* **72(11)**:2041-59.

**ASSAY PROTOCOL:**

**All samples and controls should be tested in duplicate.**

**Step 1:**

- 1) Rehydrate the microwells by adding 200  $\mu$ l of TBST buffer (1x TBS, pH 8.0, containing 0.05% Tween-20) to every well. Incubate 15 minutes at room temperature. Tap the plate onto clean paper towels to remove liquid.
- 2) Thaw **S-adenosylmethionine** on ice. Upon first thaw, briefly spin tube containing **S-adenosylmethionine** to recover full contents of the tube. Aliquot **S-adenosylmethionine** into single use aliquots and store at  $-80^{\circ}\text{C}$ . *Note: S-adenosylmethionine is sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles.*
- 3) Add 30  $\mu$ l **0.5M DTT** (not provided) to **4x PRMT5 Assay Buffer**. Prepare the master mixture: N wells  $\times$  (7.5  $\mu$ l **4X PRMT5 buffer** + 0.5  $\mu$ l **100  $\mu$ M S-adenosylmethionine** + 17  $\mu$ l water). Add 25  $\mu$ l of master mixture to all wells labeled "Positive Control," "Test Sample," and "Blank." For wells labeled "Substrate control," add 7.5  $\mu$ l **4X PRMT5 buffer** + 17.5  $\mu$ l water.
- 4) Add 5  $\mu$ l of inhibitor solution of each well designated "Test Inhibitor". For the "Positive Control," "Substrate Control," and "Blank," add 5  $\mu$ l of 10% DMSO in water (inhibitor buffer).

	Blank	Substrate Control	Positive Control	Test Sample
4X PRMT5 buffer	7.5 $\mu$ l	7.5 $\mu$ l	7.5 $\mu$ l	7.5 $\mu$ l
100 $\mu$ M S-adenosylmethionine	0.5 $\mu$ l	–	0.5 $\mu$ l	0.5 $\mu$ l
H <sub>2</sub> O	17 $\mu$ l	17.5 $\mu$ l	17 $\mu$ l	17 $\mu$ l
Test Inhibitor	–	–	–	5 $\mu$ l
10% DMSO in water (Inhibitor buffer)	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l	–
1X PRMT5 buffer	20 $\mu$ l	–	–	–
Diluted PRMT5 (5 ng/ $\mu$ l)		20 $\mu$ l	20 $\mu$ l	20 $\mu$ l
<b>Total</b>	<b>50 <math>\mu</math>l</b>	<b>50 <math>\mu</math>l</b>	<b>50 <math>\mu</math>l</b>	<b>50 <math>\mu</math>l</b>

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court West, Ste. B  
San Diego CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

- 5) Prepare **1x PRMT5 Assay Buffer** by adding 1 part of **4x PRMT5 Assay Buffer** to 3 parts water (v/v).
- 6) Thaw **PRMT5/MEP50 enzyme** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full content of the tube. Aliquot **PRMT5/MEP50 enzyme** into single use aliquots. Store remaining undiluted enzyme in aliquots at -80°C immediately. *Note: PRMT5/MEP50 enzyme is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme.*
- 7) Dilute **PRMT5/MEP50 enzyme** in **1x PRMT5 buffer** at 5-10 ng/μl (100-200 ng/20 μl). Keep diluted enzyme on ice until use. Discard any unused diluted enzyme after use. *Note: Diluted enzyme may not be stable. Dilute the enzyme immediately before use.*
- 8) Add 20 μl of **1 x PRMT5 buffer** to the wells designated "Blank."
- 9) Initiate reaction by adding 20 μl of diluted **PRMT5/MEP50 enzyme** to the wells designated "Positive Control," "Substrate Control," and "Test Sample." Incubate at room temperature for 2 hours.
- 10) Remove the supernatant from the wells and wash the strip three times with 200 μl TBST buffer. Blot dry onto clean paper towels.
- 11) Add 100 μl of **Blocking buffer** to every well. Shake on a rotating platform for 10 min. Remove supernatant as above.

#### **Step 2:**

- 1) Dilute "**Primary antibody 4-3**" 100-fold with **Blocking buffer**.
- 2) Add 100 μl per well. Incubate 1 hour at room temperature with slow shaking.
- 3) Wash plate three times with 200 μl TBST buffer and incubate in **Blocking buffer** as in steps 1-10 and 1-11.

#### **Step 3:**

- 1) Dilute "**Secondary HRP-labeled antibody 2**" 1,000-fold with **Blocking buffer**.
- 2) Add 100 μl per well. Incubate for 30 min. at room temperature with slow shaking.
- 3) Wash plate three times with 200 μl TBST buffer and incubate in **Blocking buffer** as in steps 1-10 and 1-11.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**  
Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)  
Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)


- 4) Just before use, mix on ice 50  $\mu$ l **HRP chemiluminescent substrate A** and 50  $\mu$ l **HRP chemiluminescent substrate B** and add 100  $\mu$ l per well. Discard any unused chemiluminescent reagent after use.
- 5) Immediately read sample in a luminometer or microtiter-plate capable of reading chemiluminescence. "Blank" value is subtracted from all readings.

### Reading Chemiluminescence:

Chemiluminescence is the emission of light (luminescence) which results from a chemical reaction. The detection of chemiluminescence requires no wavelength selection because the method used is emission photometry and is not emission spectrophotometry.

To properly read chemiluminescence, make sure the plate reader is set for LUMINESCENCE mode. Typical integration time is 1 second, delay after plate movement is 100 msec. Do not use a filter when measuring light emission. Typical settings for the Synergy 2 BioTek plate reader are: use the "hole" position on the filter wheel; Optics position: Top; Read type: endpoint. Sensitivity may be adjusted based on the luminescence of a control assay without enzyme (typically we set this value as 100).

### Example of Assay Results:

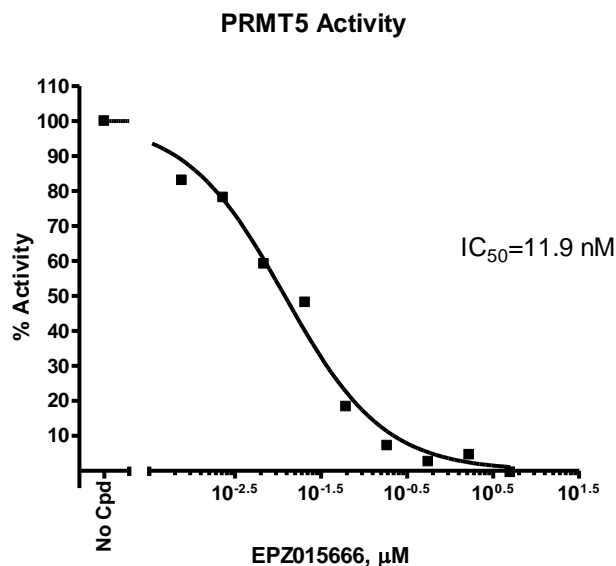

PRMT5 enzyme activity, measured using the *PRMT5 Chemiluminescent Assay Kit*, BPS Bioscience Catalog #52002L. Luminescence was measured using a Bio-Tek fluorescent microplate reader. *Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at [info@bpsbioscience.com](mailto:info@bpsbioscience.com).*

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court West, Ste. B  
San Diego CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

#### RELATED PRODUCTS

PRMT5/MEP50 (expressed in Sf9 cells)	#51048	20 µg
PRMT1 (expressed in E. coli)	#51040	50 µg
PRMT1 (expressed in Sf9 cells)	#51041	20 µg
PRMT3 (expressed in E. coli)	#51043	50 µg
PRMT4/CARM1 (expressed in HEK293)	#51047	20 µg
PRMT5/MEP50 (expressed in HEK293)	#51045	20 µg
PRMT6 (expressed in HEK293)	#51049	20 µg
PRMT7 (expressed in Sf9)	#51054	20 µg
PRMT8 (expressed in Sf9)	#51052	20 µg
PRMT9 (expressed in Sf9)	#51053	20 µg
PRMT1 Chemiluminescent Assay Kit	#52004L	96 reactions
PRMT3 Chemiluminescent Assay Kit	#52005L	96 reactions
PRMT4 Chemiluminescent Assay Kit	#52041L	96 reactions
PRMT6 Chemiluminescent Assay Kit	#52046	96 reactions
PRMT1 Homogeneous Assay Kit	#52054	384 reactions
PRMT3 Homogeneous Assay Kit	#52055	384 reactions
PRMT5 Homogeneous Assay Kit	#52052	384 reactions
PRMT6 Homogeneous Assay Kit	#52056	384 reactions
PRMT8 Homogeneous Assay Kit	#52058	384 reactions

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court West, Ste. B  
 San Diego CA 92121  
 Tel: 1.858.829.3082  
 Fax: 1.858.481.8694  
 Email: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

### TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Luminescence signal of positive control reaction is weak	PRMT5 enzyme has lost activity	Enzyme loses activity upon repeated freeze/thaw cycles. Use fresh enzyme (PRMT5, BPS Bioscience #51045). Store enzyme in single-use aliquots. Increase time of enzyme incubation. Increase enzyme concentration.
	Antibody reaction is insufficient	Increase time for primary antibody incubation. Avoid freeze/thaw cycles of antibodies.
	Incorrect settings on instruments	Record light signals at 5 second intervals. Refer to instrument instructions for settings to increase sensitivity of light detection. See "Reading Chemiluminescence" section.
	Chemiluminescent reagents mixed too soon	Chemiluminescent solution should be used within 15 minutes of mixing. Ensure both reagents are properly mixed.
Luminescent signal is erratic or varies widely among wells	Inaccurate pipetting/technique	Run duplicates of all reactions. Use a multichannel pipettor. Use master mixes to minimize errors.
	Bubbles in wells	Pipette slowly to avoid bubble formation. Tap plate lightly to disperse bubbles; be careful not to splash between wells.
Background (signal to noise ratio) is high	Insufficient washes	Increase number of washes. Increase wash volume. Increase Tween-20 concentration to 0.1% in TBST.
	Sample solvent is inhibiting the enzyme	Run negative control assay including solvent. Maintain DMSO level at <1% Increase time of enzyme incubation.
	Results are outside the linear range of the assay	Use different concentrations of enzyme (PRMT5, BPS Bioscience #51045) to create a standard curve.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)