



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court West, Suite B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet

ONE-Step™ Luciferase Assay System **Catalog #: 60690-5** **Size: 500 mL**

Description

The One-Step Luciferase Assay System is designed to be used for high-throughput, sensitive quantitation of firefly luciferase activity in mammalian cell culture. The reagent consists of two components, a Luciferase Reagent Buffer (**Component A**) and Luciferase Reagent Substrate (**Component B**). **Component A** and **Component B** are combined to form a working solution that contains all the necessary components for cell lysis and luciferase quantitation. This assay system has several features:

- Sensitive – highly sensitive detection of firefly luciferase activity.
- Stable – the signal output is stable for more than two hours, providing flexibility with regard to incubation time
- Convenient – simple one-step, homogeneous protocol.
- High-throughput – one-step homogeneous protocol minimizes handling steps to support high-throughput screening applications
- Compatibility – works well with a variety of common media containing 0-10% serum and phenol red.
- Instrumentation – does not require a luminometer with injectors.

Application

- Monitor firefly (*Photinus pyralis*) luciferase activity in cultured mammalian cells.
- High-throughput screening using luciferase reporter cell lines

Components

Component	Amount	Storage
Luciferase Reagent Buffer (Component A)	5 x 100 ml	-20°C
Luciferase Reagent Substrate, 100x (Component B)	5 x 1 ml	-20°C <i>Protect from light</i>

Each system contains sufficient reagents to perform 10,000 assays of 100 µl each in 96-well plate.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com


Stability

At least 6 months when stored as directed. Upon first thaw, store in aliquots at -20°C. The reagent may be subjected to several freeze/thaw cycles with no effect on functionality, but it is recommended that freeze/thaw cycles be avoided whenever possible.

Background

Luciferase is the general term given to a class of oxidative enzymes that catalyze reactions that give off light, a process known as bioluminescence (Fig. 1). In biology, researchers can take advantage of this reaction and use it as a readout for various biological processes. This has perhaps been exploited most in luciferase reporter cell lines where a promoter region from a gene of interest is placed immediately upstream of the coding sequence for luciferase. In this system, transcriptional activation of the gene of interest leads to a level of luciferase expression that is proportional to the level of gene activation.

The level of activation is then assessed by lysing the cells and adding the luciferase substrate, D-luciferin. Using the One-Step Luciferase Assay System from BPS Bioscience, this is accomplished in a single step as the lysis buffer has been formulated with D-luciferin. After a brief incubation period, bioluminescence is read on a luminometer.

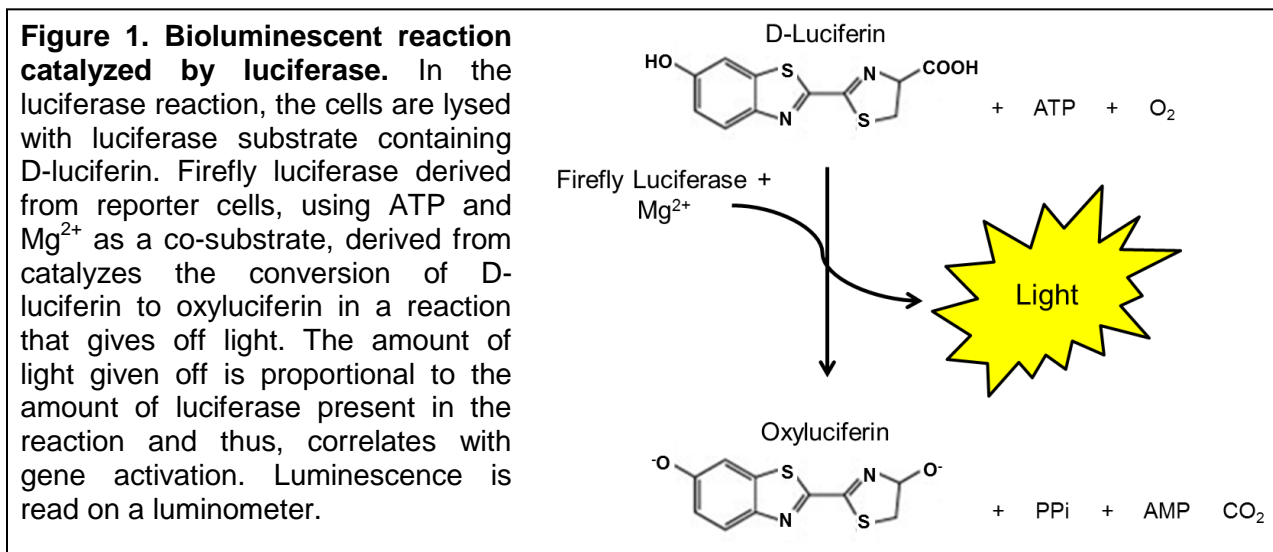

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court West, Suite B
San Diego, CA 92121

Tel: 1.858.829.3082

Fax: 1.858.481.8694

Email: info@bpsbioscience.com

Important Product Information

- The reagent has been validated in a 96-well format. Other formats will require scaling and optimization by the end-user.
- Luciferase Reagent Buffer must be at ~ room temperature (20- 25°C) before use.
- Luciferase Assay Reagent should be added to cells for at least 5 minutes before measuring luminescence to allow complete cell lysis.
- For maximal light intensity, measure samples within 1 hour of reagent addition.
- Avoid exposing to excessive heat or light during incubation.
- Different cell lines may exhibit variation in lysis ability and/or luminescence signal and may require slight optimization by the end-user.
- Luminescence signal is affected by assay conditions. Results should be compared between samples measured using the same cell type and media/serum combinations.
- To analyze multiple plates, include a common control sample in each plate and normalize the luminescence of each plate to the control contained in the same plate.
- Background luminescence is a characteristic of luminometer performance, therefore, background luminescence must be subtracted from all readings for accuracy.

Materials Required but Not Supplied

- Multiwell tissue culture plates that are compatible with luminometer being used
- Mammalian cells that express firefly luciferase
- Appropriate cell culture medium
- Laboratory platform shaker
- Luminometer

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**


Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com

General Assay Procedure



1. Thaw Luciferase Reagent Buffer (**Component A**) by placing the reagent in a room temperature water bath. Equilibrate the buffer to room temperature and mix well before use.
2. Calculate the amount of Luciferase Reagent needed for the experiment (**Component A + Component B**). Immediately prior to performing the experiment, prepare the luciferase assay working solution by diluting Luciferase Reagent Substrate (**Component B**) into Luciferase Reagent Buffer (**Component A**) at a 1:100 ratio and mix well. Avoid exposing to excessive light. *Only use enough of each component for the experiment, remaining **Component A** and **Component B** should be stored separately at -20°C.*
3. Remove multi-well plate containing mammalian cells from incubator. *Note: plates must be compatible with luminescence measurement with luminometer being used.*
4. Add equal volume of luciferase assay working solution (**Component A + Component B**) to the culture medium in each well. Example: 96-well plate with 100 µl of culture medium requires 100 µl of luciferase assay working solution per well.
5. Gently rock the plates for ≥15 minutes at room temperature. Measure firefly luminescence using a luminometer.

The signal under these conditions will be stable for more than 2 hours at room temperature. For maximal light intensity, measure samples within 1 hour of reagent addition.

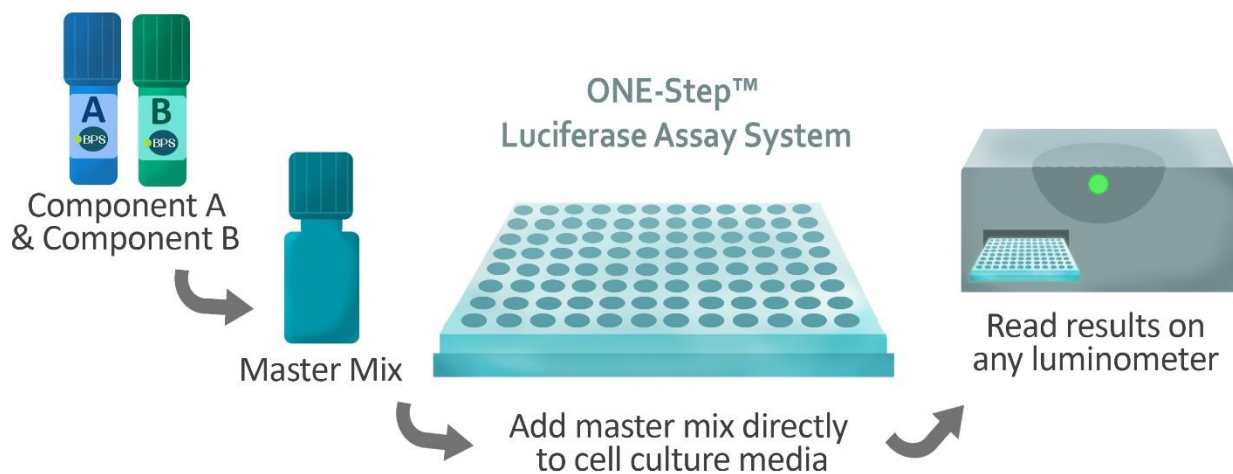

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com

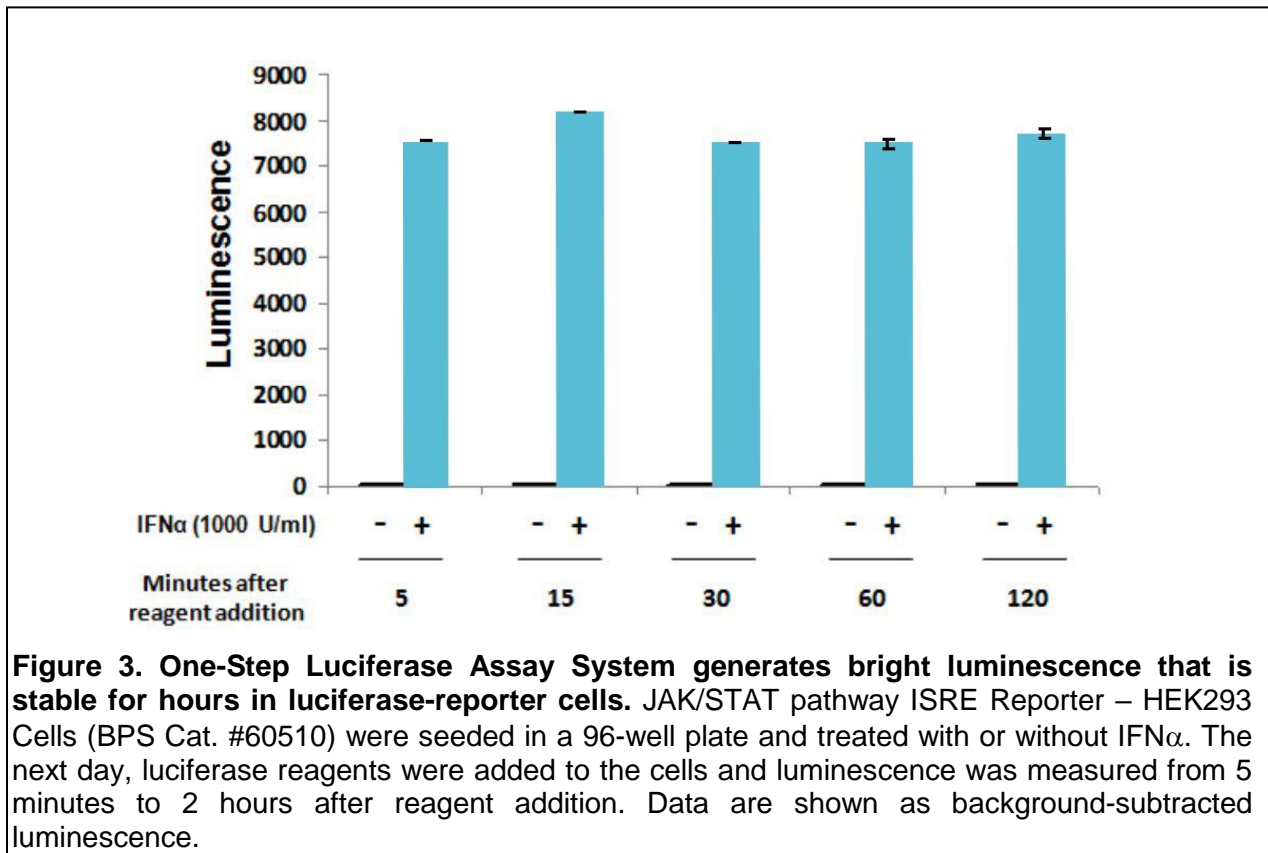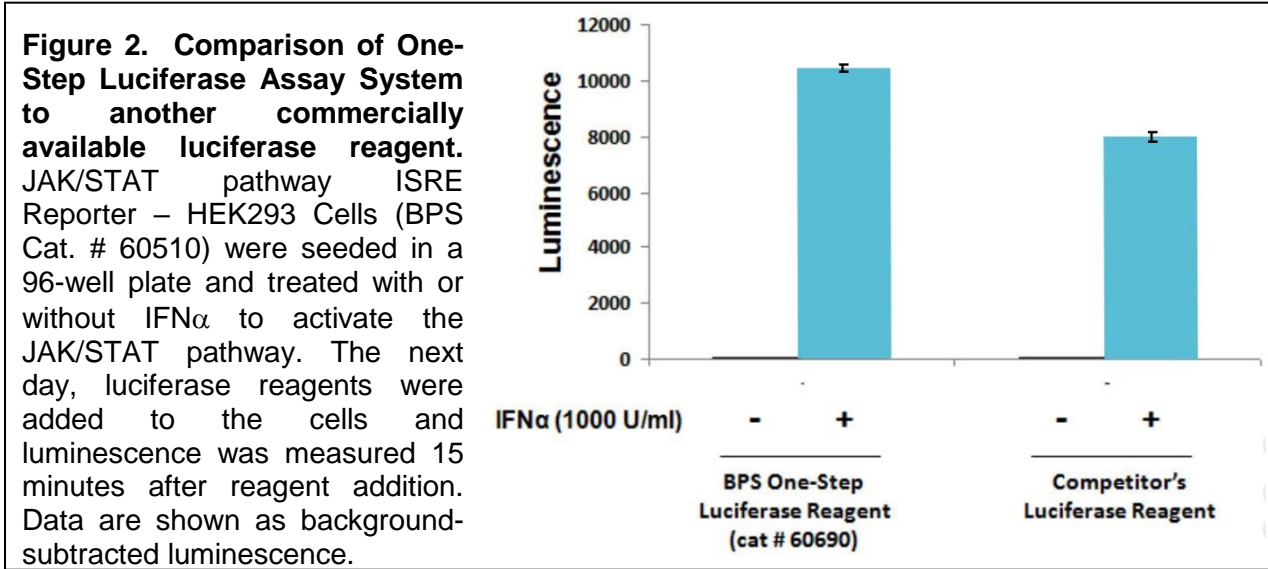

Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court West, Suite B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Related Products:

<u>Product</u>	<u>Cat. #</u>	<u>Size</u>
Myc Reporter (Luc) – HCT116 Cell Line	60520	2 vials
Myc Reporter (Luc) – HCT116 Cell Line	60406	2 vials
CRE/CREB Reporter (Luc) – HEK293 Cell Line	60515	2 vials
AP-1 Reporter – HEK293 Cell Line	60405	2 vials
ARE reporter cell line	60513	2 vials
GAL4 Reporter (Luc)-HEK293 Cell Line	60655	2 vials
NF- κ B reporter (Luc) - HEK293 Cell line	60650	2 vials
ISRE Reporter Cell Line	60510	2 vials
RAR α Reporter Cell Line	60503	2 vials
RAR β Reporter Cell Line	60603	2 vials
RAR γ Reporter Cell Line	60604	2 vials
Gli Reporter – NIH3T3 Cell Line	60409	2 vials
PDE7A/CRE Reporter - HEK293 Cell Line	60413	2 vials
SBE Reporter Kit	60654	500 rxns.
Myc Reporter Kit	60519	500 rxns.
ARE Reporter Kit	60514	500 rxns.
SBE Reporter Kit	60654	500 rxns.
SRE Reporter Kit	60511	500 rxns.
NF- κ B Reporter Kit	60614	500 rxns.
ISRE Reporter Kit	60613	500 rxns.
CRE/CREB Reporter Kit	60611	500 rxns.
FOXO Reporter Kit	60643	500 rxns.
AP1 Reporter Kit	60612	500 rxns.
Notch1/CSL Reporter Kit	60509	500 rxns.
GAL4 Reporter Kit	60522	500 rxns.
TCF/LEF Reporter Kit	60500	500 rxns.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com