



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Product Information

AccuGreen™ High Sensitivity dsDNA Quantitation Kit

For use with handheld fluorometers such as Qubit®

Kit Contents

| Component | 31066-T 100 assays | 31066 500 assays |
|--|-----------------------|-------------------------|
| AccuGreen™ Dye, 200X | 99817-T 250 uL | 99817 1.25 mL |
| DNA Quantitation Buffer, 20X Concentrate | 99979-3mL 3 mL | 99979-12.5mL 12.5 mL |
| AccuGreen™ Standard 1 (0 ng/uL) | 99819-T 1 mL | 99819 5 mL |
| AccuGreen™ Standard 2 (10 ng/uL) | 99820-T 1 mL | 99820 5 mL |

Reagents to be supplied by user

0.5 mL clear PCR tubes

Storage and Handling

Store kit at 4°C. Protect dye from light. The kit is stable for at least 12 months from date of receipt when stored as recommended. The dye color may change when stored at 4°C, but this will not affect performance. AccuGreen™ dye is stable for storage at room temperature for up to 6 months. AccuGreen™ dye is a potentially harmful chemical. Exercise universal laboratory safety precautions when handling the dye.

Spectral Properties

Ex/Em: 509/530 nm (bound to dsDNA). See Figure 1 for spectra.

Product Description

The AccuGreen™ High Sensitivity dsDNA Quantitation Kit is designed for use with handheld fluorometers such as the Qubit® fluorometers from Thermo Fisher. Unlike absorbance-based measurements, AccuGreen™ Dye is highly selective for double-stranded DNA over single-stranded DNA or RNA.

The linear range of the AccuGreen™ assay is 0.1 to 100 ng of DNA per tube. DNA samples with concentrations between 5 pg/uL and 100 ng/uL may be quantified using sample volumes between 1 and 20 uL (for example, 1 uL of 100 ng/uL is 100 ng total, and 20 uL of 5 pg/uL is 0.1 ng total, which will both fall within the linear range of the assay). If you use the most common sample volume of 10 uL, the starting sample concentration should be 10 pg/uL to 10 ng/uL.

The AccuGreen kits provide enough reagents to quantify approximately 100 samples (31066-T) or 500 samples (31066), plus the two standards. There are enough reagents for 250 reactions (31066-T) or 1,250 reactions (31066) including the standards.

Biotium also offers the AccuGreen™ High Sensitivity dsDNA Quantitation Solution (31068). While the AccuGreen™ Kit comes with a calf thymus DNA standard, the AccuGreen™ Solution does not provide a DNA standard. It is intended for those who wish to provide their own standard.

For high-sensitivity DNA quantitation using a fluorescent plate reader, we recommend using our AccuClear® or AccuBlue® dsDNA Quantitation Kits. The AccuClear® assay (31028) combines excellent sensitivity and the broadest linear range. The AccuBlue® NextGen assay (31060) provides unprecedented sensitivity and accuracy.

Protocol for reading the AccuGreen™ assay on the Qubit® Fluorometer

The first time that you use this kit, dilute the 20X buffer concentrate to 1X by adding dH₂O directly to the bottle containing the concentrated buffer. To 99979-3mL add 57 mL dH₂O; to 99979-12.5mL add 237.5 mL dH₂O. Invert the bottle several times to mix completely and mark the top of the cap to indicate that the buffer has been diluted to 1X.

This protocol describes how to measure AccuGreen™ fluorescence on a Qubit® 3.0 Fluorometer using the pre-programmed dsDNA High Sensitivity program. Instructions may vary for older Qubit® models.

Note: The linear range for this assay on the Qubit® 3.0 is 100 pg-100 ng DNA in the assay tube. However, samples even slightly below 100 pg will return the error message "Out of Range". Therefore for best results use more than 100 pg DNA per assay.

1. Warm all components to room temperature before use. AccuGreen™ Dye is provided in DMSO, which may freeze during storage at 4°C. You can place all kit components in a 37°C water bath for rapid warming; be sure to allow solutions to cool to room temperature before using.
2. Prepare the 1X working solution by diluting 200X AccuGreen™ Dye 1:200 in 1X DNA Quantitation Buffer. You will need 200 uL of working solution for each standard and sample to be tested. For example, to quantify three DNA samples, you will need 1 mL of working solution (2 standards and 3 samples).
3. For each sample and standard, pipette 190 uL of the working solution into a clear 0.5 mL PCR tube (if using the Qubit® fluorometer).
4. Into one tube, pipet 10 uL of AccuGreen™ Standard 1 (0 ng/uL).
5. Into a second tube, pipet 10 uL of AccuGreen™ Standard 2 (10 ng/uL).
6. Pipette 1-20 uL of each DNA sample to be quantified into its own tube.
7. Gently vortex or invert the tubes several times to mix, and then incubate the tubes at room temperature for at least 2 minutes.
8. Turn on the Qubit® 3.0 instrument. On the home screen select dsDNA. Choose the High Sensitivity assay.
9. Follow the prompts on the screen, and first read the AccuGreen™ Standard 1 and then the AccuGreen™ Standard 2. The program will use these values to quantify your unknown samples.
10. One at a time, measure each of your samples.
11. The data can be recorded manually or exported as a csv file.

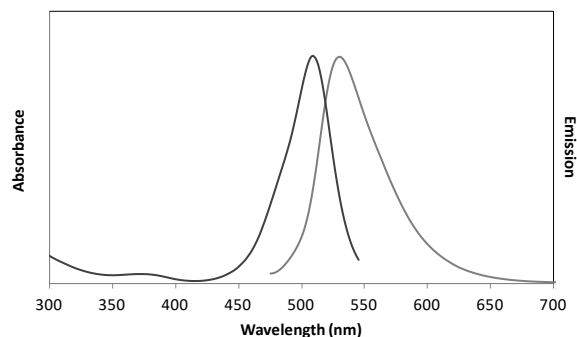


Figure 1. Absorbance and emission spectra of AccuGreen™ Dye bound to dsDNA.

Considerations for Data Analysis

Calf thymus DNA can serve as a reference for most plant and animal DNA because it is double-stranded, highly polymerized and is approximately 58% AT (42% GC). Lambda dsDNA yields similar results. You may wish to use a standard similar to your unknown samples in DNA length, structure (i.e., linear vs. circular), or GC content. For bacterial DNA, a species-specific standard may be desired because the GC content varies widely depending on the species.

The effects of common DNA contaminants such as salts, solvents, detergents and protein on the AccuGreen™ assay are listed in Table 1. The same contaminants were tested in the Qubit® HS DNA assay, and found to give essentially identical results.

Table 1. Effect of common DNA contaminants on AccuGreen™ assay signal

| Compound | Initial concentration in DNA sample | Final concentration in assay (200 uL) | Signal change (negative value indicates decrease) |
|--------------------|-------------------------------------|---------------------------------------|---|
| Sodium Chloride | 1 M | 50 mM | -13% |
| Magnesium Chloride | 100 mM | 5 mM | -21% |
| Sodium Acetate | 600 mM | 30 mM | -16% |
| Ammonium Acetate | 1 M | 50 mM | -14% |
| Ethanol | 20% | 1% | +9% |
| Phenol | 2% | 0.10% | +4% |
| SDS | 0.2% | 0.01% | -89% |
| SDS | 0.02% | 0.001% | -8% |
| Triton™ X-100 | 2% | 0.1% | +1% |
| Triton™ X-100 | 0.02% | 0.001% | +4% |
| Tween® 20 | 0.1% | 0.005% | +4% |
| CTAB | 0.01% | 0.0005% | -96% |
| BSA | 2 mg/mL | 0.1 mg/mL | no change |
| dNTPs | 2 mM | 100 uM | +2% |

Related Products

| Catalog number | Product |
|----------------|---|
| 31068 | AccuGreen™ High Sensitivity dsDNA Quantitation Solution |
| 31069 | AccuGreen™ Broad Range dsDNA Quantitation Kit |
| 31070 | AccuGreen™ Broad Range dsDNA Quantitation Solution |
| 31060 | AccuBlue® NextGen dsDNA Quantitation Kit |
| 31028 | AccuClear® Ultra High Sensitivity dsDNA Quantitation Kit with 7 DNA Standards |
| 31007 | AccuBlue® Broad Range dsDNA Quantitation Kit with 9 DNA Standards |
| 31073 | AccuBlue® Broad Range RNA Quantitation Kit |
| 41003 | GelRed® Nucleic Acid Gel Stain, 10,000X in water |
| 41005 | GelGreen® Nucleic Acid Gel Stain, 10,000X in water |
| 31041 | Forget-Me-Not™ EvaGreen® qPCR Master Mix (2-Color Tracking) |
| 31045-T | Forget-Me-Not™ EvaGreen® qPCR Master Mix (low ROX) |
| 31046-T | Forget-Me-Not™ EvaGreen® qPCR Master Mix (high ROX) |
| 31043 | Forget-Me-Not™ Universal Probe Master Mix |
| 22020 | 10X Phosphate-Buffered Saline (PBS) |
| 41024-4L | Water, Ultrapure Molecular Biology Grade |
| 31030 | DNA Gel Extraction Kit |
| CD504 | RNAstorm™ RNA Isolation Kit |
| CD501 | RNAstorm™ Kit for Isolation of RNA from FFPE Tissue Samples |
| CD502 | DNAsorm™ Kit for Isolation of RNA from FFPE Tissue Samples |

Please visit our website at www.biotium.com for information on our life science research products, including environmentally friendly EvaGreen® qPCR master mixes, fluorescent CF® dye antibody conjugates and reactive dyes, apoptosis reagents, fluorescent probes, and kits for cell biology research.

Qubit is a registered trademark of Thermo Fisher Scientific.

AccuBlue and AccuClear are registered trademarks of Biotium, Inc.

AccuGreen is a trademark of Biotium, Inc.

AccuBlue, AccuClear and AccuGreen are covered by granted and/or pending U.S. and/or international patents.

Materials from Biotium are sold for research use only, and are not intended for food, drug, household, or cosmetic use.