



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



CCL2/MCP-1 (human)
TR-FRET Biomarker Assay Kit

Powered by Bioauxilium's THUNDER™ TR-FRET Technology

Item No. 500237

www.caymanchem.com

Customer Service 800.364.9897

Technical Support 888.526.5351

1180 E. Ellsworth Rd • Ann Arbor, MI • USA

TABLE OF CONTENTS

| | | |
|------------------------------|----|-----------------------------------|
| GENERAL INFORMATION | 3 | Materials Supplied |
| | 3 | Safety Data |
| | 4 | Precautions |
| | 4 | Before You Start |
| | 4 | If You Have Problems |
| | 4 | Storage and Stability |
| | 5 | THUNDER™ General Information |
| | 6 | Materials Needed but Not Supplied |
| INTRODUCTION | 7 | Background |
| | 7 | About This Assay |
| | 8 | Principle Of This Assay |
| PRE-ASSAY PREPARATION | 9 | Reagent Preparation |
| | 13 | Sample Preparation |
| ASSAY PROTOCOL | 14 | Workflow |
| | 14 | Performing the Assay |
| | 15 | TR-FRET Plate Reader Settings |
| ANALYSIS | 16 | Calculations |
| | 17 | Performance Characteristics |
| RESOURCES | 25 | References |
| | 26 | Notes |
| | 27 | Warranty and Limitation of Remedy |

GENERAL INFORMATION

Materials Supplied

| Item Number | Item Name | 96 wells Quantity/Size | 480 wells Quantity/Size |
|-------------|--------------------------------------|---------------------------|----------------------------|
| 400290 | Europium-Labeled CCL2/MCP-1 Antibody | 1 vial/5 µl | 1 vial/25 µl |
| 400293 | Acceptor-Labeled CCL2/MCP-1 Antibody | 1 vial/20 µl | 1 vial/100 µl |
| 400275 | Assay Buffer 1 (5X) | 1 vial/1 ml | 3 vials/1 ml |

To request this kit in bulk quantities, please contact our Sales department. If any of the items listed above are damaged or missing, please contact our Customer Service department. Sales and Customer Service departments can be reached at (800) 364-9897 or (734) 971-3335. We cannot accept any returns without prior authorization.



WARNING: THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

Safety Data

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent *via* email to your institution.

Precautions

Please read these instructions carefully before beginning this assay.

Do not mix or substitute reagents or materials from other kit lots or kits. Kits are quality control tested as a set of components and performance cannot be guaranteed if utilized separately or substituted.

We cannot guarantee the performance of the product outside the conditions detailed in this kit booklet.

Before You Start

Please note the following:

ONLY white plates should be used for TR-FRET.

DO NOT modify the assay protocol or volumes.

If You Have Problems

Technical Service Contact Information

Phone: 888-526-5351 (USA and Canada only) or 734-975-3888

Email: techserv@caymanchem.com

In order for our staff to assist you quickly and efficiently, please be ready to supply the lot number of the kit (found on the outside of the box).

Storage and Stability

This kit will perform as specified if stored as directed at -80°C and used before the expiration date indicated on the outside of the box.

THUNDER™ General Information

THUNDER™ TR-FRET Biomarker Assay Kits are designed for the quantitative measurement of human cytokines in cell culture supernatants using homogeneous (no wash) TR-FRET technology.

THUNDER™ TR-FRET Biomarker Assay Kits are based on Bioauxilium's enhanced proprietary time-resolved Förster resonance energy transfer (TR-FRET) technology. THUNDER™ assays can be read on most commercially available TR-FRET-compatible plate readers (a list of suitable TR-FRET readers can be found at www.Bioauxilium.com). TR-FRET-based assays are homogeneous because they do not require any washing or separation steps. In addition, this THUNDER™ assay uses a standardized, simple, and rapid "add-incubate-measure" protocol with a two-step reagent addition. This streamlined assay protocol dramatically decreases hands-on time and provides a powerful alternative to cumbersome, error-prone and time-consuming techniques such as Western blot and ELISA.

THUNDER™ TR-FRET Biomarker Assay Kits do not require the addition of potassium fluoride to ensure signal stability. The non-toxic nature of THUNDER™ kits make them ideal candidates for multiplexing with other assay technologies using a sequential assay protocol.

THUNDER™ TR-FRET Biomarker Assay Kits contain the essential reagents necessary to carry out the measurement of human cytokines in cell culture supernatants, with the exception of the standard and the plate(s).

Materials Needed But Not Supplied

1. A plate reader equipped with a TR-FRET option
2. Adjustable pipettes; multichannel or repeating pipettor recommended
3. A source of ultrapure water, with a resistivity of 18.2 M Ω -cm and total organic carbon (TOC) levels of <10 ppb, is recommended. Pure water - glass-distilled or deionized - may not be acceptable. *NOTE: UltraPure Water is available for purchase from Cayman (Item No. 400000).*
4. Low-volume 384-well, or half-area 96-well, white plates
5. Recombinant human CCL2/MCP-1 standard, available from R&D Systems® (Catalog No. 279-MC)
6. Adhesive sealing film for plates

INTRODUCTION

Background

Chemokine (C-C motif) ligand 2 (CCL2), also known as monocyte chemoattractant protein-1 (MCP-1), is a chemokine and member of the γ chemokine family.¹⁻³ CCL2 is expressed by a variety of cells, including epithelial, myeloid, smooth muscle, and endothelial cells, as well as fibroblasts, either constitutively or after induction by inflammatory stimuli in a cell-type dependent manner.³ It binds to chemokine (C-C motif) receptor 2 (CCR2), inducing downstream signaling via the JAK2/STAT3, MAPK, and PI3K pathways to promote cell migration and phospholipase C-mediated calcium release.² CCL2 functions as a chemoattractant for monocytes, basophils, T cells, and natural killer (NK) cells and has roles in various immune processes, including basophil and mast cell degranulation and respiratory burst, myeloid-derived suppressor cell-induced immunosuppression, and macrophage lipid body biogenesis, endocytosis, phagocytosis, and efferocytosis. CCL2 levels are increased in tumor microenvironments and positively correlate with tumor cell proliferation, migration, and invasion.³ Brain levels of CCL2 are increased in patients with Parkinson's disease, Alzheimer's disease, and amyotrophic lateral sclerosis (ALS) and are associated with pathogenesis.⁴

About This Assay

This CCL2/MCP-1 (human) TR-FRET Biomarker Assay Kit uses a homogenous TR-FRET sandwich immunoassay method amenable to rapid measurement of CCL2/MCP-1 levels in cell culture supernatants. The standard curve spans a dynamic range of 13-10,000 pg/ml and a lower limit of quantification (LLOQ) of 61 pg/ml using the assay buffer included in the kit with laser excitation. The amount of reagents provided is sufficient for testing either 96 or 480 wells, depending on the size of the kit.

Principle Of This Assay

This assay is a homogenous TR-FRET-based sandwich immunoassay (Figure 1, below). Native or recombinant human CCL2/MCP-1 in cell culture supernatant samples is detected with a pair of fluorophore-labeled antibodies reactive to human CCL2/MCP-1.

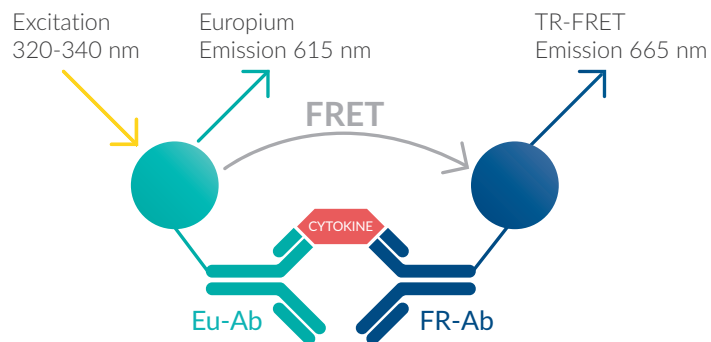


Figure 1. Schematic of the TR-FRET biomarker assay principle

The first antibody is labeled with a long-lifetime donor fluorophore (a europium chelate; Eu-Ab) and the second with a far-red acceptor fluorophore (FR-Ab). The binding of the two labeled antibodies to distinct epitopes on the target cytokine takes place in solution and brings the two dyes into close proximity. Excitation of the donor Eu chelate molecules with a flash lamp (320 or 340 nm) or laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which, in turn, emit a TR-FRET signal at 665 nm. The signal at 665 nm is proportional to the concentration of human CCL2/MCP-1 in the cell culture supernatant. Residual energy from the Eu chelate generates light at 615 nm.

TR-FRET assays exhibit very low background fluorescence levels and high signal-to-background (S/B) ratios. The data can be expressed and analyzed as either the signal at 665 nm or the 665 nm/615 nm ratio. The ratiometric measurement further increases assay reproducibility and robustness.

PRE-ASSAY PREPARATION

Reagent Preparation

The instructions described below are for testing the entire number of wells in each kit. Adjust volumes accordingly when testing fewer wells. Prepare only as much reagent as is needed on the day of the assay.

Bring all reagents to room temperature prior to use.

Centrifuge all tubes before use to improve recovery of content (2,000 xg, 10-15 seconds).

Use ultrapure water (18 M Ω -cm) to dilute the Assay Buffer 1 (5X).

NOTE: It is recommended to test all standards in triplicate and samples at least in duplicate.

1. Assay Buffer 1 (1X)

Mix the Assay Buffer 1 (5X) (Item No. 400275) end-over-end before use. The thawed Assay Buffer 1 (5X) will be stable for 1-2 weeks when stored at 4°C or for longer periods of time when stored at -80°C.

Add 1 ml (96 wells) or 3 ml (480 wells) of Assay Buffer 1 (5X) with 4 ml (96 wells) or 12 ml (480 wells) of ultrapure water. Store unused Assay Buffer 1 (1X) at 4°C; it will be stable for approximately two days.

2. CCL2/MCP-1 Standard Serial Dilutions

Recombinant human CCL2/MCP-1 is not included in the kit. The kit has been validated using recombinant human CCL2/MCP-1 from R&D Systems® (Catalog No. 279-MC). Other commercial recombinant human CCL2/MCP-1 proteins may also work but have not been tested with the kit.

Prepare standard solutions just before use and use within one hour. Do NOT store the standard solutions.

Each well requires 15 µl of serially diluted standards.

Reconstitute the analyte following the manufacturer's instructions (typically 100-200 µg/ml) to create a stock solution. Dilute the stock solution in either Assay Buffer 1 (1X) or culture medium containing 10% FBS, depending on which will be used for the standard curve dilutions, to create a 1 µg/ml CCL2/MCP-1 working stock solution. The working stock solution will be used to prepare a standard dilution series.

NOTE: If all of the working stock solution will not be used within one hour, aliquot into screw-capped polypropylene vials and store at -80°C for long-term storage. The undiluted stock solution may also be stored at -80°C.

To prepare the standard solutions for use: obtain 10 clean test tubes and label them #1-10. Alternatively, serial dilutions can be conducted using a 96-well polypropylene plate.

Prepare serial dilutions for the standard curve as indicated in Table 1, on page 11, using either Assay Buffer 1 (1X) or culture medium containing 10% FBS. Change the pipette tip between each standard dilution. Mix each tube thoroughly between each transfer.

| Tube | Volume of CCL2/MCP-1 Standard | Volume of Assay Buffer 1 (1X) OR Culture Medium + 10% FBS (µl) | CCL/MCP-1 in standard curve (g/ml in 15 µl) |
|------|--------------------------------------|--|---|
| 1 | 10 µl diluted CCL2/MCP-1 (0.1 µg/ml) | 90 | 1.0E-8 |
| 2 | 30 µl of tube 1 | 70 | 3.0E-9 |
| 3 | 30 µl of tube 2 | 60 | 1.0E-9 |
| 4 | 30 µl of tube 3 | 70 | 3.0E-10 |
| 5 | 30 µl of tube 4 | 60 | 1.0E-10 |
| 6 | 30 µl of tube 5 | 70 | 3.0E-11 |
| 7 | 30 µl of tube 6 | 60 | 1.0E-11 |
| 8 | 30 µl of tube 7 | 70 | 3.0E-12 |
| 9 | 30 µl of tube 8 | 60 | 1.0E-12 |
| 10* | 0 | 70 | 0 |

Table 1. Standard curve serial dilutions

*In order to calculate the limit of detection (LOD), three additional zero-standard tubes per plate should be prepared and assayed in triplicate.

3. Antibody Detection Mix (4X)

NOTE: Due to the low reagent volumes, the antibodies are diluted with Assay Buffer 1 (1X) directly in the vial when assaying only 96 wells.

NOTE: Each well requires 2.5 µl of Europium-Labeled CCL2/MCP-1 Antibody (Eu-Ab) (8X) and 2.5 µl of Acceptor-Labeled CCL2/MCP-1 Antibody (FR-Ab) (8X).

Prepare and mix just before use.

For 96 wells:

Eu-Ab Solution (8X): Add 255 µl of Assay Buffer 1 (1X) into the vial containing 5 µl of Europium-Labeled CCL2/MCP-1 Antibody (Eu-Ab) (Item No. 400290), and gently mix.

FR-Ab Solution (8X): Add 240 µl of Assay Buffer 1 (1X) into the vial containing 20 µl of Acceptor-Labeled CCL2/MCP-1 Antibody (FR-Ab) (Item No. 400293), and gently mix.

OR

For 480 wells:

Eu-Ab Solution (8X): Gently mix 25 µl of Eu-Ab (Item No. 400290) with 1,275 µl of Assay Buffer 1 (1X).

FR-Ab Solution (8X): Gently mix 100 µl of FR-Ab (Item No. 400293) with 1,200 µl of Assay Buffer 1 (1X).

Store unused Eu-Ab Solution (8X) and FR-Ab Solution (8X) at 4°C; they will be stable for approximately two days. If all unused Eu-Ab Solution (8X) and FR-Ab Solution (8X) will not be used within two days, aliquot and store at -80°C. Avoid repeated freeze/thaw cycles.

Sample Preparation

Each well requires 15 µl of sample (cell culture supernatant).

After collection, place the samples at 4°C and assay immediately or aliquot and store at -80°C.

Avoid repeated freeze/thaw cycles. Thaw samples at room temperature. Do NOT use a water bath to thaw samples.

Prepare serial dilutions using the same media used to culture the cells. If the analyte concentration in the sample possibly exceeds the highest point of the standard curve, prepare one or more 10-fold dilutions of the sample.

ASSAY PROTOCOL

Workflow

The THUNDER™ TR-FRET Cell Signaling Assay workflow consists of a two-step assay reagent addition, with an incubation step following each step.



Figure 2. Assay workflow

Performing the Assay

Samples and standards MUST be assayed at least in duplicate (triplicate recommended) each time the assay is performed.

When loading reagents in the low-volume 384-well microplate, change pipette tips between each standard or sample addition and after each set of reagents.

If using a multichannel pipettor, always use a new disposable reagent reservoir.

When reagents are added to the microplate, ensure the liquids are at the bottom of the well.

1. Add 15 µl of each working standard or 15 µl of the supernatant sample to the designated wells.
2. Add 2.5 µl of Eu-Ab Solution (8X) to each of the assay wells.
3. Cover the plate with adhesive plate-sealing film and incubate for 30 minutes at room temperature.
4. Add 2.5 µl of FR-Ab Solution (8X) to each of the assay wells.
5. Cover the plate with adhesive plate-sealing film and incubate for 1 hour at room temperature.
6. Gently remove the adhesive plate-sealing film. Read the plate on a TR-FRET-compatible microplate reader (see Table 2, page 15).

NOTE: The same plate can be read several times without a negative effect on the assay performance.

TR-FRET Plate Reader Settings

We recommend reading the TR-FRET assays at two wavelengths, detecting both the emission from the Eu chelate donor fluorophore at 615 nm and the acceptor fluorophore at 665 nm. Table 2, below, provides instrument settings to be used as guidelines.

| Parameter | TR-FRET-compatible Plate Reader | |
|---------------------------|---------------------------------|--------------------|
| | Flash lamp excitation | Laser excitation |
| Excitation filter | 320 nm (or 340 nm) | N/A |
| Emission filter | 615 nm (or 620 nm) | 615 nm (or 620 nm) |
| Delay time | 90 µs | 50 µs |
| Flash energy level | 100% or High | 100% |
| Number of flashes | 100-200 | 20 |
| Window (integration time) | 300 µs | 100 µs |

Table 2. Recommended TR-FRET plate reader settings

Calculations

Because the TR-FRET signal is read in a time-resolved manner, background subtraction is usually not necessary.

Do not use a value of zero for the zero standard. Use the actual background values obtained with the zero standard.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

1. TR-FRET data are typically calculated and presented ratiometrically using the following formula:

$$[(665 \text{ nm}/615 \text{ nm}) \times 1,000]$$

2. Calculate the TR-FRET ratio for each well. Alternatively, the signals at 665 nm can be used directly to analyze the data.
3. Create a standard curve by analyzing the data according to a nonlinear regression using the four-parameter logistic equation (sigmoidal dose-response curve with variable slope) and a $1/Y^2$ data weighting. For correct analysis, the highest standard point should be removed from the curve if it is present after the hook point (i.e., if it shows lower counts compared to the next standard).

Performance Characteristics

Representative Data

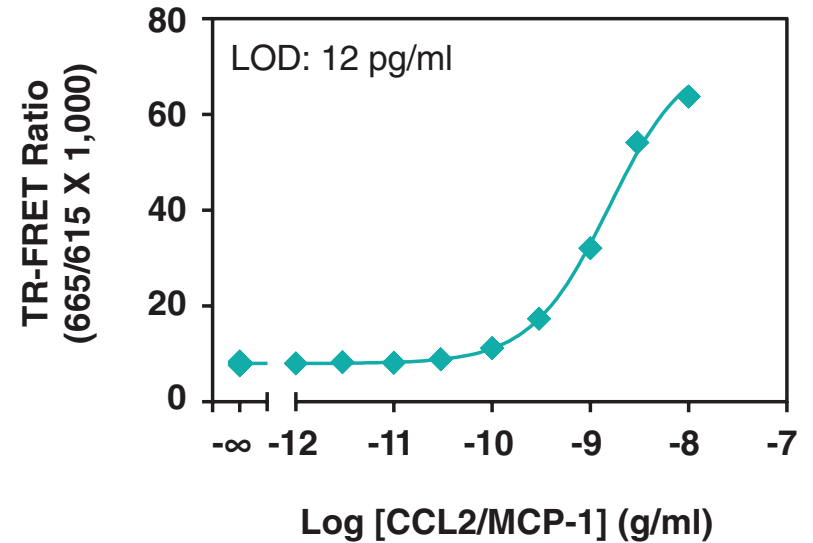


Figure 3. Typical standard curve using lamp excitation Data represent the mean \pm standard deviation of triplicate measurements (three wells) for each standard dilution. The data was generated using a low-volume 384-well white plate read with an EnVision® Multilabel Plate Reader equipped with a TR-FRET option.

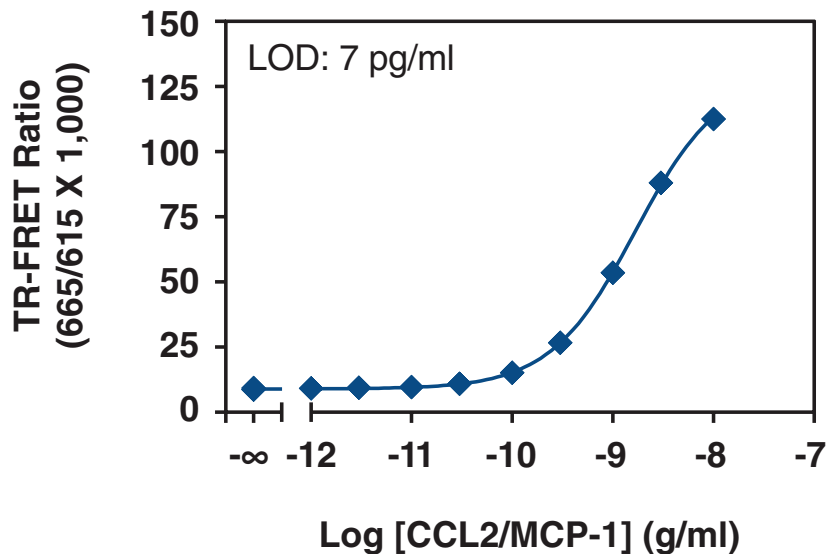


Figure 4. Typical standard curve using laser excitation Data represent the mean \pm standard deviation of triplicate measurements (three wells) for each standard dilution. The data was generated using a low-volume 384-well white plate read with an EnVision® Multilabel Plate Reader equipped with a TR-FRET option.

Sensitivity

The limit of detection (LOD) was calculated by adding two standard deviations to the mean background counts of 12 zero-standard replicates and calculating the corresponding concentration on the standard curve.

The lower limit of quantification (LLOQ) was calculated by adding 10 standard deviations to the mean background counts of 12 zero-standard replicates and calculating the corresponding concentration on the standard curve.

| | Assay Buffer 1 (1X) | | DMEM + 10% FBS | | RPMI + 10% FBS | |
|-----------------------------|---------------------|------|----------------|------|----------------|------|
| | Laser | Lamp | Laser | Lamp | Laser | Lamp |
| # of Standard Curves | 12 | 12 | 3 | 3 | 3 | 3 |
| Standard Curves LOD (pg/ml) | 13 | 20 | 18 | 20 | 11 | 16 |
| LLOQ (pg/ml) | 61 | 103 | 85 | 110 | 72 | 96 |

Table 3. Sensitivity of the CCL2/MCP-1 (human) TR-FRET Assay Kit

Precision

Intra-assay precision:

Three samples of known concentration (low, medium, high) were tested twenty times each on the same plate to assess intra-assay precision. The assay was run in a 384-well plate with the standards and samples prepared in Assay Buffer 1 (1X). The resulting low %CV (<10%) for each sample indicates good reproducibility within an assay.

| | Sample Concentration | | |
|--------------------|----------------------|--------|-------|
| | Low | Medium | High |
| # of Samples | 20 | 20 | 20 |
| Mean (pg/ml) | 353 | 1,403 | 2,908 |
| SD (pg/ml) | 27 | 57 | 223 |
| CV (%) Intra-Assay | 7.6 | 4.1 | 7.7 |

Table 4. Intra-assay precision

Inter-assay precision:

Three samples of known concentration were tested in separate assays to assess inter-assay precision. A total of seven independent experiments were performed by two operators using two different kit lots: The assay was run in 384-well plates with the standards and samples prepared in Assay Buffer 1 (1X) and tested in triplicate. The resulting low %CV (<10%) for each sample indicates good reproducibility between assays.

| | Sample Concentration | | |
|--------------------|----------------------|--------|-------|
| | Low | Medium | High |
| # of Samples | 25 | 25 | 25 |
| Mean (pg/ml) | 305 | 1,108 | 3,255 |
| SD (pg/ml) | 27 | 85 | 245 |
| CV (%) Intra-Assay | 8.9 | 7.6 | 7.5 |

Table 5. Inter-assay precision

Recovery

To assess the accuracy of the assay, spike-to-recovery experiments were conducted using Assay Buffer 1 (1X) and culture media. Each assay consisted of one standard curve and three spike concentrations (low, medium, high) assayed in triplicate in Assay Buffer 1 (1X), DMEM (supplemented with 10% FBS), or RPMI (supplemented with 10% FBS). The percent recovery is reported as the mean concentration of nine wells determined using a standard curve prepared in the corresponding diluent. The data shows recovery values within the range of 80-120% demonstrating assay accuracy.

| | Spike Level | Expected (pg/ml) | Observed (pg/ml) | % Recovery |
|---------------------|-------------|------------------|------------------|------------|
| Assay Buffer 1 (1X) | High | 3,000 | 3,008 | 100 |
| | Medium | 1,000 | 1,034 | 103 |
| | Low | 300 | 261 | 87 |
| DMEM + 10% FBS | High | 3,000 | 3,197 | 107 |
| | Medium | 1,000 | 1,000 | 100 |
| | Low | 300 | 253 | 84 |
| RPMI + 10% FBS | High | 3,000 | 3,254 | 108 |
| | Medium | 1,000 | 915 | 92 |
| | Low | 300 | 262 | 87 |

Table 6. Recovery in assay buffer or culture media

Linearity

To assess the linearity of the assay, a high concentration of CCL2/MCP-1 was spiked into culture medium supplemented with 10% FBS and diluted with the same culture medium with 10% FBS. The assay was run in a 384-well plate with standards prepared in the same type of supplemented culture medium. Each sample was measured in triplicate. Data shows recovery values within 80-120% confirming assay linearity.

| | Expected (pg/ml) | Observed (pg/ml) | % Recovery |
|----------------|------------------|------------------|------------|
| DMEM + 10% FBS | 3,125 | 3,484 | 111 |
| | 1,563 | 1,850 | 118 |
| | 781 | 924 | 118 |
| | 391 | 451 | 115 |
| | 195 | 200 | 103 |
| | 98 | 109 | 112 |
| RPMI + 10% FBS | 3,125 | 3,093 | 99 |
| | 1,563 | 1,667 | 107 |
| | 781 | 853 | 109 |
| | 391 | 404 | 103 |
| | 195 | 200 | 103 |

Table 7. Dilutional linearity in culture medium

| Assay Summary | | |
|--------------------------------------|---|-----------------|
| Assay Type | Homogeneous sandwich immunoassay | |
| Format | 384-well plate | |
| Hands-on Time | 15 minutes | |
| Incubation Time | 1.5 hours | |
| Sample Type (Volume) | Cell culture supernatant (15 µl) | |
| Specificity | Native and recombinant human CCL2/MCP-1 | |
| | Laser excitation | Lamp excitation |
| Limit of Detection (LOD) | 13 pg/ml | 20 pg/ml |
| Lower Limit of Quantification (LLOQ) | 61 pg/ml | 103 pg/ml |
| EC ₅₀ Value | 1.6 ng/ml | 1.8 ng/ml |
| Dynamic Range | 13-10,000 pg/ml | 20-10,000 pg/ml |

Table 8. Assay validation summary for the CCL2/MCP-1 (human) TR-FRET Biomarker Assay Kit

RESOURCES

References

1. Deshmane, S.L., Kremlev, S., Amini, S., *et al.* Monocyte chemoattractant protein-1 (MCP-1): An overview. *J. Interferon. Cytokine Res.* **29(6)**, 313-326 (2009).
2. Gschwandtner, M., Derler, R., and Midwood, K.S. More than just attractive: How CCL2 influences myeloid cell behavior beyond chemotaxis. *Front. Immunol.* **10**, 2759 (2009).
3. Lim, S.Y., Yuzhalin, A.E., Gordon-Weeks, A.N., *et al.* Targeting the CCL2-CCR2 signaling axis in cancer metastasis. *Oncotarget* **7(19)**, 28697-28710 (2016).
4. O'Connor, T., Borsig, L., and Heikenwalder, M. CCL2-CCR2 signaling in disease pathogenesis. *Endocr. Metab. Immune Disord. Drug Targets* **15(2)**, 105-118 (2015).

NOTES

Warranty and Limitation of Remedy

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

This document is copyrighted. All rights are reserved. This document may not, in whole or part, be copied, photocopied, reproduced, translated, or reduced to any electronic medium or machine-readable form without prior consent, in writing, from Cayman Chemical Company.

©10/26/2021, Cayman Chemical Company, Ann Arbor, MI, All rights reserved.
Printed in U.S.A.

