



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com


[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

CRF™ Anti-Polyvalent HRP Polymer

Description: CRF™ Anti-Polyvalent HRP Polymer has been developed to provide the cleanest, most consistent staining available. Developed in the research laboratories of ScyTek, the system is based on a polymerized peroxidase label that eliminates biotin and its' associated background issues from the equation. In addition, this product reduces the steps required for immunohistochemical staining by combining two steps from the traditional Biotin-Streptavidin system. CRF™ Anti-Polyvalent HRP Polymer is effective with antibodies of mouse, rat, rabbit and guinea pig.

Species of Origin:	Goat
Antigen Specificity:	Anti-Polyvalent (Mouse, Rat, Rabbit and Guinea Pig)
Preadsorbed Against:	Human
Enzyme Conjugate:	Peroxidase

Uses/Limitations: Not to be taken internally.
For In-Vitro Diagnostic use only.
Histological applications.
Do not use if reagents become cloudy.
Do not use past expiration date.
Use caution when handling reagents.
Non-Sterile.

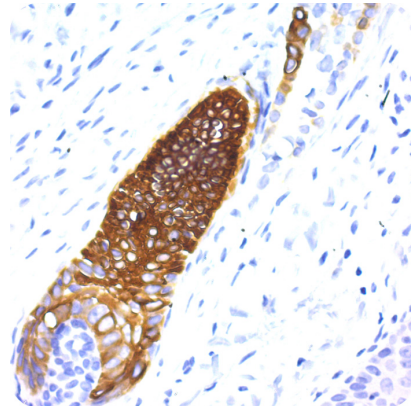

Control Tissue: Any well-fixed tissue section.
Frozen tissue section.
Cytocentrifuge preparation.

Ordering Information and Current Pricing at www.scytek.com


Availability:	<u>Item #</u>	<u>Volume</u>
	ABZ008	8 ml
	ABZ015	15 ml
	ABZ125	125 ml
	ABZ500	500 ml
	ABZ999	1000 ml

Recommended, But Not Included:


<u>Item #</u>	<u>Description</u>
CPL500	Citrate Plus
ADA500	Peroxide Block for Image
AAA500	Super Block
ACT500	DAB Chromogen/Substrate Kit (High Contrast)
HMM500	Hematoxylin, Mayer's (Lillie's Modification)
BRT500	Bluing Reagent

Storage: Store at 2-8°C.

Storage: 2° C  8° C

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.


Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands

Precautions: Avoid contact with skin and eyes.
Harmful if swallowed.
Follow all Federal, State, and local regulations regarding disposal.

Procedure:

1. Rehydrate tissue slides.
2. In a glass or plastic (Autoclavable) Coplin jar, add 5 ml of Citrate Plus (CPL) and 45 ml of deionized water.
3. Submerge slides in diluted Citrate Plus and loosely cap.
4. Add Distilled water to bottom of Autoclave or Pressure Cooker (about 1 inch deep in Pressure Cooker).
5. Place Coplin jar in Pressure Cooker or Autoclave.
6. Turn heat on and allow pressure to rise to 20-25 PSI.
7. Maintain pressure at 20-25 PSI for 5 minutes.
8. Turn off heat source and allow to cool.
9. When pressure has dropped to ambient, carefully remove lid or open door.
10. Using tongs, remove Coplin Jar and place on counter.
11. Once Coplin Jar cools to room temperature remove slides, rinse several times in buffer and proceed with staining as usual.
12. Apply Peroxide Block for Image Analysis (ADA) and incubate slide for 10-15 minutes.
13. Rinse 3 times in buffer.
14. Apply Super Block (AAA), and incubate for 5 minutes at room temperature to block nonspecific background staining. **Note:** Do not exceed 10 minutes or there may be a reduction in desired stain.
15. Rinse 3 times in buffer.
16. Apply primary antibody and incubate according to manufacturer's protocol.
17. Rinse 3 times in buffer.
18. Apply CRF™ Anti-Polyvalent HRP Polymer and incubate for 30 minutes at room temperature.
19. Rinse 3 times deionized water.


WARNING: DAB is a suspected carcinogen. Handle with care and dispose of according to all regulations.

20. Combine 50-80µl of DAB Chromogen with 1ml of DAB Substrate and apply to tissue for 5 minutes.

NOTE: Range = 50-80µl of DAB Chromogen / 1 ml of DAB Substrate. 50µl /ml produces optimal staining quality, 80µl/ml produces maximum sensitivity. Combined mixture may be used for up to six hours.


21. Rinse 1 time in deionized water.
22. Apply DAB Chromogen/Substrate mixture and incubate for a second 5-minute period.
23. Rinse 3 times in deionized water.
24. Apply Hematoxylin, Mayer's (HMM) and incubate for 1 minute.
25. Rinse 3 times in distilled water.
26. Apply Bluing Reagent (BRT) and incubate for 5-10 seconds.
27. Rinse immediately in distilled or deionized water.
28. Dehydrate slides and clear in xylene or xylene substitute.

Storage: 2° C  8° C



ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.

CE 


Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands

29. Coverslip using a permanent mounting media.

-Troubleshooting Guide-

Overstaining:

1. Concentration of the primary antibody was too high or the incubation time was too long.
2. Temperature during incubation was too high.
3. Incubation times were too long.

Non-Specific Background Staining:

1. Rinsing between steps was inadequate.
2. Tissue was allowed to dry with reagents on.
3. Folds in tissue trapped reagents.
4. Antigen migrated in tissue.
5. Excessive tissue adhesive on slides.
6. Inadequate blocking with protein block.


Weak Staining:

1. Primary antibody concentration was too low or incubation time was too short.
2. Reagents are past their expiration date.
3. Inadequate removal of wash buffer between steps, resulting in dilution of reagents.
4. Room temperature was excessively cool.
5. The primary antibody does not recognize an antigen that survives fixation and embedding in high enough amounts.
6. Excessive incubation with protein block (Super Block or normal serum).


No Staining:

1. Steps were inadvertently left out.
2. There is no antigen in the tissue.
3. The primary antibody is not of mouse, rat, rabbit or guinea pig origin.
4. Chromogenic substrate has been replaced with another that is not intended for use with peroxidase.
5. One or more components of the kit have been inactivated.

Storage: 2° C  8° C

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.


Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands