



SZABO
SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic

VectaFluor™ Duet Double Labeling Kit

DyLight™ 594 Anti-Rabbit IgG, DyLight 488 Anti-Mouse IgG



Cat. No. **DK-8828**

Storage Store reagents in original bottles at 2-8 °C. Do not freeze.


Description Instructions for immunofluorescent staining.

VectaFluor Duet Double Labeling Kit offers maximum convenience to achieve double label immunofluorescence staining. This kit can detect mouse and rabbit primary antibodies with green and red DyLight fluorescent dyes in one step and is intended for use on non-rodent specimens.

Kit Components

| Product Name | Volume |
|---|--------|
| Normal Horse Serum, 2.5% | 15 ml |
| VectaFluor Duet Reagent [DyLight 594 Anti-Rabbit IgG (Red), DyLight 488 Anti-Mouse IgG (Green)] | 15 ml |

The VectaFluor Double Labeling Kit will stain approximately 150 sections based on 100 µl per section.

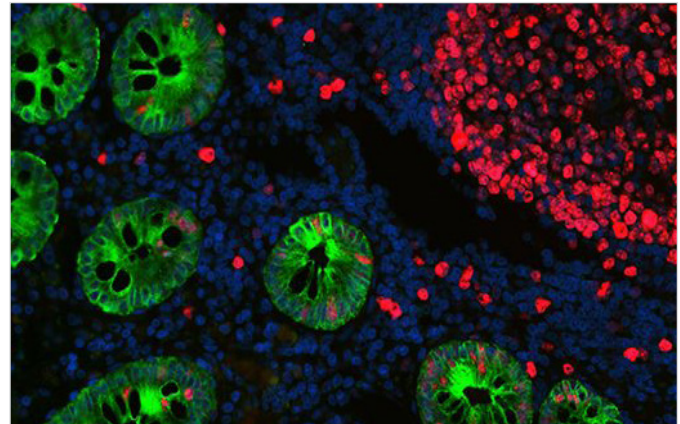

Colon: Cytokeratin (AE1/AE3) (m) and Ki67 (rb) detected simultaneously with VectaFluor Duet IF Double Labeling Kit, DyLight 488 Anti-Mouse (green)/DyLight 594 Anti-Rabbit (red). Mounted in VECTASHIELD® HardSet™ Antifade Mounting Medium.

Preparation of Working Solution

VectaFluor Double Labeling Kit reagents are ready-to-use—no mixing or titering is necessary to obtain optimal staining.

The staining procedure should be performed at room temperature (20–25°C). VectaFluor Double Labeling Kit reagents should be equilibrated to room temperature for optimal performance.

A number of different wash buffers can be used. One of the most common is 10 mM sodium phosphate, pH 7.5, 0.9% saline (PBS). 0.1% Tween® 20 detergent may be added to the wash buffer and is especially recommended for use with automated stainers.

Staining Procedure

- For paraffin sections, deparaffinize and hydrate tissue sections through xylenes or other clearing agents and graded alcohol series. For frozen sections or cell preparations fix with acetone or an appropriate fixative for the antigen under study, if required.
- If antigen unmasking is required, perform this procedure using a Vector® Antigen Unmasking Solution, Citrate-based (H-3300) or Tris-based (H-3301).
- Wash in buffer for 5 minutes.

- Incubate for 20 minutes with 2.5% Normal Horse Serum (or blocking solution of choice).
- Tip off excess serum from sections.
- Incubate with mouse primary antibody diluted in an appropriate diluent.
- Wash in buffer for 5 minutes.
- Incubate with rabbit primary antibody diluted in an appropriate diluent.
- Wash in buffer for 2 x 5 minutes.
- Incubate for 30 minutes with VectaFluor Duet Reagent.
- Wash for 2 x 5 minutes in buffer.
- Mount in a media suitable for fluorescence, such as one of the VECTASHIELD Antifade Mounting Media.

Detailed product listing, specifications, protocols and additional information are available on our website: vectorlabs.com

For research use only. Not intended for animal or human therapeutic or diagnostic use.

©2022 Vector Laboratories. All rights reserved.

HardSet, VectaFluor, VECTASHIELD, and Vector are trademarks of Vector Laboratories, Inc., DyLight is a trademark of Pierce Biotechnology Inc. Tween 20 is a trademark of Atlas Chemical Industries.
DK-8828.LBL-02266.Rev.00