



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



ZYMO RESEARCH

THE  
Epigenetics  
COMPANY™

## **EZ-96 DNA Methylation™ Kit** (Shallow-Well Format)

For high-throughput bisulfite conversion of DNA

### Highlights

- High throughput (96-well), proven procedure for bisulfite conversion of DNA.
- 96-well desulphonation and recovery of bisulfite-treated DNA.
- Recovered DNA is ideal for downstream analyses including PCR, endonuclease digestion, sequencing, microarrays, etc.

Catalog Numbers:  
D5003



Scan with your smart-phone camera to  
view the online protocol/video.



[tech@zymoresearch.com](mailto:tech@zymoresearch.com)



[www.zymoresearch.com](http://www.zymoresearch.com)



Toll Free: (888) 882-9682

# Table of Contents

---

<b>Product Contents</b> .....	<b>01</b>
<b>Introduction to DNA Methylation</b> .....	<b>02</b>
<b>Specifications</b> .....	<b>04</b>
<b>Product Description</b> .....	<b>05</b>
<b>Protocol</b> .....	<b>07</b>
Reagent Preparation .....	<b>07</b>
Sample Processing .....	<b>08</b>
<b>Appendix</b> .....	<b>10</b>
<b>Frequently Asked Questions</b> .....	<b>12</b>
<b>Ordering Information</b> .....	<b>13</b>
<b>Guarantee</b> .....	<b>15</b>

# Product Contents

EZ-96 DNA Methylation™ Kit	D5003 (2 x 96 Rxns.)	Storage Temperature
CT Conversion Reagent*	2 Bottles	Room Temp.
M-Dilution Buffer	5.2 ml	Room Temp.
M-Binding Buffer	80 ml	Room Temp.
M-Wash Buffer**	2 x 36 ml	Room Temp.
M-Desulphonation Buffer	40 ml	Room Temp.
M-Elution Buffer	8 ml	Room Temp.
Silicon-A™ Binding Plate	2 Plates	Room Temp.
Conversion Plates w/ Pierceable Cover Film	2 Plates/2 Films	Room Temp.
Collection Plates	2 Plates	Room Temp.
Elution Plates	2 Plates	Room Temp.
Instruction Manual	1	-

\* 7.5 ml water and 2.1 ml **M-Dilution Buffer** are added per bottle of **CT Conversion Reagent** and mixed prior to use.


\*\* Add 144 ml of 100% ethanol to the 36 ml **M-Wash Buffer** concentrate before use.

# Introduction to DNA Methylation

DNA methylation is a naturally occurring event in both prokaryotic and eukaryotic organisms. In prokaryotes DNA methylation provides a way to protect host DNA from digestion by restriction endonucleases that are designed to eliminate foreign DNA, and in higher eukaryotes DNA methylation functions in the regulation/control of gene expression (1). It has been demonstrated that aberrant DNA methylation is a widespread phenomenon in cancer and may be among the earliest changes to occur during oncogenesis (2). DNA methylation has also been shown to play a central role in gene imprinting, embryonic development, X-chromosome gene silencing, and cell cycle regulation. In many plants and animals, DNA methylation consists of the addition of a methyl group to the fifth carbon position of the cytosine pyrimidine ring via a methyltransferase enzyme (3). The majority of DNA methylation in mammals occurs in 5'-CpG-3' dinucleotides, but other methylation patterns do exist. In fact, about 80 percent of all 5'-CpG-3' dinucleotides in mammalian genomes are found to be methylated, whereas the majority of the twenty percent that remain unmethylated are within promoters or in the first exons of genes.

The ability to detect and quantify DNA methylation efficiently and accurately has become essential for the study of cancer, gene expression, genetic diseases, as well as many other important aspects of biology. To date, a number of methods have been developed to detect/quantify DNA methylation including: high-performance capillary electrophoresis (4) and methylation-sensitive arbitrarily primed PCR (5). However, the most common technique used today remains the bisulfite conversion method (6). This technique involves treating methylated DNA with bisulfite, which converts unmethylated cytosines into uracil. Methylated cytosines remain unchanged during the treatment. Once converted, the methylation profile of the DNA can be determined by PCR amplification followed by DNA sequencing (see next page).

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	
Original DNA with methylated C <sup>m</sup> pG	▶	G	T	T	G	C <sup>m</sup> G	C	T	C	A	C	T	G	C	C	
DNA Sequencing after CT conversion	▶	G	T	T	G	C	G	T	T	T	A	T	T	G	T	T

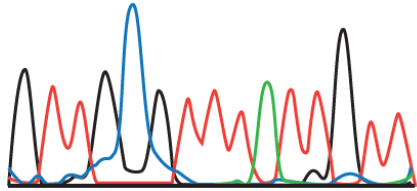

**DNA sequencing results following bisulfite treatment.** DNA with methylated CmpG at nucleotide position #5 was processed using the **EZ DNA Methylation™ Kit**. The recovered DNA was amplified by PCR and then sequenced directly. The methylated cytosine at position #5 remained intact while the unmethylated cytosines at positions #7, 9, 11, 14 and 15 were completely converted into uracil following bisulfite treatment and detected as thymine following PCR.

## References:

1. Costello JF, Plass CJ. *Med. Genet.* 2001; 38(5): 285-303.
2. Stirzaker C. *Cancer Res.* 1997; 57(11): 2229-2237.
3. Adams RL. *Bioessays.* 1995; 17(2): 139-145.
4. Fraga MF, *et al.* *Electrophoresis.* 2000; 21(14): 2990-2994.
5. Gonzalgo ML. *Cancer Res.* 1997; 57(4): 594-599.
6. Frommer M. *Proc. Natl. Acad. Sci. USA.* 1992; 89(5): 1827-1831.


# Specifications

- **DNA Input** – Samples containing 500 pg - 2 µg of DNA. For optimal results, the amount of input DNA should be from 200 to 500 ng.
- **Conversion Efficiency** – > 99% of non-methylated C residues are converted to U; > 99% protection of methylated cytosines.
- **DNA Recovery** – > 80%

# Product Description

The **EZ-96 DNA Methylation™ Kit** features a high-throughput (96-well), simplified procedure that streamlines bisulfite conversion of DNA. The kit is based on the three-step reaction that takes place between cytosine and sodium bisulfite where cytosine is converted into uracil. The product's innovative desulphonation technology eliminates otherwise cumbersome precipitations. The kit is designed to reduce template degradation, minimize DNA loss during treatment and clean-up, while ensuring complete conversion of the DNA. Purified, converted DNA is ideal for PCR amplification for downstream analyses including endonuclease digestion, sequencing, microarrays, etc.

An outline comparing the **EZ DNA Methylation™ Kit** procedure to Zymo Research's other methylation kits is shown below.

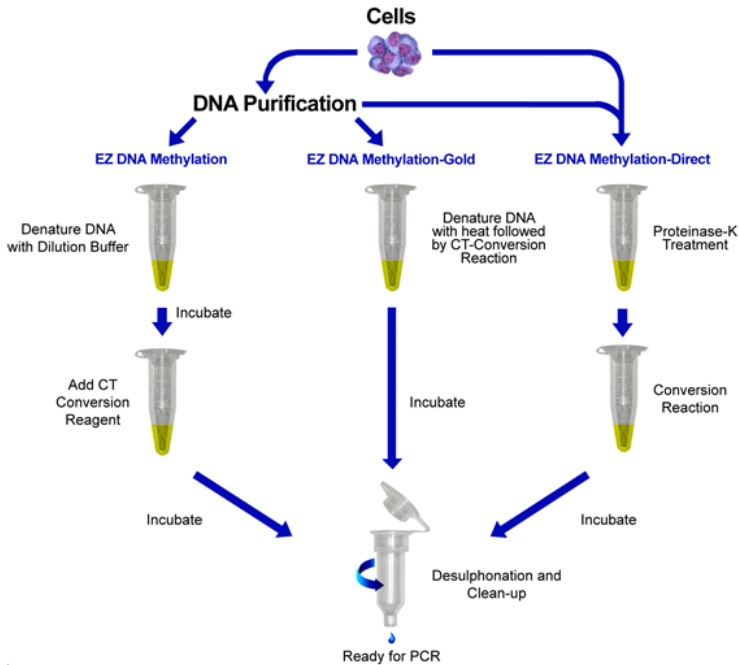

Outline of the **EZ DNA Methylation™**, **EZ DNA Methylation-Gold™** and **EZ DNA Methylation-Direct™** Kit procedures.



**Note:** Single spin-column formats are available for processing smaller numbers of samples. Also, MagPrep kits are available (p. 13) for adaptation to liquid handling robots (e.g., Tecan – Freedom EVO®) and automated sample prep.

**Selected EZ DNA Methylation™ Kit Citations:**

1. Ehrich M, et al. *Nuc. Acids Res.* 2007; 35 (5): e29
2. Kaneda M, et al. *Nature.* 2004; 429: 900-903
3. Zhang F, et al. *Proc. Natl. Acad. Sci. USA.* 2007; 104 (11): 4395-4400.
4. Oda M, et al. *Genes & Dev.* 2006; 20: 3382-3394.
5. England RPM, et al. *Nature Meth.* 2005; 2: 1-2.

# Protocol

## Reagent Preparation

### Preparation of CT Conversion Reagent

- ✓ The **CT Conversion Reagent** supplied within this kit is a solid mixture and must be prepared prior to first use. Prepare as follows:
  1. Add 7.5 ml water and 2.1 ml of **M-Dilution Buffer** to a bottle of **CT Conversion Reagent**.
  2. Mix at room temperature with frequent vortexing or shaking for 10 minutes.

**Note:** It is normal to see trace amounts of undissolved reagent in the **CT Conversion Reagent**. Each bottle of **CT Conversion Reagent** is designed for 96 separate DNA treatments.

**Storage:** The **CT Conversion Reagent** is light sensitive, so minimize its exposure to light. For best results, the **CT Conversion Reagent** should be used immediately following preparation. If not used immediately, the **CT Conversion Reagent** solution can be stored overnight at room temperature, one week at 4°C, or up to one month at -20°C. Stored **CT Conversion Reagent** solution must be warmed to 37°C, then vortexed prior to use.

### Preparation of M-Wash Buffer

- ✓ Add 144 ml of 100% ethanol to the 36 ml **M-Wash Buffer** concentrate before use.

## Sample Processing

1. Add 5  $\mu\text{l}$  of **M-Dilution Buffer** to each DNA sample in a **Conversion Plate** and adjust the total volume to 50  $\mu\text{l}$  with water. Mix each sample by pipetting up and down.

**Example:** For 14  $\mu\text{l}$  of a DNA sample add 5  $\mu\text{l}$  M-Dilution Buffer and 31  $\mu\text{l}$  water.

2. Incubate the **Conversion Plate** containing the samples at 37°C for 15 minutes.
3. After the above incubation, add 100  $\mu\text{l}$  of the prepared **CT Conversion Reagent** to each sample and mix.
4. Incubate the **Conversion Plate** in the dark at 50°C for 12-16 hours (e.g., using a thermal cycler).

The **CT Conversion Reagent** is light sensitive, so try to minimize the reaction's exposure to light whenever possible.

5. Incubate the sample at 0-4°C (e.g., on ice or using a thermal cycler) for 10 minutes. Samples may be kept at 4°C for up to 20 hours.
6. Add 400  $\mu\text{l}$  of **M-Binding Buffer** to each well of a **Silicon-A™ Binding Plate** on a **Collection Plate**.

**Note:** The capacity of each well of the Binding Plate is 600  $\mu\text{l}$ . The capacity of each well of the Collection Plate is 800  $\mu\text{l}$ . Empty the Collection Plate whenever necessary to prevent contamination of the Binding Plate contents by the flow-through.

7. Load the samples (from Step 5) into the wells of the **Silicon-A™ Binding Plate** containing the **M-Binding Buffer**. Mix by pipetting up and down.
8. Centrifuge at  $\geq 3,000 \times g$  (5,000  $\times g$  max.) for 5 minutes. Discard the flow-through.
9. Add 400  $\mu\text{l}$  of **M-Wash Buffer** to each well and centrifuge at  $\geq 3,000 \times g$  for 5 minutes.
10. Add 200  $\mu\text{l}$  of **M-Desulphonation Buffer** to each well of and let stand at room temperature (20-30°C) for 15-20 minutes. After the incubation, centrifuge at  $\geq 3,000 \times g$  for 5 minutes.

11. Add 500  $\mu\text{l}$  of **M-Wash Buffer** to each well and centrifuge at  $\geq 3,000 \times g$  for 5 minutes. Add another 500  $\mu\text{l}$  of **M-Wash Buffer** and centrifuge for 10 minutes.
12. Place the **Silicon-A™ Binding Plate** onto an **Elution Plate**. Add **30  $\mu\text{l}$  of M-Elution Buffer** directly to the binding matrix in each well. Centrifuge for 3 minutes at  $\geq 3,000 \times g$  to elute the DNA.

**Note:** Alternatively, water or TE (pH  $\geq 6.0$ ) can be used for elution if required for your experiments.

The DNA is ready for immediate analysis or can be stored at or below  $-20^{\circ}\text{C}$  for later use. For long term storage, store at or below  $-70^{\circ}\text{C}$ . We recommend using 1-4  $\mu\text{l}$  of eluted DNA for each PCR, however, up to 30  $\mu\text{l}$  can be used if necessary. The elution volume can be  $> 30 \mu\text{l}$  depending on the requirements of your experiments, but small elution volumes will yield more concentrated DNA.



Zymo Research provides primer design assistance with its Bisulfite Primer Seeker Program, available at:

[www.zymoresearch.com/tools/bisulfite-primer-seeker](http://www.zymoresearch.com/tools/bisulfite-primer-seeker)

- 4. Amount of DNA Required for Bisulfite Conversion.** The minimal amount of human or mouse genomic DNA required for bisulfite treatment and subsequent PCR amplification is 100 µg. The optimal amount of DNA per bisulfite treatment is 200 to 500 ng. Although, up to 2 µg of DNA can be processed, it should be noted that high input levels of DNA may result in incomplete bisulfite conversion for some GC-rich regions.
- 5. PCR Conditions.** Usually, 35 to 40 cycles are required for successful PCR amplification of bisulfite converted DNA. Optimal amplicon size should be between 150-300 bp; however larger amplicons (up to 1 kb) can be generated by optimizing the PCR conditions. Annealing temperatures between 55-60°C typically work well.

As most non-methylated cytosine residues are converted into uracil, the bisulfite-treated DNA usually is AT-rich and has low GC composition. Non-specific PCR amplification is relatively common with bisulfite treated DNA due to its AT-rich nature. PCR using “hot start” polymerases is strongly recommended for the amplification of bisulfite-treated DNA.

**Note:** **ZymoTaq™** is a “hot start” DNA polymerase specifically designed for the amplification of bisulfite treated DNA.

## Frequently Asked Questions

**Q: Should the input DNA be dissolved in TE, water, or some other buffer prior to its conversion?**



*A: Water, TE or modified TE buffers can be used to dissolve the DNA and do not interfere with the conversion process.*

**Q: Which Taq polymerase(s) do you recommend for PCR amplification of converted DNA?**

*A: We recommend a “hot start” DNA polymerase (e.g., ZymoTaq™ DNA Polymerase).*

**Q: Why are there two different catalog numbers for the EZ-96 DNA Methylation™ Kit?**

*A: The two different catalog numbers are used to differentiate between the binding plates that are included in the kit. Deep and shallow-well binding plates are available to accommodate most rotors and microplate carriers. Below is a comparison of the two binding plates.*

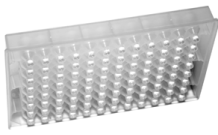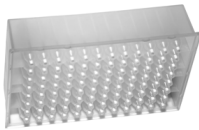		
<b>Binding Plate</b>	<b>Silicon-A™ Plate</b>	<b>Zymo-Spin™ I-96 Plate</b>
<b>Style</b>	Shallow-Well	Deep-Well
<b>Height of Binding Plate</b>	19 mm (0.75 inches)	35 mm (1.38 inches)
<b>Binding Plate/Collection Plate Assembly</b>	43 mm (1.69 inches)	60 mm (2.36 inches)
<b>Binding Cap./Minimum Elution Volume</b>	5 µg/30 µl	5 µg/15 µl
<b>Catalog Numbers</b>	<b>D5003</b>	<b>D5004</b>

# Ordering Information

Product Description	Catalog No.	Size
EZ DNA Methylation™ Kit	D5001 D5002	50 Rxns. 200 Rxns.
EZ-96 DNA Methylation™ Kit (Shallow-Well)	D5003	2 x 96 Rxns.
EZ-96 DNA Methylation™ Kit (Deep-Well)	D5004	2 x 96 Rxns.
EZ-96 DNA Methylation™ MagPrep*	D5040 D5041	4 x 96 Rxns. 8 x 96 Rxns.

\* MagPrep kits are adaptable to liquid handling robots (e.g., Tecan – Freedom EVO®) making them ideal for automated sample prep.

Individual Kit Components	Catalog No.	Amount
CT Conversion Reagent	D5001-1 D5003-1	1 Tube 1 Bottle
M-Dilution Buffer	D5001-2 D5002-2	1.3 ml 5.2 ml
M-Binding Buffer	D5005-3 D5006-3 D5040-3	30 ml 125 ml 250 ml
M-Wash Buffer	D5001-4 D5002-4 D5007-4 D5040-4	6 ml 24 ml 36 ml 72 ml
M-Desulphonation Buffer	D5001-5 D5002-5 D5040-5	10 ml 40 ml 80 ml
M-Elution Buffer	D5001-6 D5002-6 D5007-6 D5041-6	1 ml 4 ml 8 ml 40 ml
Zymo-Spin™ IC Columns (capped)	C1004-50 C1004-250	50 Pack 250 Pack
Collection Tubes	C1001-50 C1001-500 C1001-1000	50 Pack 500 Pack 1,000 Pack



<b>MagBinding Beads</b>	D4100-5-3 D4100-5-8 D4100-5-16	6 ml 8 ml 16 ml
<b>Zymo-Spin™ I-96 Binding Plates</b>	C2004	2 Plates
<b>Silicon-A™ Binding Plates</b>	C2001	2 Plates
<b>Conversion Plates w/ Pierceable Cover Film</b>	C2005	2 Plates /Films
<b>Collection Plates</b>	C2002	2 Plates
<b>Elution Plates</b>	C2003	2 Plates



**100% satisfaction guarantee on all Zymo Research products,  
or your money back.**

Zymo Research is committed to simplifying your research with quality products and services. If you are dissatisfied with this product for any reason, please call 1(888) 882-9682.

---

Integrity of kit components is guaranteed for up to one year from date of purchase. Reagents are routinely tested on a lot-to-lot basis to ensure they provide the highest performance and reliability.

This product is for research use only and should only be used by trained professionals. It is not for use in diagnostic procedures. Some reagents included with this kit are irritants. Wear protective gloves and eye protection. Follow the safety guidelines and rules enacted by your research institution or facility.

Use of Methylation Specific PCR (MSP) is protected by US Patents 5,786,146 & 6,017,704 & 6,200,756 & 6,265,171 and International Patent WO 97/46705. No license under these patents to use the MSP process is conveyed expressly or by implication to the purchaser by the purchase of this product.

™ Trademarks of Zymo Research Corporation



ZYMO RESEARCH

The **BEAUTY** of **SCIENCE** is to Make Things **SIMPLE**®



[tech@zymoresearch.com](mailto:tech@zymoresearch.com)



[www.zymoresearch.com](http://www.zymoresearch.com)



Toll Free: (888) 882-9682