

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH


Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


MCAK Antibody

Rabbit Polyclo					
Antigen Affinity Purified			Protein ID	NP_006836.1	
Catalog No. A300–807A			GenelD	11004	BETHYL
Lot No.	A300-80	17A-2 Laboratories, inc			
APPLICATIONS		WB, IP			
SPECIES REACTIVITY		Human			
AMOUNT		100 µl			
CONCENTRATION		1000 μg/ml			
STORAGE/SHELF LIFE		2 - 8° C / 1 year from date of receipt			
PHYSICAL STATE		Liquid			
BUFFER		Tris-citrate/phosphate buffer, pH 7 to 8 containing 0.09% Sodium Azide			
ISOTYPE		IgG			
ORIGIN		USA			
PRODUCTION PROCEDURES		Antibody was affinity purified using an epitope specific to MCAK immobilized on solid support.			
		The epitope recognized by A300-807A maps to a region between residue 550 and 600 of human Mitotic Centromere-Associated Kinesin using the numbering given in entry NP_006836.1 (GeneID 11004).			
		Immunoglobulin concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.			
APPLICATIONS		Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.			
		Western Blot	1:2	,000 - 1:10,000	
		Immunoprecipit	ation 2 –	5 µg/mg lysate	
APPLICATION N	(Western blot of immunoprecipitates performed using Normal Pig Serum (Cat. No. S100–020), Goat anti-Rabbit Light Chain HRP Conjugate (Cat. No. A120–113P) and 4–8% SDS–PAGE (link to IP-western blot protocol in Additional Info section below).			
	Ŋ	Western blot of l	ysates perfor	med using standard western b	lot reagents and 4-8% SDS-PAGE.
ADDITIONAL IN	NFO	https://www.bethyl.com/product/A300-807A			
				S, a current list of citations, ar s://www.bethyl.com/content/	nd other product specific information. protocol_IP_WB

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer Date: June 21, 2019

For in vitro laboratory use only. Not for any clinical, therapeutic or diagnostic use in humans or animals. Not for human or animal consumption. This product may not be resold or modified for resale without the prior written approval of Bethyl Laboratories, Inc. The information provided in this data sheet is believed to be correct but does not purport to be all-inclusive and is intended to be used as a guide. Bethyl Laboratories, Inc. shall not be liable or responsible in any way for use of either this information or the material supplied. Disposal of hazardous material may be subject to federal, state or local laws or regulations.

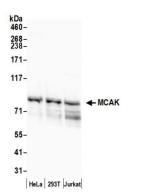

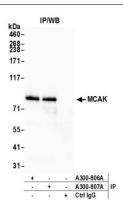

Detection of human MCAK by western blot. Samples: Whole cell lysate (50 μ g) from HeLa, HEK293T, and Jurkat cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-MCAK antibody A300-807A (lot A300-807A-2) used for WB at 0.1 μ g/ml. Detection: Chemiluminescence with an exposure time of 30 seconds. Detection of human MCAK by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1.0 mg per IP reaction; 20% of IP loaded) from HeLa cells prepared using NETN lysis buffer. *Antibodies:* Affinity purified rabbit anti-MCAK antibody A300-807A (lot A300-807A-2) used for IP at 3 µg per reaction. MCAK was also immunoprecipitated by rabbit anti-MCAK antibody A300-806A. For blotting immunoprecipitated MCAK, A300-807A was used at 1 µg/ml. *Detection:*

Chemiluminescence with an exposure time of 10 seconds.

Bethyl Laboratories, Inc. • 25043 West FM 1097 • Montgomery, TX 77356 • 800.338.9579 • 936.597.6111 • 866.597.6105 (FAX) • www.bethyl.com • technical@bethyl.com