

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


ATRX Antibody

| Rabbit Polyclo | nal | | | | | the second second |
|---------------------------|-------|---|------------|---|--|-----------------------------|
| Antigen Affinity Purified | | | Protein ID | NP_000480.2 | | |
| Catalog No. A301–045A–T | | GenelD | 546 | | | |
| Lot No. | A301- | 045A-T-4 | | | | BEINYL LABORATORIES, INC |
| APPLICATIONS | | WB, IP, IHC | | | | |
| SPECIES REACTIVITY | | Human | | | | |
| PRESUMED REACTIVITY | | Based on 100% sequence identity, this antibody is predicted to react with Chimpanzee | | | | |
| AMOUNT | | 20 µl (2 blots) | | | | |
| CONCENTRATION | | 100 µg/ml | | | | |
| STORAGE/SHELF LIFE | | 2 - 8° C / 1 year from date of receipt | | | | |
| PHYSICAL STATE | | Liquid | | | | |
| BUFFER | | Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide | | | | |
| ISOTYPE | | IgG | | | | |
| ORIGIN | | USA | | | | |
| PRODUCTION PROCEDURES | | Antibody was affinity purified using an epitope specific to ATRX immobilized on solid support. | | | | |
| | | The epitope recognized by A301-045A-T maps to a region between residue 2442 and the C- terminus (residue 2492) of human alpha thalassemia/mental retardation syndrome X-linked using the numbering given in entry NP_000480.2 (GeneID 546). | | | | |
| APPLICATIONS | | Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use. | | | | |
| | | Western Blot | 1:1 | 1000 | | |
| | | Immunoprecipit | us | e in immunoprecipitation | hin A301-045A-T has been qualified for n; however, we recommend using the his antibody found as product A301-045A. | |
| | | Immunohistoch | , | 50 to 1:200. Epitope retri commended for FFPE tiss | | fer pH6.0 is |
| APPLICATION NOTES | | Western blot of immunoprecipitates performed using Normal Pig Serum (Cat. No. S100–020), Goat anti–Rabbit Light Chain HRP Conjugate (Cat. No. A120–113P) and 4–8% SDS–PAGE (link to IP–western blot protocol in Additional Info section below). | | | | |
| ADDITIONAL INFO | | Western blot of lysates performed using standard western blot reagents and 4-8% SDS-PAGE. https://www.bethyl.com/product/A301-045A-T Use the link above to view SDS, a current list of citations, and other product specific information. IP-western blot protocol https://www.bethyl.com/content/protocol_IP_WB | | | | |

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer Date: October 30, 2019

For in vitro laboratory use only. Not for any clinical, therapeutic or diagnostic use in humans or animals. Not for human or animal consumption. This product may not be resold or modified for resale without the prior written approval of Bethyl Laboratories, Inc. The information provided in this data sheet is believed to be correct but does not purport to be all-inclusive and is intended to be used as a guide. Bethyl Laboratories, Inc. shall not be liable or responsible in any way for use of either this information or the material supplied. Disposal of hazardous material may be subject to federal, state or local laws or regulations.

2

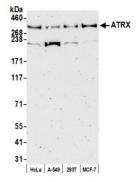

Detection of human ATRX by western blot. Samples: Whole cell lysate (50 μ g) from HeLa, A-549, HEK293T, and MCF-7 cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-ATRX antibody A301-045A-T (lot A301-045A-T-4) used for WB at 1:1000. Detection: Chemiluminescence with an exposure time of 3 minutes.

Bethyl Laboratories, Inc. • 25043 West FM 1097 • Montgomery, TX 77356 • 800.338.9579 • 936.597.6111 • 866.597.6105 (FAX) • www.bethyl.com • technical@bethyl.com