



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien


T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

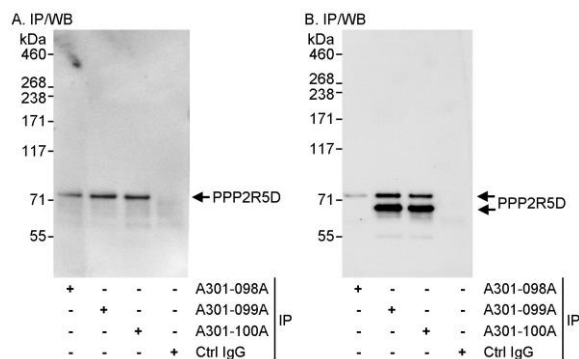

Detection of human PPP2R5D by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1 mg for IP, 20% of IP loaded) from HeLa cells. *Antibody:* Affinity purified rabbit anti-PPP2R5D antibody A301-099A-T used for IP at 3 µg/mg lysate. PPP2R5D was also immunoprecipitated by rabbit anti-PPP2R5D antibodies A301-098A and A301-100A, which recognize other epitopes. For blotting immunoprecipitated PPP2R5D, A301-098A (A) and A301-100A (B) were used at 0.1 µg/ml. *Detection:* Chemiluminescence with exposure times of 10 seconds (A) and 30 seconds (B).