



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien



T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

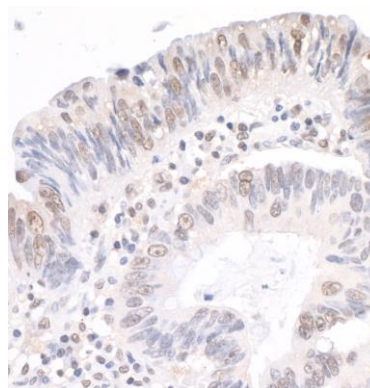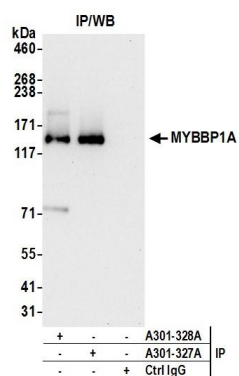

Detection of human MYBBP1A by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1.0 mg per IP reaction; 20% of IP loaded) from HeLa cells prepared using NETN lysis buffer. *Antibodies:* Affinity purified rabbit anti-MYBBP1A antibody A301-327A-T (lot A301-327A-T-2) used for IP at 3 μ g per reaction. MYBBP1A was also immunoprecipitated by rabbit anti-MYBBP1A antibody A301-328. For blotting immunoprecipitated MYBBP1A, A301-328A was used at 0.4 μ g/ml. *Detection:* Chemiluminescence with an exposure time of 10 seconds.

Detection of human MYBBP1A by immunohistochemistry. *Sample:* FFPE section of human colon carcinoma. *Antibody:* Affinity purified rabbit anti-MYBBP1A (Cat. No. A301-327A Lot 2) used at a dilution of 1:1,000 (1 μ g/ml). *Detection:* DAB