

## Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

## SZABO-SCANDIC HandelsgmbH


Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com

## KANK2/SIP Antibody

| ,  |       |   |            |                            | ( ) m                              |
|--|-------|---|------------|----------------------------|------------------------------------|
| Rabbit Polyclonal<br>Antigen Affinity Purified |       |   | Protein ID | NP_056308.2                |                                    |
| Catalog No.                                    | A301- | 408A  | GenelD     | 25959                      | DETLIVI                            |
| Lot No.  | A301- | 408A-1  |            |                            | LABORATORIES, INC                  |
| APPLICATIONS                                   |       | IP  |            |                            |                                    |
| SPECIES REACTIVITY                             |       | Human   |            |                            |                                    |
| PRESUMED REACTIVITY                            |       | Based on 100% sequence identity, this antibody is predicted to react with Mouse and Bovine  |            |                            |                                    |
| AMOUNT   |       | 100 µl  |            |                            |                                    |
| CONCENTRATION                                  |       | 1000 μg/ml  |            |                            |                                    |
| STORAGE/SHELF LIFE                             |       | 2 – 8° C / 1 year from date of receipt  |            |                            |                                    |
| PHYSICAL STATE                                 |       | Liquid  |            |                            |                                    |
| BUFFER   |       | Tris-citrate/phosphate buffer, pH 7 to 8 containing 0.09% Sodium Azide  |            |                            |                                    |
| ISOTYPE  |       | IgG   |            |                            |                                    |
| ORIGIN   |       | USA   |            |                            |                                    |
| PRODUCTION<br>PROCEDURES                       |       | Antibody was affinity purified using an epitope specific to KANK2/SIP immobilized on solid support.   |            |                            |                                    |
|  |       | The epitope recognized by A301-408A maps to a region between residue 625 and 675 of human KN motif and ankyrin repeat domains 2 (SRC-1 interacting protein) using the numbering given in entry NP_056308.2 (GeneID 25959).                            |            |                            |                                    |
|  |       | Immunoglobulin concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.  |            |                            |                                    |
| APPLICATIONS                                   |       | Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.   |            |                            |                                    |
|  |       | Western Blot  | No         | ot recommended. Use rabbit | anti-KANK2/SIP antibody A301-407A. |
|  |       | Immunoprecipi   | tation 2 - | – 5 µg/mg lysate           |                                    |
| APPLICATION NOTES                              |       | Western blot of immunoprecipitates performed using Normal Pig Serum (Cat. No. S100–020),<br>Goat anti–Rabbit Light Chain HRP Conjugate (Cat. No. A120–113P) and 4–8% SDS–PAGE<br>(link to IP–western blot protocol in Additional Info section below). |            |                            |                                    |
| ADDITIONAL INFO                                |       | https://www.bethyl.com/product/A301-408A  |            |                            |                                    |
|  |       | Use the link above to view SDS, a current list of citations, and other product specific information.<br>IP-western blot protocol: https://www.bethyl.com/content/protocol_IP_WB   |            |                            |                                    |

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer Date: June 21, 2019

For in vitro laboratory use only. Not for any clinical, therapeutic or diagnostic use in humans or animals. Not for human or animal consumption. This product may not be resold or modified for resale without the prior written approval of Bethyl Laboratories, Inc. The information provided in this data sheet is believed to be correct but does not purport to be all-inclusive and is intended to be used as a guide. Bethyl Laboratories, Inc. shall not be liable or responsible in any way for use of either this information or the material supplied. Disposal of hazardous material may be subject to federal, state or local laws or regulations.

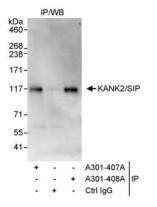

Detection of human KANK2/SIP by western blot of immunoprecipitates. *Samples:* Whole cell lysate (1 mg for IP, 20% of IP loaded) from HeLa cells. *Antibodies:* Affinity purified rabbit anti-KANK2/SIP antibody A301-408A used for IP at 3 µg/mg lysate. KANK2/SIP was also immunoprecipitated by rabbit anti-KANK2/SIP antibody A301-407A, which recognizes an upstream epitope. For blotting immunoprecipitated KANK2/SIP, A301-407A was used at 1 µg/ml. *Detection:* Chemiluminescence with an exposure time of 30 seconds.