

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand


SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


TAB182 Antibody							
Rabbit Polyclonal Antigen Affinity Purified			Protein ID	NP_203754.2			
Catalog No. A301-		439A	GenelD	85456		DETLIVI	
Lot No.	A301-4	439A-2	19A-2				
APPLICATIONS		WB, IP					
SPECIES REACTIVITY		Human					
AMOUNT		100 µl					
CONCENTRATION		200 µg/ml					
STORAGE/SHELF LIFE		2 - 8° C / 1 year from date of receipt					
PHYSICAL STATE		Liquid					
BUFFER		Tris-buffered Saline containing 0.1% BSA and 0.09% Sodium Azide					
ISOTYPE		IgG					
ORIGIN		USA					
PRODUCTION PROCEDURES		Antibody was affinity purified using an epitope specific to TAB182 immobilized on solid support.					
		The epitope recognized by A301-439A maps to a region between residue 1679 and 1729 of human tankyrase 1-binding protein of 182 kDa using the numbering given in entry NP_203754.2 (GeneID 85456).					
		Immunoglobulin concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.					
APPLICATIONS		Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.					
		Western Blot	1:2	2,000 - 1:10,000			
		Immunoprecipi	itation 2 –	10 µg/mg lysate			
APPLICATION NOTES		Western blot of immunoprecipitates performed using Normal Pig Serum (Cat. No. S100–020), Goat anti-Rabbit Light Chain HRP Conjugate (Cat. No. A120–113P) and 4–8% SDS–PAGE (link to IP-western blot protocol in Additional Info section below).					
Western blot of lysates performed using standard western blot reagents ar					and 4-8% SDS-PAGE.		
ADDITIONAL INFO		https://www.bethyl.com/product/A301–439A					
		Use the link above to view SDS, a current list of citations, and other product specific information. IP-western blot protocol: https://www.bethyl.com/content/protocol_IP_WB					

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc. Eric McIntush, PhD | Chief Scientific Officer Date: June 21, 2019

For in vitro laboratory use only. Not for any clinical, therapeutic or diagnostic use in humans or animals. Not for human or animal consumption. This product may not be resold or modified for resale without the prior written approval of Bethyl Laboratories, Inc. The information provided in this data sheet is believed to be correct but does not purport to be all-inclusive and is intended to be used as a guide. Bethyl Laboratories, Inc. shall not be liable or responsible in any way for use of either this information or the material supplied. Disposal of hazardous material may be subject to federal, state or local laws or regulations.

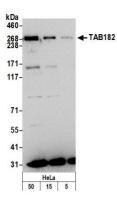

Detection of human TAB182 by western blot. Samples: Whole cell lysate (5, 15 and 50 μ g) from HeLa cells prepared using NETN lysis buffer. Antibody: Affinity purified rabbit anti-TAB182 antibody A301-439A (lot A301-439A-2) used for WB at 0.1 μ g/ml. Detection: Chemiluminescence with an exposure time of 3 minutes.

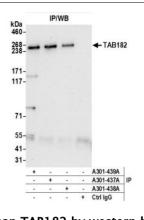

Detection of human TAB182 by western blot of immunoprecipitates. *Samples:* Whole cell lysate (0.5 or 1.0 mg per IP reaction; 20% of IP loaded) from HeLa cells prepared using NETN lysis buffer. *Antibodies:* Affinity purified rabbit anti-TAB182 antibody A301-439A (lot A301-439A-2) used for IP at 6 µg per reaction. TAB182 was also immunoprecipitated by rabbit anti-TAB182 antibodies A301-437A and A301-438A. For blotting immunoprecipitated TAB182, A301-439A was used at 1 µg/ml. *Detection:* Chemiluminescence with an exposure time of 3 minutes.