



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Goat IgG Heavy and Light Chain Cross-Adsorbed Antibody

Rabbit Polyclonal Conjugate DyLight® 800

Antigen Affinity Purified

Catalog No. A50-200D8

Lot No. A50-200D8-9

APPLICATIONS WB, IHC, ICC, Flow Cyt, IF

SPECIES REACTIVITY Goat. Minimum reactivity to chicken, horse, human, mouse, pig and rat

AMOUNT 1 ml

CONCENTRATION 0.5 mg/ml

STORAGE/SHELF LIFE 2 – 8°C / 1 year from date of receipt

PHYSICAL STATE Liquid

BUFFER Phosphate Buffered Saline (PBS) containing 0.09% Sodium Azide

FLUOROPHORE/PROTEIN 4.4

ISOTYPE IgG

ORIGIN USA

PRODUCTION PROCEDURES Antiserum was cross adsorbed using chicken, horse, human, mouse, pig and rat immunosorbents to remove cross reactive antibodies. The antibody to goat IgG was isolated by affinity chromatography using antigen coupled to agarose beads and conjugated to DyLight® 800.

Antibody concentration was determined by extinction coefficient: absorbance at 280 nm of 1.4 equals 1.0 mg of IgG.

By immunoelectrophoresis and ELISA this antibody reacts specifically with goat IgG and with light chains common to other goat immunoglobulins. No antibody was detected against non-immunoglobulin serum proteins. Less than 1% cross reactivity to chicken, horse, human, mouse, pig and rat IgG was detected.

This antibody may cross react with IgG from other species.

APPLICATIONS Centrifuge tube to remove product from lid. Optimal working dilutions should be determined experimentally by the investigator. Prepare working dilution immediately before use.

Western Blot 1:1,000 – 1:20,000. 5% non-fat dry milk in PBST or TBST is recommended for blocking and incubation of antibodies. BSA is not recommended.

Immunohistochemistry 1:50 – 1:500

Immunocytochemistry 1:50 – 1:500

Flow Cytometry 1:50 – 1:200

Immunofluorescence 1:50 – 1:500

APPLICATION NOTES Not all listed applications have been specifically tested by our laboratory.

DyLight® 800 is excited at 770 (in PBS) and emits at 794 (in PBS).


DyLight® is a trademark of Thermo Fisher Scientific Inc. and its subsidiaries.

ADDITIONAL INFO <https://www.fortislife.com/p/A50-200D8>

Use the link above to view SDS, a current list of citations, and other product specific information.

This document certifies that this product has met all of the quality control standards defined by Bethyl Laboratories, Inc.
Michael Spencer, PhD

Date: January 19, 2023

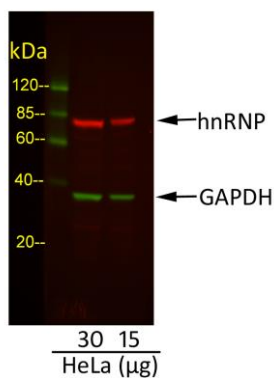

Detection of GAPDH and hnRNP in HeLa Whole Cell Lysate. *Primary Antibodies:* cocktail of goat anti-GAPDH A303-878A (A303-878A-1) and mouse anti-hnRNP A500-011A (A500-011A-1) at 1 µg/ml each. *Secondary Antibodies:* cocktail of Dylight® 800-conjugated rabbit anti-goat A50-200D8 (A50-200D8-1) (blue) and Dylight® 680-conjugated rabbit anti-mouse A90-317D6 (A90-317D6-1) (red) at 0.5 µg/ml each. *Acquisition:* Syngene G:Box, 47 seconds (green) and 25 seconds (red).