



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

TECHNICAL DATA SHEET



bioauxilium

BETTER TOOLS. REAL DISCOVERIES.

THUNDER™ Total MEK1 TR-FRET Cell Signaling Assay Kit

CATALOG NUMBERS KIT-MEKIT-100 (100 tests)
KIT-MEKIT-500 (500 tests)

Store at -80°C
For research use only.
Not for use in diagnostic procedures.

PRODUCT DESCRIPTION

This assay kit measures intracellular levels of **total MEK1** protein in cell lysates using a simple, rapid and sensitive immunoassay based on the homogeneous (no-wash) THUNDER™ TR-FRET technology. The kit is compatible with both adherent and suspension cells.

SPECIFICITY

This assay kit contains two specific and selective antibodies that recognize **total** (both phosphorylated and unphosphorylated) **MEK1**.

SPECIES REACTIVITY

Human (Swiss-Prot Acc.: Q02750; Entrez-Gene Id: 5604).

Other species should be tested on a case-by-case basis.

TR-FRET ASSAY PRINCIPLE

The **Total MEK1** assay kit is a homogeneous time-resolved Förster resonance energy transfer (TR-FRET) sandwich immunoassay (Figure 1). The THUNDER™ Cell Signaling assay workflow consists of 3 steps (Figure 2). Following cell treatment, cells are first lysed with the specific Lysis Buffer provided in the kit. Then **Total MEK1** in the cell lysates is detected with a pair of fluorophore-labeled antibodies in a simple "add-incubate-measure" format (single-step reagent addition; no wash steps). One antibody is labeled with a donor fluorophore (Europium chelate; Eu-Ab1) and the second with a far-red acceptor fluorophore (FR-Ab2). The binding of the two labeled antibodies to distinct epitopes on the target protein takes place in solution and brings the two dyes into close proximity. Excitation of the donor Europium chelate molecules with a flash lamp (320 or 340 nm) or a laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which in turn emit a TR-FRET signal at 665 nm. Residual energy from the Eu chelate generates light at 615 nm. The signal at 665 nm is proportional to the concentration of **Total MEK1** in the cell lysate. Data can be expressed as either the signal at 665 nm or the 665 nm/615 nm ratio.

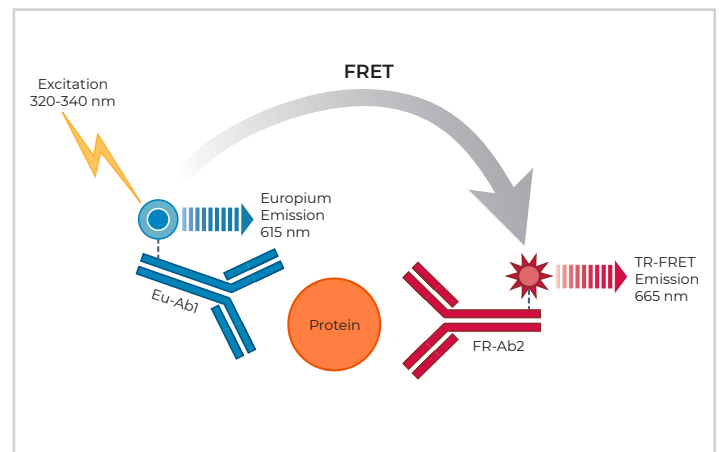


Figure 1 Schematic representation of the TR-FRET cell signaling assay principle.

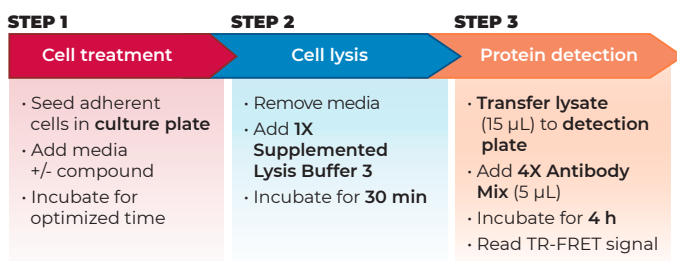


Figure 2 Assay workflow using the 2-plate (transfer) protocol.

KIT COMPONENTS

	100 points*	500 points*
Eu-labeled total-MEK1 antibody (Eu-Ab1)	5 µL	25 µL
Acceptor-labeled total-MEK1 antibody (FR-Ab2)	20 µL	100 µL
Lysis Buffer 3 (5X)	1 mL	5 mL
Detection Buffer (10X)	50 µL	250 µL
Positive control cell lysate	100 µL	500 µL


* The number of assay points is based on an assay volume of 20 µL in half-area 96-well or low-volume 384-well assay plates using the kit components at the recommended concentrations (refer to the User Manual).

VALIDATION DATA

This assay kit has been validated for the relative quantification of total MEK1 in HeLa cell lysates using the 2-plate assay protocol.


- Adherent cells were cultured overnight in a 96-well tissue culture plate (DMEM +10% FBS).
- Following cell treatment, the media was removed and cells were lysed with the 1X **Lysis Buffer 3** (50 μ L) supplemented with the phosphatase inhibitors sodium fluoride (1 mM) and sodium orthovanadate (2 mM).
- Following a **30-min** incubation at room temperature (RT) on an orbital shaker (400 rpm), lysates (15 μ L) were then transferred to a 384-well assay plate followed by addition of the labeled antibodies Eu-Ab1 and FR-Ab2 (5 μ L) for detection of total MEK1.
- The plate was incubated at RT for **4 hours** and the TR-FRET signal was recorded at 665 and 615 nm (EnVision[®]; lamp excitation).

STIMULATION OF PHOSPHO-MEK1 (S218/S222) IN HELA CELLS

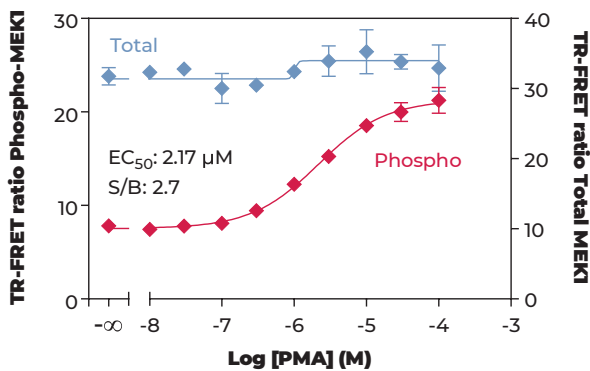


HeLa cells (25,000 cells/well; in triplicate) were incubated with serial dilutions of PMA for 5 min at 37°C. Data show that treatment of HeLa cells with PMA stimulates phosphorylation of MEK1 at S218/S222, but does not affect the levels of total MEK1.

INHIBITION OF PHOSPHO-MEK1 (S218/S222) IN HELA CELLS

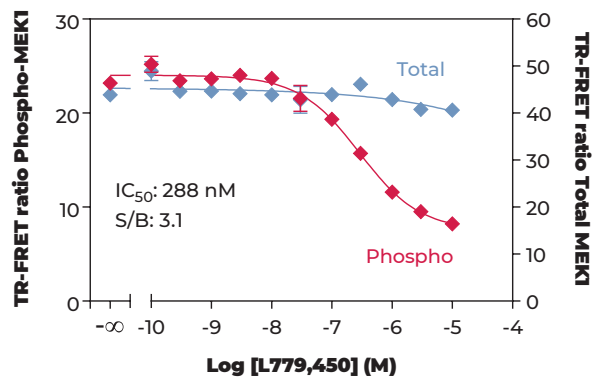

HeLa cells (25,000 cells/well; in triplicate) were incubated with serial dilutions of the inhibitor L779,450 for 15 min at 37°C. Cells were then stimulated with 30 μ M PMA for 5 min at 37°C. Data show that treatment of HeLa cells with L779,450 inhibits phosphorylation of MEK1 at S218/S222 by PMA, but does not affect the levels of total MEK1.

HELA CONTROL LYSATE TITRATION (QC TEST)

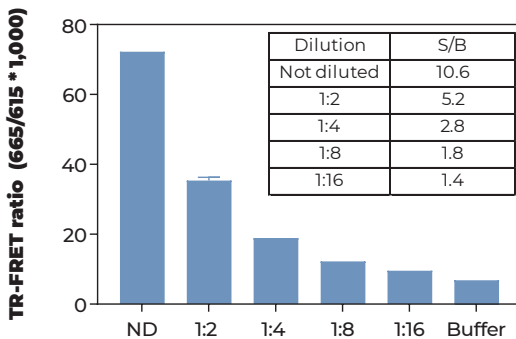

Quality Control: the Total MEK1 assay kit is routinely tested against PMA-treated HeLa lysates. HeLa cells were cultured in a T175 flask to 90% confluence and stimulated with 30 μ M PMA for 5 min at 37°C. Following cell lysis using 4 mL of 1X Lysis Buffer 3, lysates were serially diluted with 1X Lysis Buffer 3 and tested in triplicate. Data show a linear relationship between lysate dilutions and TR-FRET ratio values.

