



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien


T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

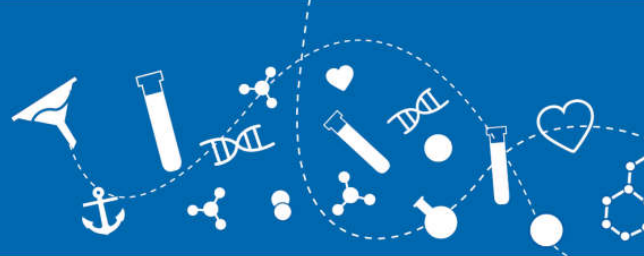

## Galectin-8, His tag, mouse recombinant (rmLGALS8-His)

**Catalog No:** 97604  
**Lot No:** XXXXX  
**Source:** *E. coli*  
**Synonyms:** Galectin-8, Gal-8, LGALS-8, A1326142, D13Ert524e, 1200015E08Rik

### Background

LGALS8 is a prostate-specific antigen that is solely overexpressed in malignant tumors and thus is a supplementary specific identifier of malignancies. LGALS8 is part of the galectin gene family which facilitates both cell-cell and cell matrix interactions in a method parallel to the selectin subgroup of C-type lectins.

### Description

LGALS8 mouse recombinant produced *E. coli* is a single polypeptide chain containing 339 amino acids (1-316) and having a molecular mass of 38 kDa. LGALS8 is fused to a 23 amino acid His Tag at N-terminus and purified by proprietary chromatographic techniques.

### Physical Appearance

Sterile filtered colorless solution.

### Formulation

The LGALS8 solution (0.5 mg/ml) contains 20 mM Tris-HCl buffer (pH 8.0), 0.15 M NaCl, 20% glycerol and 1 mM DTT.

### Stability

Store at 4°C if entire vial will be used within 2-4 weeks. Store frozen at -20°C for longer periods of time. For long term storage it is recommended to add a carrier protein (0.1% HSA or BSA). Avoid multiple freeze-thaw cycles.

### Purity

Greater than 90.0% as determined by SDS-PAGE.

### Amino Acid Sequence

MGSSHHHHHH SSSLVPRGSH MGSMLSLNLL QNIIYNPIIP YVGTITEQLK PGSLIVIRGH VPKDSERFQV DFQLGNSLKP  
RADVAHFHNP RFKRSSCIVC NTLTQEKWGW EEITYDMPFR KEKSFEIVFM VLKNKFQVAV NGRHVLLYAH RISPEQIDTV  
GIYGKVNIIHS IGFRFSSDLQ SMETSALGLT QINRENIQKP GKLQLSLPFE ARLNASMGPG RTVVIKGEVN TNARSFNVDL  
VAGKTRDIAL HLNPRNLNVKA FVRNSFLQDA WGEERNITC FPFSSGMYFE MIIYCDVREF KVAINGVHSL EYKHRFKDLS  
SIDTLSVDGD IRLLDVRSW

### Activity

The ED50 for this effect is 3.9 - 15.63 µg/ml. Measured by its ability to agglutinate human red blood cells.

### Usage

**This product is offered by Biomol for research purposes only. Not for diagnostic purposes or human use. It may not be resold or used to manufacture commercial products without written approval of Biomol GmbH.**