



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com


[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Product Information

Description:	Recombinant human PKN3 also known as PRK3 (protein kinase N3), full length. This construct contains an N-terminal GST-tag. The recombinant protein was affinity purified and kinase active.
Species:	Human
Construct:	PKN3 (GST-Full Length)
Concentration:	0.05 mg/ml
Expression System:	Sf9
Purity:	70%
Format:	Aqueous buffer solution.
Formulated In:	50 mM Tris-HCl, pH 7.5, 150 mM NaCl, 10 mM glutathione, 0.1 mM EDTA, 0.25 mM DTT, 0.1 mM PMSF, 25% glycerol
MW:	130 kDa
Genbank Accession:	NM_013355
Stability:	At least 6 months at -80°C.
Storage:	-80°C
Instructions for Use:	Thaw on ice and gently mix prior to use. DO NOT VORTEX. Perform a quick spin before opening. Aliquot into small volumes and flash freeze for long term storage. Avoid multiple freeze/thaw cycles.
Specific Activity:	2.4 pmol/min/μg
Assay Conditions:	<p>Kinase activity was measured using the ADP-Glo™ Kinase Assay Kit (Promega; Cat# V9101) which quantifies the amount of ADP produced. The ADP-Glo™ Reagent is added to terminate the reaction and deplete the remaining ATP. The Kinase Detection Reagent is then added to convert ADP to ATP and to measure the newly synthesized ATP using a luciferase reaction.</p> <p>PKN3/PRK3 activity was measured by using RSK substrate (KRRRLSSLRA) diluted in water to a final concentration of 1 mg/ml. Reaction was initiated by mixing increasing amounts of the PKN3/PRK3 with 25 μM ATP in 40 mM Tris-HCl, pH 7.4, 20 mM MgCl₂, 0.1 mg/ml BSA prepared with 250 μM DTT and the 20 μg/ml substrate.</p> <p>After a 40-minute incubation at 37°C, the reaction was terminated by addition of the AMP-Glo™ Reagent followed by a subsequent 40-minute incubation at room temperature. Kinase Detection Reagent was then added and incubated for another 30 minutes. Detection of luminescence was measured using the Luminescence Module Protocol on GloMax®-Multi Microplate reader. The corrected activity (RLU) was calculated by removing the blank value for each sample divided by the (specific activity of ADP in RLU/pmol)*(Reaction time in min)*(Enzyme amount in μg or mg). The blank was determined from a “no kinase” sample by replacing the enzyme working solution with an equal volume of Kinase Dilution Buffer IX (1X).</p>
Applications:	Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.

Quality Control Data

4-20% SDS-PAGE Coomassie Staining

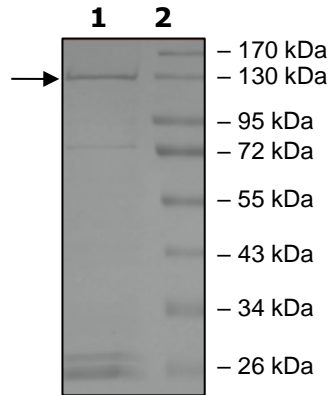

Specific Activity

