



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Ct. West, Ste. B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet

Protein Tyrosine Kinase Substrate (poly-Glu,Tyr 4:1)

Ready to use PTK substrate (50X)

Catalog #: 40217 Lot #: 191001

Size: 100 μ l @ 10 mg/ml (50X)

Formulation: Proprietary buffer for peptide reconstitution

Stability: At least 6 months at -20°C .

References:

1. Varkondi E. *et al.*, *J of Receptors and Signal Transduction*, 2005; **25(1)**; 45-56.

Description:


Poly-(Glu,Tyr 4:1)peptide is reconstituted in BPS proprietary buffer. 0.2 mg/ml at final concentration is enough for a majority of PTKs.

Application:

Useful as a substrate for protein tyrosine kinase inhibitor screening.

Quality Assurance

PTK Substrate, poly (Glu, Tyr 4:1)

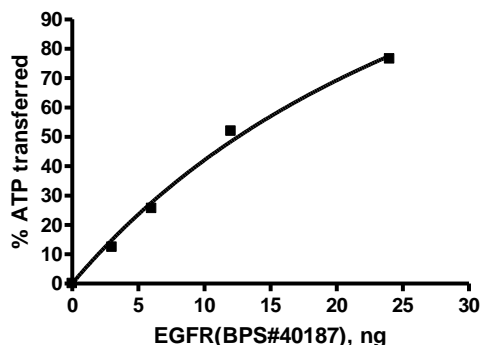

Enzyme reaction was performed at 30°C for 40 min utilizing 20 μM ATP and PTK substrate (0.2 mg/ml at final). Luminescence signal was generated using Kinase-Glo max reagent (Promega).

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com