



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: info@bpsbioscience.com

**Data Sheet**  
**SMARCA2 TR-FRET Assay Kit**  
Catalog # 40342  
Size: 384 reactions

**DESCRIPTION:** The SMARCA2 TR-FRET Assay Kit is designed to measure the inhibition of SMARCA2 (BRM) binding to its substrate in a homogeneous 384 reaction format. This FRET-based assay requires no time-consuming washing steps, making it especially suitable for high throughput screening applications. The assay procedure is straightforward and simple; a sample containing terbium-labeled donor, dye-labeled acceptor, SMARCA2, substrate, and an inhibitor is incubated for 120 minutes. Then, the fluorescence intensity is measured using a fluorescence reader.

**COMPONENTS:**

Catalog #	Component	Amount	Storage	
31141	SMARCA2	10 µg	-80°C	<b>(Avoid freeze/thaw cycles!)</b>
	Bromodomain Ligand 2	50 µl	-80°C	
	Tb-labeled donor	2 x 10 µl	-20°C	
	Dye-labeled acceptor	2 x 10 µl	-20°C	
33012	3x BRD TR-FRET Assay Buffer 1	4 ml	-20°C	
	White, Nonbinding, low volume, microtiter plate	1	Room temp.	

**MATERIALS REQUIRED BUT NOT SUPPLIED:**

Fluorescent microplate reader capable of measuring Time Resolved Fluorescence Resonance Energy Transfer (TR-FRET)  
Adjustable micropipettor and sterile tips

**APPLICATIONS:** Great for screening small molecular inhibitors for drug discovery and HTS applications.

**STABILITY:** At least 6 months from date of receipt when stored as directed.

**REFERENCE:** Schiaffino-Ortega, S., *et al.*, 2014. *J. Hematol. Oncol.* 7:81.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**  
Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: info@bpsbioscience.com

## ASSAY PROTOCOL:

**All samples and controls should be tested in duplicate.**

- 1) Dilute one part **3x BRD TR-FRET Assay Buffer 1** with 2 parts distilled water (3-fold dilution) to make **1x BRD TR-FRET Assay Buffer 1**. Make only a sufficient quantity needed for the assay; store remaining stock solution in aliquots at  $-20^{\circ}\text{C}$ .
- 2) Dilute **Tb-labeled donor** and **dye-labeled acceptor** 100-fold in **1x BRD TR-FRET Assay Buffer 1**. Make only sufficient quantities needed for the assay; store remaining stock solution in aliquots at  $-20^{\circ}\text{C}$ .
- 3) Add 5  $\mu\text{l}$  of diluted **Tb-labeled donor**, and 5  $\mu\text{l}$  of diluted **dye-labeled acceptor** to each well designated "Test Inhibitor", "Negative Control", and "Positive Control".
- 4) Add 2  $\mu\text{l}$  of inhibitor solution to each well designated "Test Inhibitor". Add 2  $\mu\text{l}$  of the same solution without inhibitor (inhibitor buffer) to the wells labeled "Negative Control", and "Positive Control". *Note: Keep DMSO concentration below 0.5%.*

	Negative Control*	Positive Control	Test Inhibitor
Tb-labeled donor	5 $\mu\text{l}$	5 $\mu\text{l}$	5 $\mu\text{l}$
Dye-labeled acceptor	5 $\mu\text{l}$	5 $\mu\text{l}$	5 $\mu\text{l}$
Test Inhibitor	–	–	2 $\mu\text{l}$
Inhibitor Buffer (no inhibitor)	2 $\mu\text{l}$	2 $\mu\text{l}$	–
Bromodomain Ligand 2	–	5 $\mu\text{l}$	5 $\mu\text{l}$
1x BRD TR-FRET Assay Buffer 1	5 $\mu\text{l}$	–	–
SMARCA2 (0.83 ng/ $\mu\text{l}$ )	3 $\mu\text{l}$	3 $\mu\text{l}$	3 $\mu\text{l}$
<b>Total</b>	<b>20 <math>\mu\text{l}</math></b>	<b>20 <math>\mu\text{l}</math></b>	<b>20 <math>\mu\text{l}</math></b>

- 5) Thaw **Bromodomain Ligand 2** on ice. Upon first thaw, briefly spin tube containing enzyme to recover the full contents of the tube. Aliquot each ligand into single-use aliquots. Store remaining undiluted ligand at  $-80^{\circ}\text{C}$  immediately. *Note: each ligand is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots.*
- 6) Individually dilute **Bromodomain Ligand 2** 40-fold in **1x BRD TR-FRET Assay Buffer 1**. Add 5  $\mu\text{l}$  of diluted **Bromodomain Ligand 2** to each well designated as "Positive Control" and "Test Inhibitor". Add 5  $\mu\text{l}$  of **1x BRD TR-FRET Assay Buffer 1** to the wells labeled as "Negative Control".

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**  
Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
Tel: 1.858.829.3082  
Fax: 1.858.481.8694  
Email: info@bpsbioscience.com

- 7) **Thaw SMARCA2** protein on ice. Upon first thaw, briefly spin tube containing protein to recover the full contents of the tube. Aliquot **SMARCA2** protein into single-use aliquots. Store remaining undiluted **SMARCA2** in aliquots at  $-80^{\circ}\text{C}$  immediately. *Note: SMARCA2 is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted protein.*
- 8) Dilute **SMARCA2** in **1x BRD TR-FRET Assay Buffer 1** to  $0.83\text{ ng}/\mu\text{l}$  ( $2.5\text{ ng}/\text{reaction}$ ). Initiate reaction by adding  $3\ \mu\text{l}$  of diluted **SMARCA2** to wells designated for the "Negative Control" "Positive Control", and "Test Inhibitor". Discard any remaining diluted **SMARCA2** protein after use.
- 9) Incubate at room temperature for 2 hours.
- 10) Read the fluorescent intensity in a microtiter-plate reader capable of TR-FRET.

#### Instrument Settings

Reading Mode	Time Resolved
Excitation Wavelength	$340\pm 20\text{ nm}$
Emission Wavelength	$620\pm 10\text{ nm}$
Lag Time	$60\ \mu\text{s}$
Integration Time	$500\ \mu\text{s}$
Excitation Wavelength	$340\pm 20\text{ nm}$
Emission Wavelength	$665\pm 10\text{ nm}$
Lag Time	$60\ \mu\text{s}$
Integration Time	$500\ \mu\text{s}$

#### CALCULATING RESULTS:

Two sequential measurements should be conducted. Tb-donor emission should be measured at  $620\text{ nm}$  followed by dye-acceptor emission at  $665\text{ nm}$ . Data analysis is performed using the TR-FRET ratio ( $665\text{ nm}$  emission/ $620\text{ nm}$  emission).

When percentage activity is calculated, the FRET value from the negative control can be set as zero percent activity and the FRET value from the positive control can be set as one hundred percent activity.


$$\% \text{ Activity} = \frac{\text{FRET}_s - \text{FRET}_{neg}}{\text{FRET}_p - \text{FRET}_{neg}} \times 100\%$$

Where  $\text{FRET}_s$  = Sample FRET,  $\text{FRET}_{Neg}$  = Negative control FRET, and  $\text{FRET}_p$  = Positive control FRET.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**  
Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**EXAMPLE OF ASSAY RESULTS:**

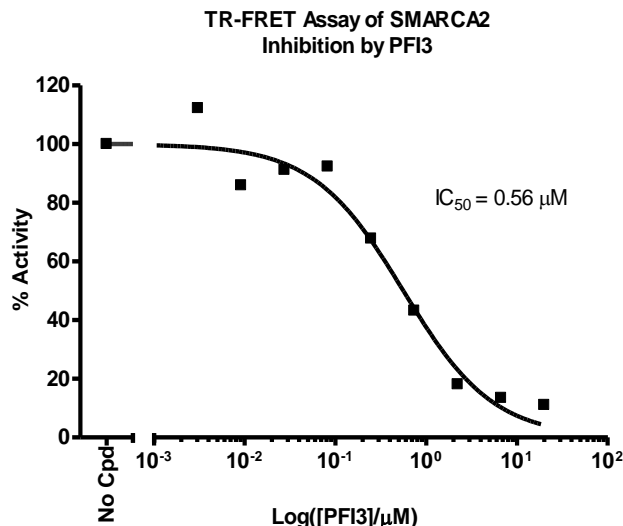

Inhibition of SMARCA2 by PFI3 (BPS Cat. # 27332), measured using the *SMARCA2 TR-FRET Assay Kit*, BPS Bioscience # 40342. Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

**RELATED PRODUCTS:**

<u>Product</u>	<u>Catalog #</u>	<u>Size</u>
SMARCA2 (1375 – 1511), GST-tag	31141	100 µg
SMARCA2 (1375 – 1511), His-tag	31111	100 µg
BRG1 (SMARCA4) (1480 – 1603), GST-tag	31132	100 µg
BRG1 (SMARCA4) (1480 – 1603), His-tag	31102	100 µg
SMARCA2, His-Tag, Biotin-Labeled	31129	50 µg
Bromodomain Ligand 2	33003	0.5 ml
SMARCA2 Inhibitor Screening Assay Kit	32610	384 rxns.
SMARCA4/BRG1 Inhibitor Screening Assay Kit	32609	384 rxns.

Note: Tb-labeled donor and dye-labeled acceptor are products of Cisbio Bioassays.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**  
 Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)