



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

MLK1

Human, recombinant, N-terminal GST-tag

Catalog #: 41102

Lot #: xxxxxx

Conc.: 0.1 mg/ml

Formulated in: 50 mM Tris-HCl, pH 7.5, 150 mM NaCl, 10 mM glutathione, 0.1 mM EDTA, 0.25 mM DTT, 0.1 mM PMSF, 25% glycerol.

Stability: >6 months at -80°C . Avoid freeze/thaw cycles.

References:

1. Dorrow, D.S., *et al.*, *Europ. J. Biochem.* 1993. **213**:701-710.
2. Durkin, J.T., *et al.*, *Biochemistry.* 2004. **43(51)**:16348-16355.


Description: Human MLK1, also known as Mixed Lineage Kinase 1, Mitogen-activated protein kinase kinase kinase 9, MAP3K9, MEKK9, and PRKE1, GenBank Accession No. NM_033141, a.a. 1-433, with N-terminal GST-tag, MW = 77 kDa, expressed in Sf9 cells via a baculovirus expression system.

Specific Activity: 116 pmole/min/ μg

Assay Conditions: Varying concentrations of MLK1 were mixed with 5 μg of MBP protein substrate and 5 μl of 250 μM [^{33}P]-ATP was added for a total volume of 25 μl . Mixture was incubated for 15 min at 30°C . Reaction was stopped by placing 20 μl of reaction onto phosphocellulose P81 paper. Paper was washed 3 times with 1% phosphoric acid solution for 10 min each and air dried. Radioactivity was measured from the P81 paper in the presence of scintillation fluid using a scintillation counter.

Application: Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.

Quality Assurance

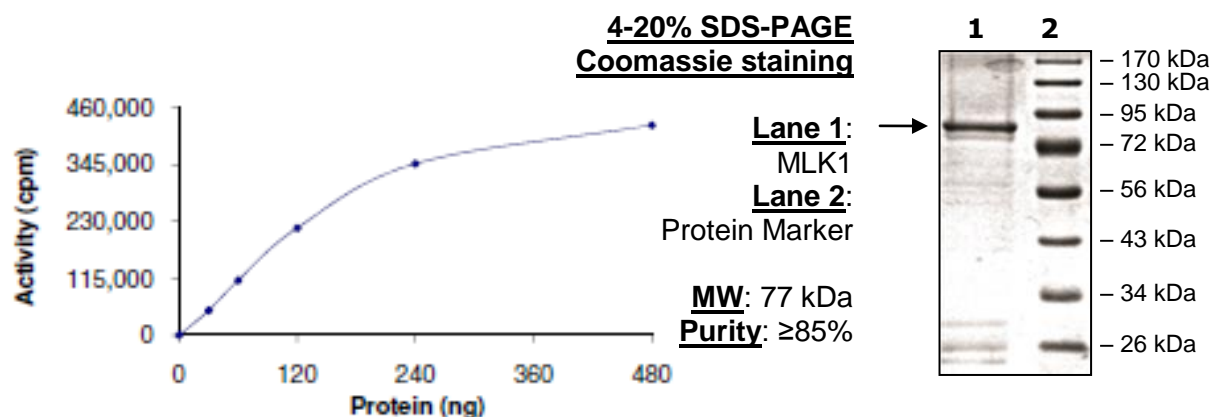

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com