



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6044 Cornerstone Court W, Ste E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet
***Hsp90 β* Assay Kit**
Catalog # 50299
Size: 384 reactions

DESCRIPTION: Hsp90 β is a molecular chaperone with essential functions in maintaining transformation. Inhibition of Hsp90 β function has been shown to play a role in tumorigenesis and disease progression. The *Hsp90 β Assay Kit* is designed for identification of Hsp90 β inhibitors using fluorescence polarization. The assay is based on the competition of fluorescently labeled geldanamycin, an HSP90 inhibitor, for binding to purified recombinant Hsp90 β .

The *Hsp90 β Assay Kit* comes in a convenient 384-well format, with enough purified Hsp90 β enzyme, FITC-labeled geldanamycin, and Hsp90 β assay buffer for 400 enzyme reactions. The key to the *Hsp90 β Assay Kit* is the fluorescently labeled geldanamycin. Using this kit, only one simple step on a microtiter plate is required for Hsp90 β reactions. The FITC-labeled geldanamycin is incubated with a sample containing Hsp90 β enzyme to produce a change in fluorescent polarization that can then be measured using a fluorescence reader.

COMPONENTS:

Catalog #	Component	Amount	Storage	
50292	Hsp90 β recombinant enzyme	140 μ g	-80°C	Avoid freeze/ thaw cycles!
	FITC-labeled geldanamycin (2.5 μ M)	2 x 30 μ l	-80°C	
50311	5x Hsp90 assay buffer 1	2 x 4 ml	-20°C	
	Black, low binding microtiter plate	1	Room temp.	

MATERIALS OR INSTRUMENTS REQUIRED BUT NOT SUPPLIED:

40 mM DTT
2 mg/ml BSA (bovine serum albumin)

APPLICATIONS: Great for studying enzyme kinetics and screening small molecular inhibitors for drug discovery and HTS applications.

STABILITY: Up to 1 year when stored as recommended.

REFERENCES:

1. Kim J, *et al.*, *Biomol. Screening* 2004; **9(5)**: 375-381.
2. Howes R, *et al.*, *Anal. Biochem.* 2006; **350**:202-213.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Court W, Ste E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

ASSAY PROTOCOL:

Immediately prior to assay:

- 1) Thaw **FITC-labeled geldanamycin** on ice. Upon first thaw, briefly spin tube containing FITC-labeled geldanamycin to recover full content of the tube. Aliquot into single use aliquots. Store remaining FITC-labeled geldanamycin in aliquots at -80°C immediately. *Note: FITC-labeled geldanamycin is very sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles.*
- 2) Thaw **Hsp90 β** on ice. Upon first thaw, briefly spin tube containing Hsp90 β to recover full content of the tube. Aliquot Hsp90 β into single use aliquots. Store remaining Hsp90 β in aliquots at -80°C immediately. *Note: Hsp90 β is very sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles.*

Step 1:

All samples and controls should be tested in duplicate.

- 1) Dilute FITC-labeled geldanamycin (2.5 μM stock) 25-fold with 1x Hsp90 assay buffer to make a 100 nM solution. (Make only sufficient quantity needed for the assay; store remaining 2.5 μM stock solution in aliquots at -80°C .)
- 2) Dilute Hsp90 β in 1x Hsp90 assay buffer to 35 ng/ μl (350 ng/reaction)*. Aliquot any remaining enzyme and store undiluted at -70°C . Keep diluted enzyme on ice. Discard any remaining diluted enzyme after use. **Note: optimal enzyme concentration may vary with the specific activity of the enzyme.*
- 3) Prepare the master mixture: N wells x (7.5 μl **5x Hsp90 assay buffer 1** + 2.5 μl **40 mM DTT** + 2.5 μl **2 mg/ml BSA** + 20 μl **H₂O**). Add 32.5 μl of master mixture to all wells.

	Blank	Enzyme Positive Control	Enzyme Negative Control	Test Inhibitor
5x Hsp90 assay buffer 1	7.5 μl	7.5 μl	7.5 μl	7.5 μl
40 mM DTT	2.5 μl	2.5 μl	2.5 μl	2.5 μl
2 mg/ml BSA	2.5 μl	2.5 μl	2.5 μl	2.5 μl
H ₂ O	20 μl	20 μl	20 μl	20 μl
FITC-Labeled geldanamycin (100 nM)	–	2.5 μl	2.5 μl	2.5 μl
Inhibitor	–	–	–	5 μl
Inhibitor Buffer (no inhibitor)	5 μl	5 μl	5 μl	–
1x HSP90 assay buffer	12.5 μl	–	10 μl	–
Hsp90 β (35 ng/ μl)	–	10 μl	–	10 μl
Total	50 μl	50 μl	50 μl	50 μl

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Court W, Ste E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

- 4) Add 2.5 μ l of diluted **FITC-labeled geldanamycin** (100 nM) to each well designated "Enzyme Positive Control", "Enzyme Negative Control", and "Test Inhibitor."
- 5) Add 5 μ l of **Inhibitor** to each well designated "Test Inhibitor." For the, "Blank", "Enzyme Positive Control" and "Enzyme Negative Control", add 5 μ l of the same solution without Inhibitor (**Inhibitor Buffer**).
- 6) Add 10 μ l of **1x HSP90 assay buffer** to the well designated "Enzyme Negative Control". Add 12.5 μ l **1x Hsp90 assay buffer** to the wells designated "Blank".
- 7) Initiate reaction by adding 10 μ l of **diluted Hsp90 β** (35 ng/ μ l), prepared as described above, to each well designated "Enzyme Positive Control" and "Test Inhibitor." Incubate at room temperature for 2 – 3 hours with slow shaking.

Step 2:

Read fluorescent polarization of the sample in a microtiter-plate reader capable of excitation at wavelengths ranging from 475-495 nm and detection of emitted light ranging from 518-538 nm. Blank value is subtracted from all other values.

CALCULATING RESULTS:

Definition of Fluorescence Polarization:

$$P = \frac{I_{\parallel} - I_{\perp}}{I_{\parallel} + I_{\perp}}$$

Where I_{\parallel} = Intensity with polarizers parallel and I_{\perp} = Intensity with polarizers perpendicular. Most instruments display fluorescence polarization in units of mP.

$$mP = \left(\frac{I_{\parallel} - I_{\perp}}{I_{\parallel} + I_{\perp}} \right) \times 1000$$

The equation above assumes that light is transmitted equally well through both parallel and perpendicular oriented polarizers. In practice, this is generally not true and a correction must be made to measure the absolute polarization state of the molecule. This correction factor is called the "G Factor".

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com




6044 Cornerstone Court W, Ste E
 San Diego, CA 92121
 Tel: 1.858.829.3082
 Fax: 1.858.481.8694
 Email: info@bpsbioscience.com

$$mP = \left(\frac{I_{II} - G(I_1)}{I_{II} + G(I_1)} \right) \times 1000 \quad \text{OR} \quad mP = \left(\frac{G(I_{II}) - I_1}{G(I_{II}) + I_1} \right) \times 1000$$

The G-factor is instrument-dependent and may vary slightly depending upon instrument and conditions. Please check the manual of your instrument to obtain the information about the establishment of the G-factor.

EXAMPLE OF ASSAY RESULTS:

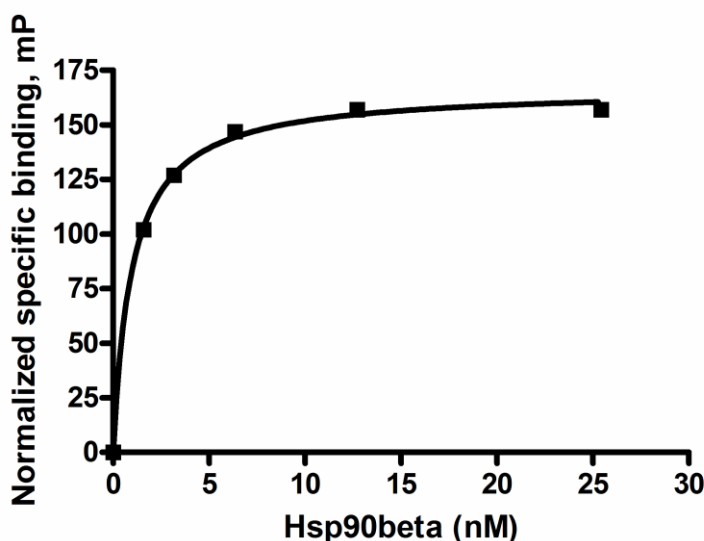

Binding of FITC-geldanamycin to HSP90β, measured using the Hsp90β Assay Kit, BPS Bioscience # 50299. Fluorescence was measured at λ_{ex} 485nm, λ_{em} 530 nm using a Bio-Tek fluorescent microplate reader. *Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at info@bpsbioscience.com*

RELATED PRODUCTS:

<u>Product</u>	<u>Cat. #</u>	<u>Size</u>
Hsp90α recombinant enzyme	50290	200 µg
Hsp90β recombinant enzyme	50292	200 µg
Aha1 recombinant enzyme	50291	200 µg
Geldanamycin inhibitor	27008	5 mg
MS-275 (Entinostat) inhibitor	27011	25 mg
Hsp90α Assay Kit (96 well)	50293	96 rxns
Hsp90α Assay Kit (384 well)	50298	384 rxns
Hsp90β Assay Kit (96 well)	50294	96 rxns

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com