



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6044 Cornerstone Ct. West, Ste. E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet
SMYD3 Homogeneous Assay Kit
Catalog # 52063
Size: 384 reactions

DESCRIPTION: The *SMYD3 Homogeneous Assay Kit* is designed to measure SMYD3 activity for screening and profiling applications. SMYD3 is a histone methyltransferase that promotes di- and trimethylation of lysine4 of histone H3. It also exhibits methylation activity toward a lysine residue of MAPK/ERK Kinase Kinase 2 (MAP3K2). The *SMYD3 Homogeneous Assay Kit* comes in a convenient AlphaLISA[®] format, with GST-tagged SMYD3 substrate, primary antibody, methylation assay buffer, and purified SMYD3 for 384 enzyme reactions. The key to the *SMYD3 Homogeneous Assay Kit* is a highly specific antibody that recognizes methylated substrate. With this kit, only three simple steps on a microtiter plate are required for methyltransferase detection. First, a sample containing SMYD3 enzyme is incubated with the substrate for three hours. Next, acceptor beads and primary antibody are added, then donor beads, followed by reading the Alpha-counts.

COMPONENTS:

Catalog #	Component	Amount	Storage	
51110	SMYD3/HSP90	20 µg	-80°C	Avoid freeze/ thaw cycles!
52120	100 µM S-adenosylmethionine	2 x 250 µl	-80°C	
52140O	Primary antibody 15	5 µl	-80°C	
	SMYD3 substrate	200 µl	-80°C	
	4x SMYD3 assay buffer*	3 x 1 ml	-20°C	
	4x Detection buffer 3	2 ml	-20°C	

*Add 5 µl of 0.5M DTT to each 1 ml tube before use.

MATERIALS OR INSTRUMENTS REQUIRED BUT NOT SUPPLIED:

AlphaLISA[®] anti-rabbit IgG acceptor beads, 5 mg/ml (PerkinElmer #AL104C)
AlphaScreen[®] Glutathione donor beads, 5 mg/ml (PerkinElmer #6765300)
Optiplate-384 (PerkinElmer #6007290)

APPLICATIONS: Great for studying enzyme kinetics and HTS applications.

CONTRAINDICATIONS: Green and blue dyes that absorb light in the AlphaScreen[®] signal emission range (520-620 nm), such as Trypan Blue. Avoid the use of the potent singlet oxygen quenchers such as sodium azide (NaN₃) or metal ions (Fe²⁺, Fe³⁺, Cu²⁺, Zn²⁺ and Ni²⁺). The presence of >1% RPMI 1640 culture medium leads to a signal reduction due to the presence of excess biotin and iron in this medium. MEM, which lacks these components, does not affect AlphaScreen[®] assays.

STABILITY: At least one year from date of receipt when stored as directed.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Ct. West, Ste. E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

- REFERENCES:**
1. Dillon, S.C., *et al.* 2005. *Genome Biology* **6**:227.
 2. Silva, F.P., *et al.* 2008. *Oncogene* **27**:2686.

ASSAY PROTOCOL:

All samples and controls should be tested in duplicate.

Step 1:

- 1) Add 5 μ l of 0.5M DTT to each 1 ml tube before use. Prepare the master mixture: N wells x (2.0 μ l **4x SMYD3 Assay buffer** + 1.0 μ l **S-adenosylmethionine** (100 μ M) + 0.5 μ l **SMYD3 substrate** + 1.5 μ l water). Add 5 μ l to wells designated "Positive Control", "Test Sample", and "Blank". To wells labeled "Substrate Control", add 2 μ l **4x SMYD3 Assay buffer** + 0.5 μ l **SMYD3 substrate** + 2.5 μ l water.
- 2) Add 3 μ l of Inhibitor solution of each well labeled as "Test Inhibitor". For the "Positive Control", "Substrate Control" and "Blank", add 3 μ l of the same solution without inhibitor (Inhibitor buffer).
- 3) Prepare **1x SMYD3 Assay buffer** by adding 1 part of **4x SMYD3 Assay buffer** to 3 parts water (v/v). Prepare only enough for the assay; store any remaining **4x SMYD3 Assay buffer** at -20°C.
- 4) Thaw **SMYD3** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full content of the tube. Aliquot **SMYD3** enzyme into single use aliquots. Store remaining undiluted enzyme in aliquots at -80°C. *Note: SMYD3 is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme. Depending on the nature of inhibitor, pre-incubation with the enzyme may be necessary.*
- 5) To the wells designated as "Blank", add 2 μ l of **1x SMYD3 Assay buffer**.
- 6) Dilute **SMYD3** in **1x SMYD3 Assay buffer** at 13-25 ng/ μ l (26-50 ng/2 μ l). Keep diluted enzyme on ice until use. Discard any unused diluted enzyme after use.
- 7) Initiate reaction by adding 2 μ l of diluted **SMYD3** enzyme to the wells designated "Positive Control", "Substrate Control", and "Test Sample". Incubate at room temperature for 3-4 hours.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Ct. West, Ste. E
 San Diego, CA 92121
 Tel: 1.858.829.3082
 Fax: 1.858.481.8694
 Email: info@bpsbioscience.com

	Positive Control	Test Sample	Substrate Control	Blank
4x SMYD3 Assay buffer	2 µl	2 µl	2 µl	2 µl
100 µM S-adenosylmethionine	1 µl	1 µl	-	1 µl
SMYD3 substrate	0.5 µl	0.5 µl	0.5 µl	0.5 µl
H ₂ O	1.5 µl	1.5 µl	2.5 µl	1.5 µl
Test Inhibitor/Activator	-	3 µl	-	-
Inhibitor Buffer (no inhibitor)	3 µl	-	3 µl	3 µl
1x SMYD3 Assay buffer	-	-	-	2 µl
SMYD3 (13-25 ng/µl)	2 µl	2 µl	2 µl	-
Total	10 µl	10 µl	10 µl	10 µl

Protect your samples from direct exposure to light for steps 2 and 3!

Step 2:

- 1) Prepare **1x Detection buffer 3** by adding 1 part **4x Detection buffer 3** to 3 parts water (v/v). Prepare only enough for the assay; store any remaining **4x Detection buffer 3** at -20°C.
- 2) Dilute anti-Rabbit Acceptor beads 1:250-fold with **1x Detection buffer 3**. Add 5 µl per well. Shake plate briefly.
- 3) Dilute **Primary antibody 15** 400-fold with **1x Detection buffer 3**. Add 5 µl per well. Shake plate. Incubate 30 min at room temperature.
(Alternatively, dilute anti-Rabbit Acceptor beads (1:500) and Primary antibody 4 (1:800) with 1x Detection buffer 3 in one step. Add 10 µL of acceptor beads/antibody mixture per well.)


Step 3:

- 1) Dilute Glutathione donor beads 125-fold with **1x Detection buffer 3**. Add 10 µl per well. Incubate for 30-60 min. at room temperature.
- 2) Read Alpha-counts.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com

Example of Assay Results:

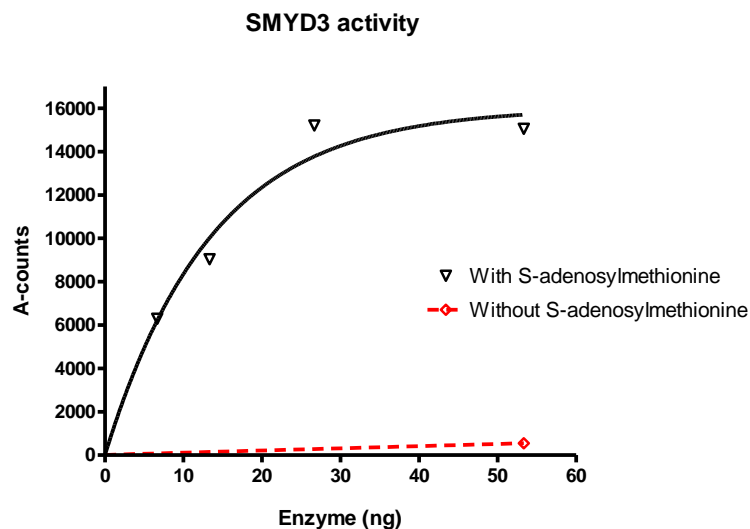

SMYD3 enzyme activity, measured using the *SMYD3 Homogeneous Assay Kit*, BPS Bioscience Cat. #.... Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at info@bpsbioscience.com

RELATED PRODUCTS:

SMYD3 (full length), FLAG-tag	#51031	20 µg
SMYD3 (35-end), GST-tag	#51015	20 µg
SMYD2 (KMT3C), FLAG-tag	#51014	20 µg
SMYD2 (KMT3C) Chemiluminescent Assay Kit	#52055	96 reactions

AlphaScreen® and AlphaLISA® are registered trademarks of PerkinElmer, Inc.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com