



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

PDE3B

Human, recombinant, N-terminal GST tag

Catalog #: 60031

Lot#: 110215G2 **Conc.:** 0.067 mg/ml

Formulated in: 40 mM Tris, pH 8.0, 110 mM NaCl, 2.2 mM KCl, 3 mM DTT, 0.04% Tween-20, 20% glycerol

Stability: >6 months at -80 °C

References:

1. Raymond, D.R. *et al.*, *Cell. Signal.* **19** (12), 2507-2518 (2007).
2. Palmer, D. *et al.*, *J. Biol. Chem.* **282** (13), 9411-9419 (2007).


Description: Human PDE3B (GenBank Accession No. NM_000922), (amino acids 592-end), with N-terminal GST tag, MW = 86 kDa, expressed in a Baculovirus infected Sf9 cell expression system.

Specific Activity: 1270 pmol/min/μg.

One unit is defined as the amount of enzyme that will convert 1 pmole of 3', 5'-cAMP to 5'-AMP per min at 37°C. Assay condition: 10 mM Tris-HCl, pH 7.4, 10 mM MgCl₂, 1 mM MnCl₂, 200 μM cAMP, and 0.2 ng/μl PDE3B at 37°C for 20 min. Quantified by 5'-nucleotidase cleaving the 5'-AMP product and releasing the phosphate group which is detected by Malachite Green Reagent.

Application: Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.

Quality Assurance

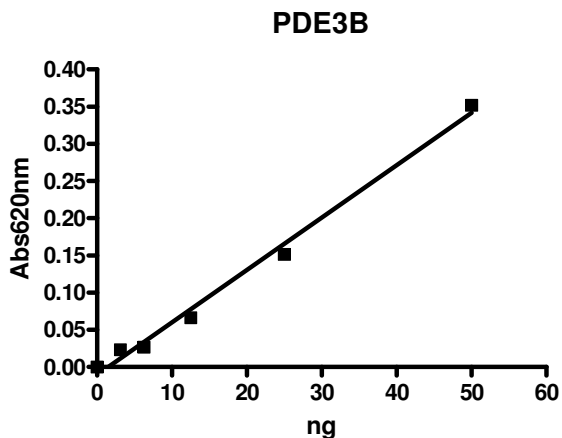


10% SDS-PAGE Coomassie staining

Lane 1:
4 μg PDE3B

Lane 2:
Protein Marker
BioLabs (#P7708L)

MW: 86 kDa
Purity: ≥90%

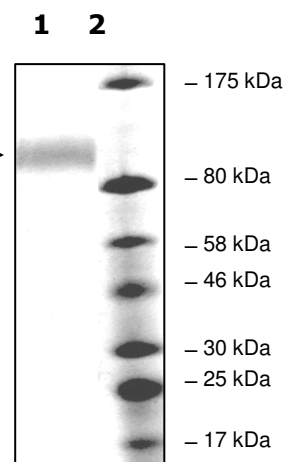

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com