



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6044 Cornerstone Court W, Ste E  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

## Data Sheet

### ***PDE1B-HEK293 Recombinant Cell line***

**Cat. # 60402**

#### **Product description**

Recombinant HEK293 cell line expressing human PDE1B (phosphodiesterase 1B, accession number NM\_000924).

#### **Format**

Each vial contains  $1 \times 10^6$  cells in 1 ml of 10% DMSO.

#### **Introduction**

PDE1B plays a role in signal transduction by regulating the intracellular concentration of cyclic nucleotides. PDE1B hydrolyzes both cAMP and cGMP to nucleoside 5'-monophosphate. Although it prefers cGMP at low substrate levels it has a  $V_{max}$  that is approximately equal for both cGMP and cAMP (Polli, JW and Koncaid, RL, 1994; Yu J. *et al.*, 1997).

cAMP is an important signal carrier necessary for the proper biological response of cells to hormones and other extracellular signals. It is required for cell communication in the hypothalamus/pituitary gland axis and for the feedback control of hormones (Alewijne A.E. *et al.*, 1997).

#### **Culture conditions**

Cells should be grown at 37°C with 7% CO<sub>2</sub> using MEM (Hyclone) medium supplemented with 10% FBS, 1% non-essential amino acid, 1 mM Na-pyruvate, plus 400 µg/ml of Geneticin (Invitrogen) to ensure the recombinant expression is maintained. PDE1B-HEK293 cells should exhibit a typical cell division time of 24 hours.

It is recommended to quickly thaw the frozen cells from liquid nitrogen in a 37°C water bath, transfer to a tube containing 10ml growth medium, spin down cells, resuspend cells and transfer to T25 flask. Cells should be split before they reach complete confluency.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [\*\*info@bpsbioscience.com\*\*](mailto:info@bpsbioscience.com)

Please visit our website at: [\*\*www.bpsbioscience.com\*\*](http://www.bpsbioscience.com)

### Functional validation


N-terminal His-tagged human PDE1B has been stably expressed in HEK293 cell line and its expression was confirmed by Western blotting.

The function of PDE1B was characterized by cAMP detection assay.

In general, forskolin is commonly used to raise levels of cyclic AMP (cAMP). Forskolin resensitizes cell receptors by activating the enzyme adenylyl cyclase and increasing the intracellular levels of cAMP.

When 293 cells are activated by forskolin, cAMP levels are upregulated in parental cells whereas cells overexpressing PDE1B show reduction in cAMP accumulation.

The data provided below represents effect of PDE1B overexpression on cAMP accumulation after forskolin stimulation in 293 cells.

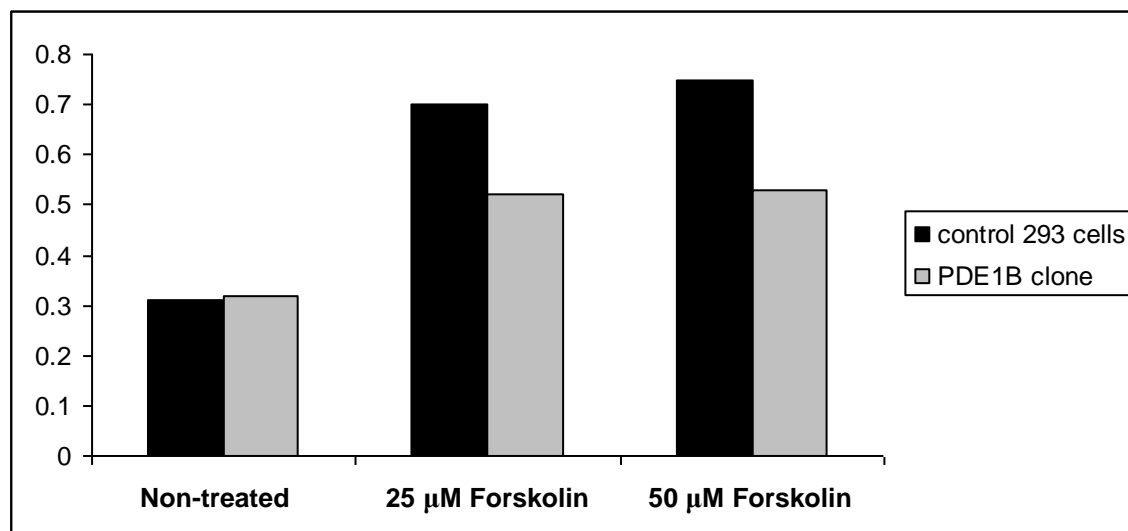

**Figure 1. PDE1B overexpression reduces cAMP accumulation after forskolin stimulation.**

Cells were incubated in PBS buffer for 15 min and consequently stimulated with 0, 25, or 50 μM forskolin for 20 min. cAMP production was measured by fluorescent quenching detection method using a 96-well assay kit (Mediomics).

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6044 Cornerstone Court W, Ste E  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

### Vector and sequence

Human His-PDE1B cloned in pcDNA3.1 (accession number NM\_000924)

His-PDE1B sequence:

MHHHHHHELSPRSPPEMLEESDCPSPLELKSAPSKMMWIKLRSLRLRYMVKQLENGEIN  
IEELKKNLEYTASLLEAVYIDETRQILDTEDELQELRSDAVPSEVRDWLASTFTQQARAK  
GRRAEKPKFRSIVHAVQAGIFVERMFRRTYTSVGPTYSTAVLNCLKNLDLWCFDVFS  
LNQAADDHALRTIVFELLTRHNLISRFKIPTVFLMSFLDALETGYGKYKNPYHNQIHAAD  
VTQTVHCFLLRRTGMVHCLSEIELLAIIFAAAIHDYEHTGTTNSFHIQTKSECAIVYNDRSV  
LENHHISSVFRLMQDDEMNFIFNLTKDEFVELRALVIEMVLATDMSCHFQQVKTMKTAL  
QQLERIDKPKALSLLLHAADISHPTKQWLVHSRWTKALMEEFFRQGDKEAELGLPFSP  
LCDRTSTLVAQSQIGFIDFIVEPTFSVLTDVAEKSVQPLADEDSSKSKNQPSFQWRQPSL  
DVEVGDPNPDVVSFRSTWVKRIQENKQKWKERAASGITNQMSIDELSPCEEEAPPSP  
AEDEHNQNGNLD.

### References

- 1) Polli, JW and Kincaid RL. (1994). "Expression of a calmodulin-dependent phosphodiesterase isoform (PDE1B1) correlates with brain regions having extensive dopaminergic innervation." *J. Neurosc.* **14**: 1251-1261.
- 2) Yu, J, et al. (1997). "Identification and Characterization of a Human Calmodulin-Stimulated Phosphodiesterase PDE1B1" *Cell Signaling* **9** (7): 519-529.
- 3) Bender A.T. *et al.*, (2005). "Selective up-regulation of PDE1B2 upon monocyte-to-macrophage differentiation." *PNAS* **102** (2): 497-502.
- 4) Alewijnse, AE *et al.* (1997). "Modulation of forskolin-mediated adenylyl cyclase activation by constitutively active G<sub>s</sub>-coupled receptors". *FEBS Letters* **419** (2-3): 171-174.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6044 Cornerstone Court W, Ste E  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

## License Disclosure

Purchase of this cell line grants you with a 10-year license to use this cell line in your immediate laboratory, for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make the cell line available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit use of this cell line in humans or for therapeutic or drug use. The license does not permit modification of the cell line in any way. Inappropriate use or distribution of this cell line will result in revocation of the license and result in an immediate cease of sales and distribution of BPS products to your laboratory. BPS does not warrant the suitability of the cell line for any particular use, and does not accept any liability in connection with the handling or use of the cell line. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees; contact [sales@bpsbioscience.com](mailto:sales@bpsbioscience.com) for details. Publications using this cell line should reference BPS Bioscience, Inc., San Diego.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)