



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

B7-1 (CD80), Fc Fusion

Human, recombinant, Fc fusion protein

Catalog #: 71125

Lot #: 150310

Conc.: 2.29 mg/ml

Formulated in: 8 mM Phosphate, pH 7.4, 110 mM NaCl, 2.2 mM KCl, and 20% glycerol.

Stability: At least 6 months at -80°C . Avoid freeze/thaw cycles. Protein may be diluted to $\geq 100 \mu\text{g/ml}$ in PBS + glycerol and stored at -80°C .

References:


1. Byrum, J.N., *et al.*, *J. Immunol.* 2013; **191(6)**: 3073-3081
2. Legat, A., *et al.*, *Front. Immunol.* 2013; **19(4)**: 455.

Description: Human secreted B7-1-Fc fusion protein, also known as CD80. GenBank Accession No. NM_005191, a.a. 35-242 expressed in a HEK293 cell expression system. MW = 50.5 kDa (monomer). This protein runs at a higher MW by SDS-PAGE due to glycosylation.

Assay Description: B7-1 binding to PD-L1 was determined by coating 100 ng of B7-1 per well and adding the indicated amount of PD-L1-biotin (BPS Cat. #71105) and incubating for 2 hr at RT. Strep-HRP was added and allowed to incubate at room 1 hr at RT, followed by HPR chemiluminescent substrate addition and plate reading.

Application: Useful for the study of protein binding studies and screening inhibitors.

Quality Assurance

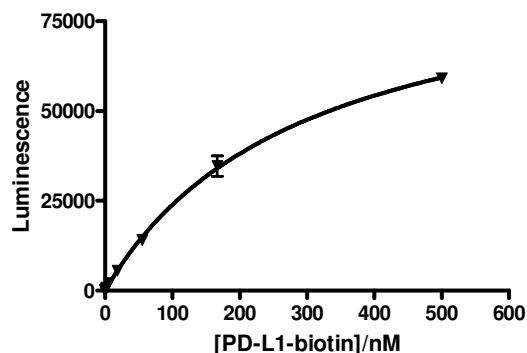


**4-20% SDS-PAGE
 Coomassie staining**

Lane 1:
 3 μg B7-1

Lane 2:
 Protein Marker

MW: 50.5 kDa + glycans

Purity: $\geq 90\%$

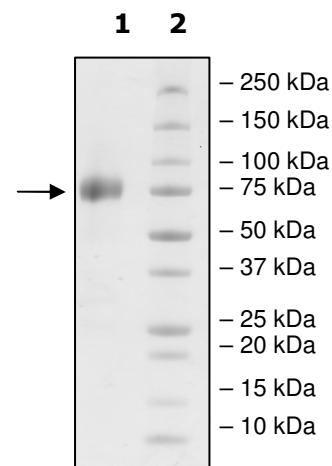

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**


Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Court W, Ste E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Gel Filtration Curve:

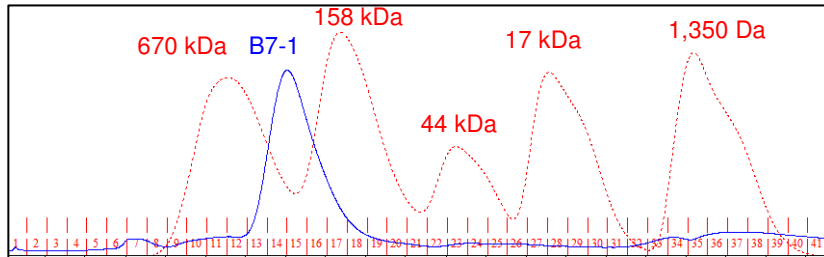

<5% Aggregation

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com