



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

PCSK9, His-tagged

Human, Recombinant, C-terminal His- tag

Catalog #: 71204-1

Lot #: 180206-2 **Conc.:** 1.45 mg/ml

Formulated in: 8 mM Phosphate, pH 7.4, 110 mM NaCl, 2.2 mM KCl, and 20% glycerol.

Stability: At least 6 months at -80°C . Avoid freeze/thaw cycles. Storing diluted protein is not recommended, if necessary, use carrier protein (BSA 0.1 – 0.5%).

References:

1. Abifadel, M., *et al.*, *Nat Genet.* 2003 Jun;**34**(2):154-6.
2. Abifadel, M., *et al.*, *Hum Mutat.* 2009 Apr;**30**(4):520-9.


Description: Human proprotein convertase subtilisin/kexin type 9 (PCSK9), also known as FH3, HCHOLA3, and PC9, GenBank Accession No. NM_174936, a.a. 31-692(end), with C-terminal His-tag, MW=71.8 kDa (calculated), expressed in a HEK293 cell expression system. PCSK9 is autocleaved to the ~14 kDa prodomain (a.a. 31-152) and the ~58 kDa mature form (a.a. 153-692(end)), which run at higher apparent MW by SDS-PAGE.

Assay Conditions: Coated plate with 100 ng/well of PCSK9-his (#71204) in 50 μl (in PBS) (overnight at 4 deg). Initiated binding reaction with addition of various concentrations of LDLR-Biotin in 50 mM MES pH6, 50 mM NaCl, 0.2 mM CaCl_2 , 0.5% Tween-20, and 0.1% BSA (AKA PCSK9-LDLR TR-FRET assay buffer). Incubated for 2h at RT. Detected binding with Strep-HRP generated luminescence.

Application: Useful for studying protein binding and screening small molecules.

Quality Assurance

PCSK9:LDLR[B]

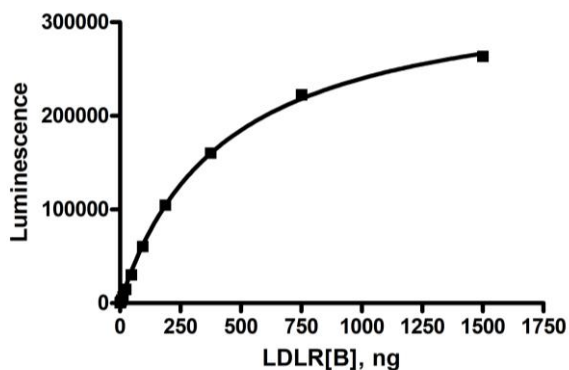


**4-20% SDS-PAGE
Coomassie staining**

Lane 1:
5 μg PCSK9

Lane 2:
Protein Marker

MW:
Pro: 14kDa
Mature: 58kDa

Purity: $\geq 90\%$

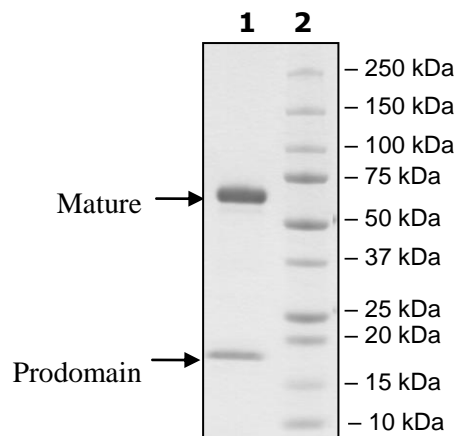

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Court W, Ste E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Protein Sequence:

Green: Pro-peptide: aa 31-152

Black: mature PCSK9 aa 153-692(end)

Red: His tag

QEDEDGDYEELVLALRSEEDGLAEAPEHGTTATFHRCADPWRLPGTYVVVLKEETHL
SQSERTARRLQAQAARRGYLTKILHVFHGLLPGLVKMSGDLLELALKLPHVDYIEEDS
SVFAQSIPWNLERITPPRYRADEYQPPDGGSLVEVYLLDTSIQSDHREIEGRVMVTD
NVPEEDGTRFHRQASKCDSHGTHLAGVVSGRDAGVAKGASMRSLRVLNQCQKGT
SGTLIGLEFIRKSQLVQPVGPLVLLPLAGGYSRVLNAACQLARAGVVLVTAAGNFR
DDACLYSPASAPEVITVGATNAQDQPVTLGTLGTNFGRCVDLFAPGEDIIGASSDCST
CFVSQSGTSQAAAHVAGIAAMMLSAPELTLAELRQRLIHFSKDVINEAWFPEDQR
VLTPNLVAALPPSTHGAGWQLFCRTVWSAHS GPTRMATAVARCAPDEELLSCSSFS
RSGKRRGERMEAQGGKLVCRAHNAFGGEGVYAIARCCLLPQANCSVHTAPPAEAS
MGTRVHCHQQGHVLTGCSHWEVEDLGTHKPPVLRPRGQPNQCVGHREASIHASC
CHAPGLECKVKEHGIPAPQEQT VACEEGWTLTGCSALPGTSHVLGAYAVDNTCVV
RSRDVSTTGSTSEGAVTAVAICCRSRHLAQASQELQH HHHHHH

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com