



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

PD-L1 (K124A), Fc fusion, Avi-Tag, HiP™

Human, Recombinant, Fc fusion protein
with C-terminal Avi-tag™

Catalog #: 71993

Lot #: 160927

Conc.: 1.02 mg/ml

Formulated in: 8 mM Phosphate, pH 7.4,
110 mM NaCl, 2.2 mM KCl, and 20%
glycerol.

Stability: At least 6 months at -80°C . Avoid
freeze/thaw cycles. Protein may be diluted
to $\geq 100 \mu\text{g/ml}$ in PBS + glycerol and stored
at -80°C .

References:

1. Freeman, G.J., *et al.*, *J Exp Med.* 2000
Oct 2;**192(7)**:1027-34.
2. Dong, H., *et al.*, *Nat Med.* 1999
Dec;**5(12)**:1365-9.

Description: Human secreted Programmed
Cell Death 1 Ligand 1 (PD-L1)-Fc fusion
protein, also known as CD274 and B7
homolog 1 (B7-H1) with a K124A mutation,
GenBank Accession No. NM_014143,
amino acids 19-239, fused at the C-
terminus with the Fc portion of human IgG1,
with C-terminal Avi-tag™ expressed in a
HEK293 cell expression system. MW = 54
kDa. This protein runs at a higher MW by
SDS-PAGE due to glycosylation.

Application:

Useful for the study of protein-protein
interaction, epitope mapping, and *in vitro*
and *in vivo* protein function.

Quality Assurance


**4-20% SDS-PAGE
Coomassie staining**

Lane 1:
3 μg PD-L1(K124a)

Lane 2:
Protein Marker

MW: 54k Da+glycans

Purity: $\geq 90\%$

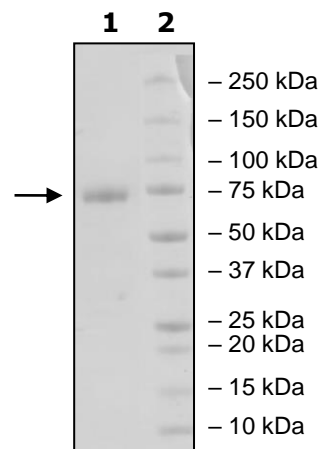

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**


Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Gel Filtration Curve

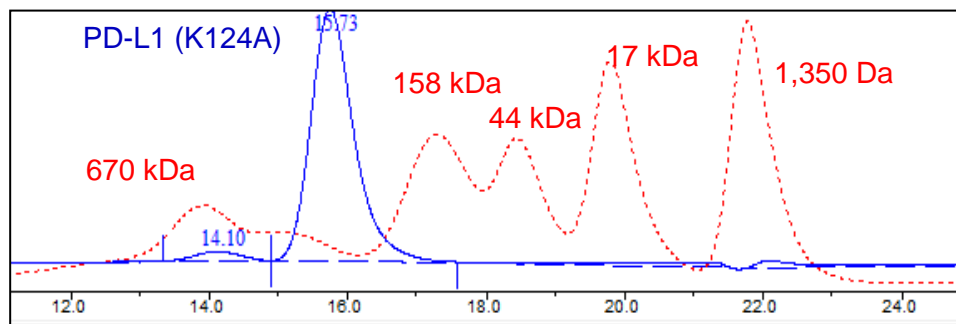

< 10% Aggregation

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com