



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.202.1401
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

Data Sheet
ARG1 Inhibitor Screening Assay Kit
Catalog #72048-1
Size: 96 reactions

BACKGROUND: The *ARG1 Inhibitor Screening Assay Kit* is designed to measure ARG1 enzyme inhibition. The kit comes in a convenient format, with enough reaction solution and enzyme to perform a total of 96 reactions. The *ARG1 Inhibitor Screening Assay Kit* is simple to use. The test inhibitor and enzyme are added to a solution containing thioarginine substrate. After a room temperature incubation, activity is determined by measuring the absorbance of the reaction product at $\lambda = 410$ nm.

DESCRIPTION: Arginase enzymes convert arginine to ornithine through hydrolysis. Two known isoforms of arginase exist, ARG1 and ARG2. These enzymes are involved in the regulation of a variety of immunological responses and are a major target in immunotherapy. ARG1/2 are overexpressed in myeloid-derived suppressor cells (MDSCs) and tumor-associated macrophages (TAMs). Overexpression of ARG1/2 results in depleted levels of arginine both intracellularly and extracellularly. As arginine levels are depleted in the microenvironment, immune cells are starved of this amino acid and the function of key immunological activators become impaired; T cell proliferation is inhibited, regulatory T cells become activated to inhibit CD4+ T cells, and immunosuppressants have increased longevity. Depleted arginine also results in the release of reactive nitrogen species and reactive oxygen species from TAMs and MDSCs. This reactive species cause T cell apoptosis and the activation and growth of antigen presenting cells.

COMPONENTS:

Catalog #	Component	Amount	Storage	
71658	ARG1	25 μ g	-80 °C	<i>Avoid freeze/ thaw cycles!</i>
	10X ARG Assay Buffer	1 ml	-80 °C	
	Thioarginine	1 mg	-80 °C	
	Detection Reagent	3 mg	-80 °C	
	96-well transparent plate	1	Room Temp	

MATERIALS OR INSTRUMENTS REQUIRED BUT NOT SUPPLIED:

Spectrophotometer capable of measuring absorbance at $\lambda = 410$ nm.
Ethanol (200 proof, 100%)

APPLICATIONS: Useful for the study of ARG1 enzymology, screening inhibitors, and selectivity profiling.

CONTRAINDICATIONS: DMSO > 1%, strong acids or bases, ionic detergents, high salt

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**
Or you can Email us at: **info@bpsbioscience.com**



6042 Cornerstone Court W, Ste B
San Diego, CA 92121
Tel: 1.858.202.1401
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

STABILITY: Up to 6 months from date of receipt, when stored as recommended.

REFERENCES:

Sedbrook, J.C., *et al.* *PNAS*. 1999. **96(3)**: 1140-1145.

Woll, P.J., *et al.* *PNAS*. 1998. **85(6)**: 1859-1863.

ASSAY PROTOCOL:

All samples and controls should be tested in duplicate.

ARG Reaction Solution preparation

1. Dilute **10X ARG assay buffer** with water (1:10) to **1X ARG Assay Buffer**.
2. Solubilize 1 mg of **Thioarginine** with 367 μ l **1X ARG Assay Buffer**. Keep on ice until use.
3. Solubilize 3 mg of **Detection Reagent** with 350 μ l ethanol. Ensure complete solubilization by vortexing for approximately 1 minute.
4. Prepare **ARG Reaction Solution** by adding 250 μ l solubilized **Thioarginine** and 56 μ l solubilized **Detection Reagent** to 9.7 ml **1X ARG Assay Buffer**. Keep on ice until use.

Note: Due to instability of the substrate and detection reagent it is recommended to initiate reaction within 1 hour of reaction solution preparation. These materials cannot be refrozen and thawed for later use.

ARG1 Assay:

1. Prepare **ARG1 Reaction Solution** (see above) and aliquot 90 μ l into each well of the assay plate.
2. Add 5 μ l of **inhibitor solution** (including no more than 20% DMSO) to each well designated "Test Inhibitor". For the wells labeled "Positive Control" and "Blank", add 5 μ l of the same solution without inhibitor (inhibitor buffer). Note: Keep the final DMSO concentration \leq 1%.
3. Thaw **ARG1** on ice. Upon first thaw, briefly spin tube containing enzyme to recover full contents of the tube. Aliquot **ARG1** into single use aliquots. Store remaining undiluted enzyme in aliquots at -80°C . Note: **ARG1** is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted enzyme.
4. Dilute **ARG1** in **1X ARG Assay Buffer** to 40 ng/ μ l. Keep diluted protein on ice until use. Discard any unused diluted protein after use.
5. Add 5 μ l of **1X ARG Assay Buffer** to the well designated "Blank".
6. Initiate reaction by adding 5 μ l of diluted **ARG1** (40 ng/ μ l) prepared as described above to the wells labeled "Positive Control", and "Test Inhibitor". Incubate at room temperature for up to 30 minutes.
7. Measure absorbance at $\lambda = 410$ nm. Subtract "Blank" value from all other values. Note: If test compounds absorb at 410 nm, it is recommended to read the plate at time 0 as well as the final timepoint at 30 min. The time 0 measurement can be subtracted from the final reading to account for compound absorbance.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.


To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com



6042 Cornerstone Court W, Ste B
 San Diego, CA 92121
 Tel: 1.858.202.1401
 Fax: 1.858.481.8694
 Email: info@bpsbioscience.com

	Positive Control	Test Inhibitor	Blank
ARG Reaction Solution	90 µl	90 µl	90 µl
Test Inhibitor	-	5 µl	-
Inhibitor buffer (no inhibitor)	5 µl	-	5 µl
1X ARG Assay Buffer	-	-	5 µl
ARG1 (40 ng/µl)	5 µl	5 µl	-
Total	100 µl	100 µl	100 µl

Example of assay results:

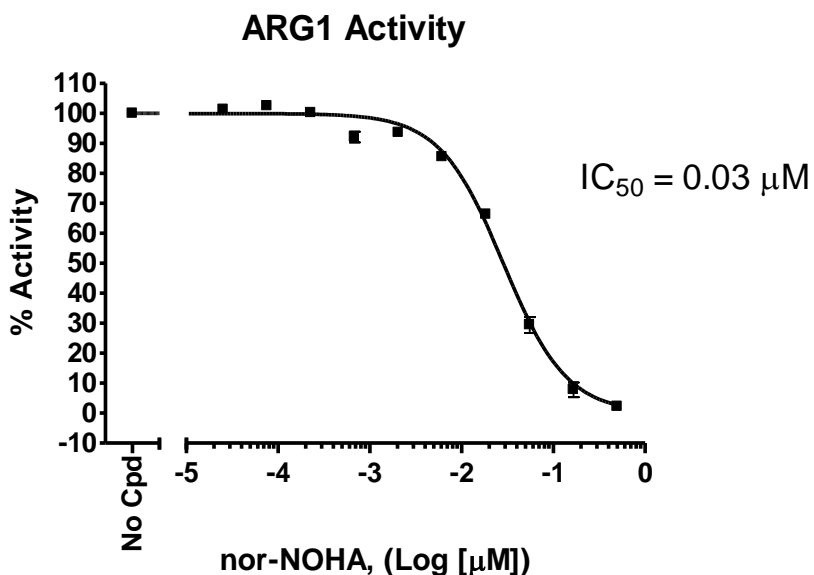

ARG1 inhibition by nor-NOHA (Cayman, 10006861) measured using the *ARG1 Inhibitor Screening Assay Kit*, BPS Bioscience, Catalog #72048-1. Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at info@bpsbioscience.com.

RELATED PRODUCTS:

<u>Product</u>	<u>Catalog#</u>	<u>Size</u>
ARG1	71658	50 µg
ARG2	71659	50 µg
ARG1 Inhibitor Screening Assay Kit	72048-2	384 reactions
ARG2 Inhibitor Screening Assay Kit	72043-2	384 reactions

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.202.1401** Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com