



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Description

The P70S6K (S6K1) Kinase Assay Kit is designed to measure P70S6K (S6K1) kinase activity for screening and profiling applications using ADP-Glo® as a detection reagent. The assay kit comes in a convenient 96-well format, with enough purified recombinant P70S6K (S6K1) kinase, kinase substrate, ATP and Kinase Assay Buffer for 100 enzyme reactions.

Background

p70S6K is responsible for the phosphorylation of the 40S ribosomal protein S6 and is ubiquitously expressed in human adult tissues. In intact cells p70S6K is activated by serum and this activation is inhibited by wortmannin and rapamycin. p70S6k activity is cell cycle-dependent and increases 20-fold in G1 cells released from G0. p70S6K activation requires sequential phosphorylation of serine residues in the putative auto-inhibitory pseudo-substrate domain, of threonine 389 phosphorylated by mTOR, and of threonine 229 phosphorylated by phosphoinositide-dependent kinase 1 (PDK-1).

Application(s)

Study enzyme kinetics and screen small molecular inhibitors for drug discovery and high throughput (HTS) applications.

Supplied Materials

Catalog #	Name	Amount	Storage
40062	P70S6K (S6K1), GST-Tag*	15 µg	-80°C
79334	Kinase Buffer 1 (5x)	1.5 ml	-20°C
79686	ATP (500 µM)	50 µl	-20°C
	S6K Substrate (1 mg/ml)	500 µl	-20°C
79696	White 96-well plate	1	Room Temperature

*The concentration of the protein is lot-specific and will be indicated on the tube

Materials Required but Not Supplied

Name	Catalog #
ADP-Glo® Kinase Assay	Promega #V6930
Microplate reader capable of reading luminescence	
Adjustable micropipettor and sterile tips	
30°C incubator	

Storage Conditions

This assay kit will perform optimally for up to 6 months from date of receipt when the materials are stored as directed.

Safety

This product is for research purposes only and not for human or therapeutic use. This product should be considered hazardous and is harmful by inhalation, in contact with skin, eyes, clothing, and if swallowed. If contact occurs, wash thoroughly.

Assay Principle

The **ADP-Glo™ Kinase Assay (Promega #V6930)** quantifies the amount of ADP produced by a kinase upon phosphorylation of a substrate. First, addition of the ADP-Glo™ reagent terminates the reaction and quenches the remaining ATP. Second, addition of the Kinase Detection reagent converts the produced ADP to ATP. The new ATP is quantified by a luciferase reaction. The luminescent signal correlates with the amount of ADP generated by the kinase and is linear to 1 mM ATP.

Contraindications

The final concentration of DMSO in the assay should not exceed 1%.

Assay Protocol

All samples and controls should be tested in duplicate.

1. Thaw **5x Kinase Assay Buffer 1**, (**500 μM**) **ATP**, and **S6K Substrate (1 mg/ml)**.
2. Prepare 3 ml of **1x Kinase Assay Buffer 1** by mixing 600 μl of **5x Kinase Assay Buffer 1** with 2,400 μl water.

Note: Three (3 ml) of 1x Kinase Assay Buffer 1 is sufficient for 100 reactions.

3. Prepare the **Master Mix** (12.5 μl/well): N wells x (6 μl of **1x Kinase Assay Buffer 1** + 0.5 μl of **ATP (500 μM)** + 5 μl of **S6K1 Substrate (1 mg/ml)** + 1 μl of distilled water. Add 12.5 μl to every well.
4. Prepare the **Test Inhibitor** (2.5 μl/well): for a titration, prepare serial dilutions at concentrations 10-fold higher than the desired final concentrations. The final volume of the reaction is 25 μl.

If the Test Inhibitor is water-soluble:

- 4.1 Prepare serial dilutions in the **1x Kinase Assay Buffer 1**, 10-fold more concentrated than the desired final concentrations.
- 4.2 For the positive and negative controls, use **1x Kinase Assay Buffer 1** (Diluent Solution).

Or

If the Test inhibitor is soluble in DMSO:

- 4.1 Prepare the test inhibitor at 100-fold the highest desired concentration in DMSO, then dilute the inhibitor 10-fold in **1x Kinase Assay Buffer 1** to prepare the highest concentration of the 10-fold intermediate dilutions. The concentration of DMSO is now 10%.
- 4.2 Prepare serial dilutions of the Test Inhibitor at 10-fold the desired final concentrations using 10% DMSO in **1x Kinase Assay Buffer 1** to keep the concentration of DMSO constant.
- 4.3 For positive and negative controls, prepare 10% DMSO in water (vol/vol) so that all wells contain the same amount of DMSO (Diluent Solution).

Note: The final concentration of DMSO should not exceed 1%.

5. Add 2.5 μl of **Test Inhibitor** to each well labeled "Test Inhibitor." For the "Positive Control" and "Blank," add 2.5 μl of **Diluent Solution** (either Kinase Assay Buffer or 10% DMSO in Kinase Assay Buffer, as described above).
6. To the wells designated as "Blank," add 10 μl of **1x Kinase Assay Buffer 1**.

7. Thaw **P70S6K (S6K1) kinase** on ice. Briefly spin the tube to recover its full contents. Dilute the protein kinase (10 µl/well) to 15 ng/µl using **1x Kinase Assay Buffer 1**.

Note: The concentration of protein is lot-specific and is indicated on the tube. Verify the initial concentration and dilute accordingly.

Note: This kinase is particularly sensitive to freeze/thaw cycles. Avoid multiple freeze/thaw cycles. Do not re-use the thawed protein and do not re-use the diluted kinase.

8. Initiate the reaction by adding 10 µl of diluted Kinase to the wells designated "Positive Control" and "Test Inhibitor."

Component	Blank	Positive Control	Test Inhibitor
Master Mix	12.5 µl	12.5 µl	12.5 µl
Test Inhibitor	-	-	2.5 µl
Diluent Solution	2.5 µl	2.5 µl	-
1x Kinase Assay Buffer 1	10 µl	-	-
P70S6K (S6K1) (15 ng/µl)	-	10 µl	10 µl
Total	25 µl	25 µl	25 µl

9. Incubate at 30°C for 45 minutes.
10. During the incubation, thaw the ADP-Glo™ reagent. At the end of the 45-minute reaction, add 25 µl of ADP-Glo™ reagent to each well. Cover the plate with aluminum foil and incubate at room temperature for 45 minutes.
11. Thaw the Kinase Detection Reagent. At the end of the 45-minute incubation, add 50 µl of Kinase Detection reagent to each well. Cover the plate with aluminum foil and incubate at room temperature for another 45 minutes.
12. Immediately read in a luminometer or a microplate reader capable of reading luminescence. The "Blank" value is subtracted from all other readings.

Reading Luminescence

Luminescence is the emission of light resulting from a chemical reaction. The detection of luminescence requires no wavelength selection because the method used is emission photometry and not emission spectrophotometry.

To properly read luminescence, make sure the plate reader is set for LUMINESCENCE mode. Typical integration time is 1 second, delay after plate movement is 100 msec. Do not use a filter when measuring light emission. Typical settings for the Synergy 2 BioTek plate reader: use the "hole" position on the filter wheel; Optics position: Top; Read type: endpoint. Sensitivity may be adjusted based on the luminescence of a control assay without enzyme (typically we set this value as 100).

Example Results

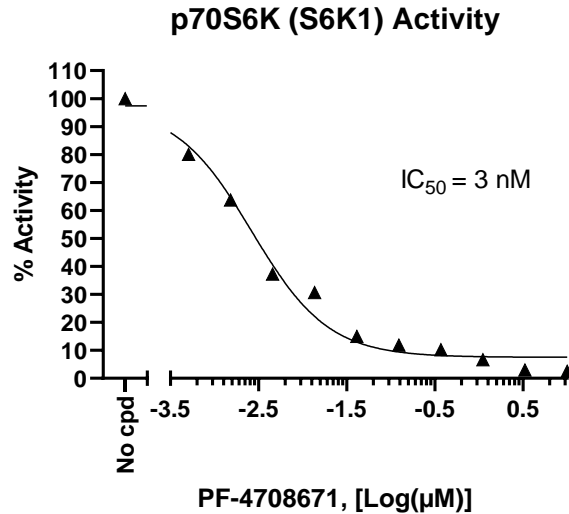


Figure 1: Inhibition of P70S6K (S6K1) kinase Activity by PF-4708671 (SelleckChem #S2163).
 The inhibition of P70S6K (S6K1) kinase activity was measured in the presence of increasing inhibitor concentrations using the P70S6K (S6K1) Kinase Assay Kit (BPS Bioscience #78806). The “Blank” value was subtracted from all other values. Results are expressed as the percent of control (kinase activity in the absence of inhibitor, set at 100%).

For lot-specific information, please contact BPS Bioscience, Inc. at support@bpsbioscience.com.

Troubleshooting Guide

Visit bpsbioscience.com/assay-kits-faq for detailed troubleshooting instructions. For all further questions, please email support@bpsbioscience.com.

Related Products

<i>Products</i>	<i>Catalog #</i>	<i>Size</i>
AXL Kinase Assay Kit	79711	96 reactions
BTK Assay Kit	79568	96 reactions
EPHA4 Kinase Assay Kit	78592	96 reactions
FGFR2 Assay Kit	79804	96 reactions