



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Data Sheet

### CD47, Fc fusion, Avi-Tag

Human, Recombinant, Fc fusion, C-Terminal  
Avi-tag

**Catalog #:** 79051

**Lot #:** 170517

**Conc.:** 0.8 mg/ml

**Formulated in:** 8 mM Phosphate, pH 7.4,  
110 mM NaCl, 2.2 mM KCl, and 20%  
glycerol

**Stability:** At least 6 months at  $-80^{\circ}\text{C}$ . *Avoid  
freeze/thaw cycles. Protein may be diluted to  
 $\geq 100 \mu\text{g/ml}$  in PBS + glycerol and stored at  
 $-80^{\circ}\text{C}$ .*

**Description:**

Human CD47, also known as leukocyte surface antigen CD47, Rh-related antigen, integrin-associated protein (IAP), antigenic surface determinant protein OA3, and MER6. GenBank Accession No. NM\_001777, a.a. 19-139 fused at the C-terminus to the Fc portion of human IgG1 and C-terminal Avi-tag expressed in a HEK293 cell expression system. MW = 42 kDa. This protein runs at a higher MW by SDS-PAGE due to glycosylation.

**References:**

1. Vernon-Wilson, E.F., *et al.*, *Eur. J. Immunol.* 2000; **30(8)**: 2130-2137.
2. Majeti, R. *et al.*, *Cell.* 2009; **138(2)**: 286-299.

---

### **Quality Assurance**


**4-20% SDS-PAGE  
Coomassie staining**

**Lane 1:**  
4 $\mu\text{g}$  CD47

**Lane 2:**  
Protein Marker

**MW:** 42 kDa+glycans

**Purity:**  $\geq 90\%$

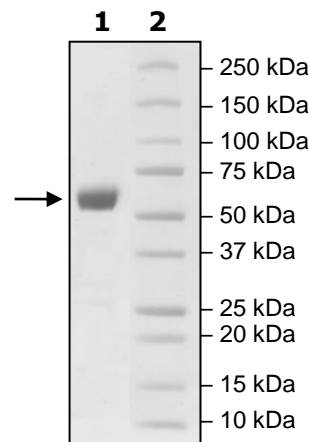

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**


Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

### Gel Filtration Curve

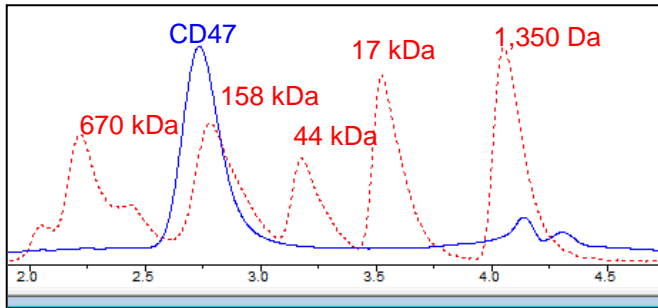

<10% Aggregation

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**  
Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)  
Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)