



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

FTO, His-Tag

Human, Recombinant, N-Terminal His-Tag

Catalog #: 79306

Lot #: 171012-A

Conc.: 1.55 mg/ml

Formulated in: 40 mM Tris-HCl pH 8.0, 110 mM NaCl, 2.2 mM KCl, 3 mM DTT and 20% glycerol.

Stability: At least 6 months at -80°C . Avoid freeze/thaw cycles. Protein may be diluted to $\geq 100 \mu\text{g/ml}$ in PBS + glycerol and remain stable for 6 months at -80°C .

References:


1. Frayling, T., *et al.*, *Science*. 2007; **316(5826)**: 889-894.
2. Dina, C., *et al.*, *Nature Genetics*. 2007; **39(6)**: 724-726.

Description: Human FTO, also known as Fat mass and obesity-associated protein, ALKBH9; BMIQ14; GDFD; KIAA1752, and alpha-ketoglutarate-dependent dioxygenase, GenBank Accession No. NM_001080432, a.a 32-505 (end), with N-terminal His-Tag and expressed in a *E. coli* cell expression system. MW = 56 kDa.

Specific Activity: 0.003 pmol/min/ μg
 Assay Conditions: Assay done using ELISA technology. 50 μl reaction mix containing HEPES-based buffer and 0-4 μg FTO is added to the wells coated with the m6A substrate. Incubate for 2 hr. at 37°C . Add antibody against N6-methyladenosine, incubate 1 hr. Then, add secondary HRP-labeled antibody and incubate 30 min. Finally, add HRP chemiluminescent substrates and read luminescence.

Applications: Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.


Quality Assurance

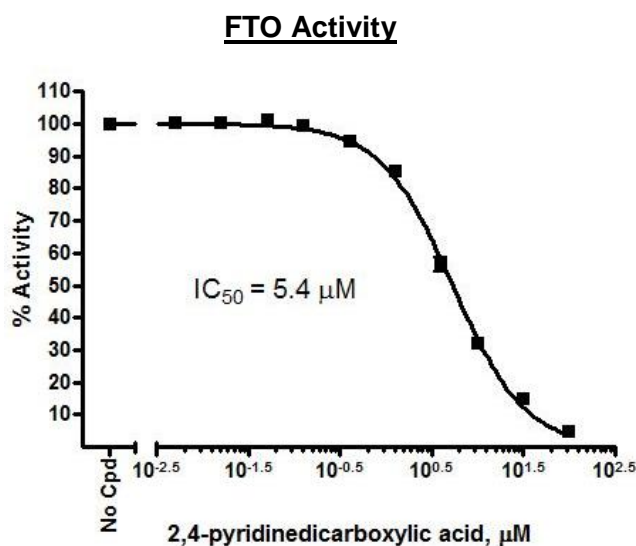

**4-20% SDS-PAGE
Coomassie staining**

Lane 1:
2 μg FTO
Lane 2:
Protein Marker

MW: 56 kDa
Purity: $\geq 90\%$

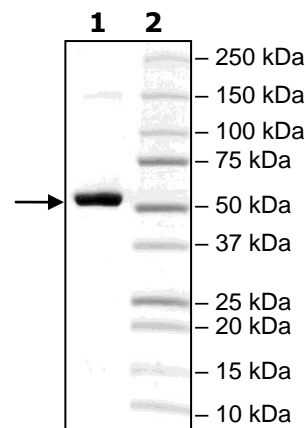

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com