



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

CERK, FLAG-Tag (Rat)

Rat, Recombinant, C-terminal FLAG-tag

Catalog #: 79414

Lot#: 180301-1 **Conc.:** 0.26 mg/ml

Formulated in: 40 mM Tris-HCl, pH 8.0, 110 mM NaCl, 2.2 mM KCl, 20% glycerol and 100ug/ml flag peptide

Stability: At least 6 months at -80°C . Avoid freeze/thaw cycles. Protein may be diluted to $\geq 100 \mu\text{g/ml}$ in PBS + glycerol and stored at -80°C .

References:

1. Sugiura, M. *et al.*, *JBC*. **277**, 23294-23300 (2002)
2. Majjalieh, S.M., *et al.*, *JBC*. **264**, 14354-14360 (1989)


Description: Rat CERK, also known as ceramide kinase, GenBank Accession No. NM_00113486, 1-531 (end) with C-terminal Flag tag, MW= 61 kDa, expressed in baculovirus expression system in Sf9 insect cells.

Specific Activity: ≥ 20 pmole/min/ μg

Assay conditions: CERK activity was measured at 37°C for 25 min in a buffer (20 mM HEPES pH7.2, 80 mM KCl, 3 mM CaCl_2 , 10 mM MgCl_2 , 2 mM DTT, 5 μM ATP) using a μ -Octylglucoside mixed micellar assay with 100 μM Ceramide (BPS# 40618). ATP reduction was detected using Kinase-Glo[®] Luminescent Kinase Assay Platform (Promega, Madison, WI).

Application: Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.

Quality Assurance

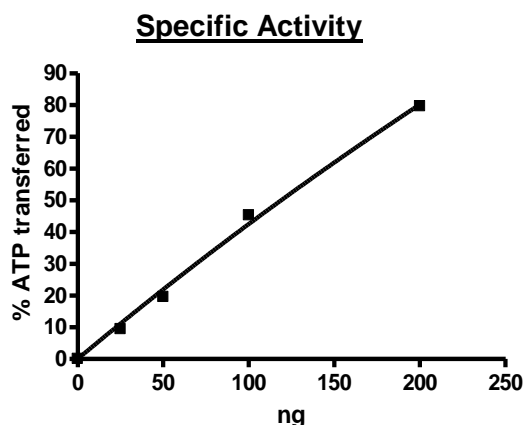


4-20% SDS-PAGE
Coomassie staining

Lane 1:
 Rat CERK

Lane 2:
 Protein Marker

Purity: 65%

MW: 61 kDa

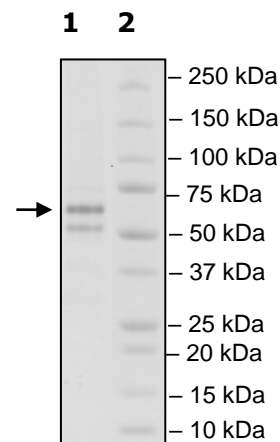

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com