



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

**Description**

The Chemi-Verse™ CDK9/CyclinT Kinase Assay Kit is designed to measure CDK9 (cyclin dependent kinase 9)/Cyclin T kinase activity for screening and profiling applications using ADP-Glo™ as a detection reagent. The assay kit comes in a convenient 384-well format, with enough purified recombinant CDK9/Cyclin T kinase, kinase substrate, ATP, and kinase assay buffer for 400 enzyme reactions.

**Background**

Cyclin-dependent kinase 9 (CDK9) is the catalytic subunit of the positive transcription elongation factor b (P-TEFb) complex, which phosphorylates the C-terminal domain of RNA polymerase II, a key player in the production of mature RNA. P-TEFb also includes cyclin T. CDK9 translocates to the nucleus once it gets auto phosphorylated, where it can associate with cyclin T. Hyperactivation of cyclins is known to result in cancer and the development of resistance to therapy, and CDK9-cyclin T is involved in the transcription of oncogenes. The inhibition or degradation of specific CDK-cyclin complexes could prove beneficial in cancer treatment, but the development of specific inhibitors has been difficult. Recently a small molecule degrader specific for CDK9-cyclin T was identified, LL-K9-3, which resulted in lower levels of androgen receptor (AR) and cMyc in 22RV1 cells. LL-K9-3 thus seemed more effective than SNS032, its parental CDK9 inhibitor, and the CDK9 PROTAC Thal-SNS032. Further research into new and more specific molecules will prove beneficial for cancer therapy.

**Applications**

Study enzyme kinetics and screen small molecule inhibitors for drug discovery and high throughput screening (HTS) applications.

**Supplied Materials**

Catalog #	Name	Amount	Storage
40307	CDK9/CyclinT, GST-Tag*	2 x 10 µg	-80°C
79334	5x Kinase Buffer 1	2 x 1.5 ml	-20°C
79686	500 µM ATP	200 µl	-20°C
79604	5x CDK Substrate Peptide 2	2 x 1 ml	-20°C
79969	White 384-well plate	1	Room Temperature

\*The concentration of the protein is lot-specific and will be indicated on the tube.

**Materials Required but Not Supplied**

Name	Ordering Information
ADP-Glo™ Kinase Assay	Promega #V6930
DTT (Dithiothreitol), 1M, optional	
Microplate reader capable of reading luminescence	
Adjustable micropipettor and sterile tips	
30°C incubator	

**Storage Conditions**

This assay kit will perform optimally for up to **6 months** from date of receipt when the materials are stored as directed.

**Safety**

This product is for research purposes only and not for human or therapeutic use. This product should be considered hazardous and is harmful by inhalation, in contact with skin, eyes, clothing, and if swallowed. If contact occurs, wash thoroughly.

**Assay Principle**

The **ADP-Glo™ Kinase Assay (Promega #V6930)** quantifies the amount of ADP produced by a kinase upon phosphorylation of a substrate. First, addition of the ADP-Glo™ reagent terminates the reaction and quenches the remaining ATP. Second, the addition of the Kinase Detection reagent converts the produced ADP to ATP. The newly generated ATP is quantified by a luciferase reaction. The luminescent signal correlates with the amount of ADP generated by the kinase and is linear to 1 mM ATP.

**Contraindications**

The final concentration of DMSO in the assay should not exceed 1%.

**Assay Protocol**

- All samples and controls should be tested in duplicate.
- The assay should include “Blank”, “Positive Control” and “Test Inhibitor” conditions.
- We recommend maintaining the diluted protein on ice during use.
- For detailed information on protein handling please refer to [Protein FAQs \(bpsbioscience.com\)](http://bpsbioscience.com).
- We recommend using Dinaciclib as internal control. If not running a dose response curve for the control inhibitor, we recommend running the control inhibitor at 0.1X, 1X and 10X the IC<sub>50</sub> value shown in the validation data below.

1. Thaw **5x Kinase Assay Buffer 1**, **500 μM ATP**, and **5x CDK Substrate Peptide 2**.

*Optional: If desired, make **5x Kinase Assay Buffer 1** with 10 mM DTT.*

2. Prepare 3 ml of **1x Kinase Assay Buffer 1** by mixing 600 μl of **5x Kinase Assay Buffer 1** with 2,400 μl of distilled water.

*Note: Three (3 ml) of **1x Kinase Assay Buffer 1** is sufficient for 100 reactions.*

3. Prepare a **Master Mix** (12.5 μl/well): N wells x (3 μl of 5x Kinase Assay Buffer 1 + 0.5 μl of 500 μM ATP + 5 μl of 5x CDK Substrate Peptide 2 + 4 μl of distilled water).
4. Add 12.5 μl of Master Mix to every well.
5. Prepare the **Test Inhibitor** (2.5 μl/well): for a titration prepare serial dilutions at concentrations 10-fold higher than the desired final concentrations. The final volume of the reaction is 25 μl.

5.1 If the Test Inhibitor is water-soluble: Prepare serial dilutions in 1x Kinase Assay Buffer 1, 10-fold more concentrated than the desired final concentrations.

For the positive and negative controls, use 1x Kinase Assay Buffer 1 (Diluent Solution).

**OR**

5.2 If the Test inhibitor is soluble in DMSO: Prepare the test inhibitor at 100-fold the highest desired concentration in 100% DMSO, then dilute the inhibitor 10-fold in 1x Kinase Assay Buffer 1 to prepare the highest concentration of the 10-fold intermediate dilutions. The concentration of DMSO is now 10%.

Prepare serial dilutions of the Test Inhibitor at 10-fold the desired final concentrations using 10% DMSO in 1x Kinase Assay Buffer 1 to keep the concentration of DMSO constant.

For positive and negative controls, prepare 10% DMSO in 1x Kinase Assay Buffer 1 (vol/vol) so that all wells contain the same amount of DMSO (Diluent Solution).

*Note: The final concentration of DMSO should not exceed 1%.*

6. Add 2.5 µl of Test Inhibitor to each well labeled "Test Inhibitor".
7. Add 2.5 µl of Diluent Solution to the "Positive Control" and "Blank" wells.
8. Add 10 µl of 1x Kinase Assay Buffer 1 to the "Blank" wells.
9. Thaw **CDK9/CyclinT Kinase** on ice. Briefly spin the tube to recover its full content.
10. Dilute the protein kinase (10 µl/well) to 5 ng/µl with 1x Kinase Assay Buffer 1.
11. Initiate the reaction by adding 10 µl of diluted Kinase to the wells designated "Positive Control" and "Test Inhibitor".

<b>Component</b>	<b>Blank</b>	<b>Positive Control</b>	<b>Test Inhibitor</b>
Master Mix	12.5 µl	12.5 µl	12.5 µl
Test Inhibitor	-	-	2.5 µl
Diluent Solution	2.5 µl	2.5 µl	-
1x Kinase Assay Buffer 1	10 µl	-	-
Diluted CDK9/CyclinT (5 ng/µl)	-	10 µl	10 µl
<b>Total</b>	<b>25 µl</b>	<b>25 µl</b>	<b>25 µl</b>

12. Incubate at 30°C for 45 minutes.
13. Thaw the ADP-Glo™ reagent.
14. At the end of the 45 minute reaction, add 25 µl of ADP-Glo™ reagent to each well.
15. Cover the plate with aluminum foil and incubate at Room Temperature (RT) for 45 minutes.
16. Thaw the Kinase Detection Reagent.


17. Add 50 µl of Kinase Detection reagent to each well.
18. Cover the plate with aluminum foil and incubate at RT for another 45 minutes.
19. Immediately read in a luminometer or a microplate reader capable of reading luminescence.
20. The “Blank” value is subtracted from all other readings.

### Reading Luminescence

Luminescence is the emission of light resulting from a chemical reaction. The detection of luminescence requires no wavelength selection because the method used is emission photometry and not emission spectrophotometry.

To properly read luminescence, make sure the plate reader is set for LUMINESCENCE mode. Typical integration time is 1 second, delay after plate movement is 100 msec. Do not use a filter when measuring light emission. Typical settings for the Synergy 2 BioTek plate reader: use the “hole” position on the filter wheel; Optics position: Top; Read type: endpoint. Sensitivity may be adjusted based on the luminescence of a control assay without enzyme (typically we set this value as 100).

### Example Results

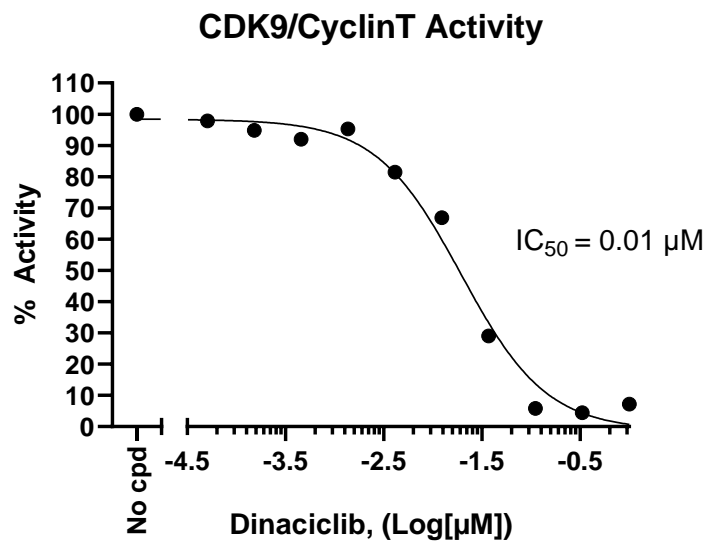

*Figure 1: Inhibition of CDK9/CyclinT kinase activity by Dinaciclib.*

CDK9/CyclinT kinase activity was measured in the presence of increasing concentrations of Dinaciclib (SelleckChem S2768). The “Blank” value was subtracted from all other values. Results are expressed as the percent of control (kinase activity in the absence of inhibitor, set at 100%).

*Data shown is representative. For lot-specific information, please contact BPS Bioscience, Inc. at [support@bpsbioscience.com](mailto:support@bpsbioscience.com).*

### Troubleshooting Guide

Visit [bpsbioscience.com/assay-kits-faq](https://bpsbioscience.com/assay-kits-faq) for detailed troubleshooting instructions. For all further questions, please email [support@bpsbioscience.com](mailto:support@bpsbioscience.com)

**References**

Napolitano G, *et al.*, 2003 *J Cell Physiol.* 197(1):1-7.

Li J., *et al.*, 2022 *J Med Chem.* 65(16):11034-11057.

**Related Products**

<i>Products</i>	<i>Catalog #</i>	<i>Size</i>
CDK9/CyclinK, GST-tag Recombinant	40106	10 µg
Chemi-Verse™ CDK9/CyclinK Kinase Assay Kit	82237	96 reactions
Chemi-Verse™ CDK3/CyclinE1 Kinase Assay Kit	78884	96 reactions
Chemi-Verse™ CDK17/CyclinY Kinase Assay Kit	78885	96 reactions
Chemi-Verse™ CDK8/CyclinC Kinase Assay Kit	78886	96 reactions
Chemi-Verse™ CDK2/CyclinE1 Kinase Assay Kit	82228	96 reactions

*Version 010524*