



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

## **Data Sheet**

### **TRPC7-HEK293 Recombinant Cell line Cat #: 90030**

#### **Description**

Recombinant HEK293 cell line expressing human TRPC7 (transient receptor potential cation channel, subfamily C, member 7, accession number NM\_020389).

#### **Host cell**

HEK293 cells

#### **Format**

Each vial contains  $\sim 1.5 \times 10^6$  cells in 1 ml of 10% DMSO.

#### **Storage**

Immediately upon receipt, store in liquid nitrogen.

#### **Mycoplasma testing**

The cell line has been screened using the PCR-based Venor<sup>®</sup>GeM Mycoplasma Detection kit (Sigma-Aldrich) to confirm the absence of Mycoplasma species.

#### **Introduction**

TRPC7 channel belongs to the transient receptor potential channel (TRP) superfamily which is divided into seven main subfamilies. TRP channels share a similar architecture of six-transmembrane domains and function as tetramers. Most channels are permeable to calcium and monovalent cations.

Specifically, TRPC7 channel belongs to a canonical TRP (TRPC) subfamily that is involved in calcium entry. TRPC7 is activated via diacylglycerol (DAG) after receptor-mediated activation of PLC.

#### **Applications**

- Monitor TRPC7 calcium channel activity
- Screen for activators or inhibitors of TRPC7 calcium channel

#### **Functional validation**

N<sup>′</sup>- terminal FLAG tagged human TRPC7 channel has been stably expressed in HEK293 cell line and its expression was confirmed by Western blotting.

The function of TRPC7 was characterized by calcium assay. In TRPC7-HEK293 cells, addition of methacholine, an activator of phospholipase C-linked muscarinic receptors, significantly activates TRPC7, resulting Ca<sup>2+</sup> entry in the cells.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [\*\*info@bpsbioscience.com\*\*](mailto:info@bpsbioscience.com)


Please visit our website at: [\*\*www.bpsbioscience.com\*\*](http://www.bpsbioscience.com)

These data show the stable expression of TRPC7 channel in HEK293 cells.

**Figure 1. Methacholine-activated Ca<sup>2+</sup> entry in TRPC7 expressed HEK293 cells.**

A) TRPC7-HEK293 cells; B) WT-HEK293 cells.

Cells were loaded with fluorescent Ca<sup>2+</sup> indicator Fluo-8 and incubated in the absence of extracellular Ca<sup>2+</sup>. Cells were treated with 300 μM methacholine (MeCh) (black line), followed by 2 mM of Ca<sup>2+</sup> was added where indicated. Cells without methacholine (pink line) were included as a control. 5 μM Gd<sup>3+</sup> was maintained throughout the experiment, to block endogenous store depletion-induced Ca<sup>2+</sup> entry. The calcium influx was measured by Fluo-8 fluorescence (excitation 485±20 nm and emission 528±20 nm).

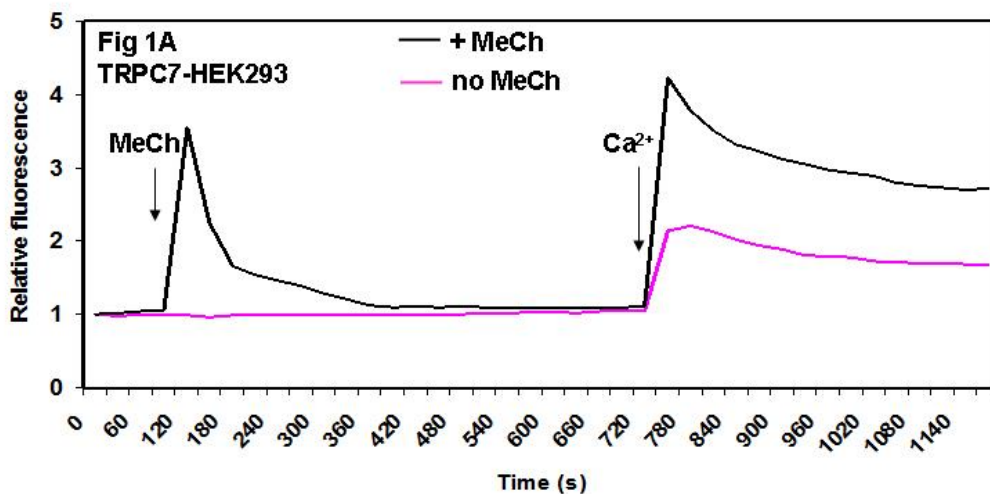


OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)

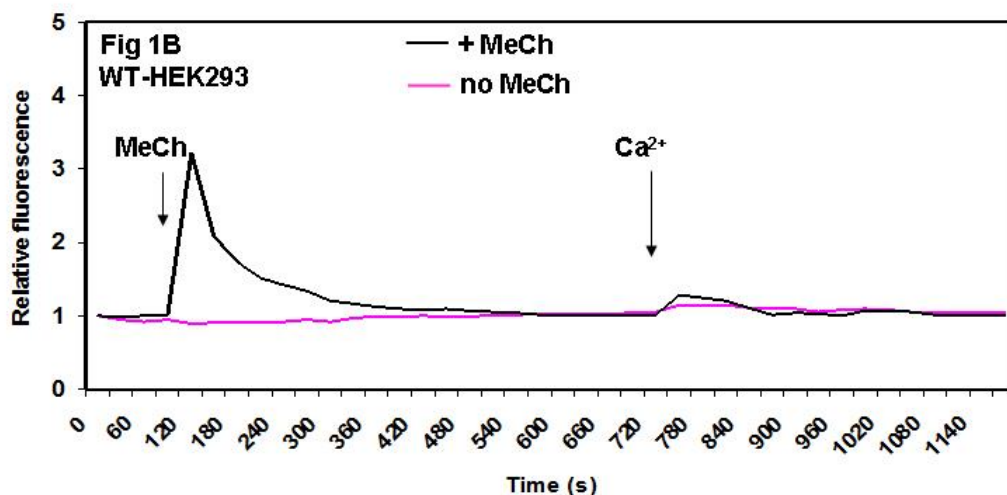

### Culture Conditions

**Thaw Medium 1 (BPS Cat. #60187):** MEM medium (Hyclone #SH30024.01) supplemented with 10% FBS (Invitrogen #26140-079), 1% non-essential amino acids (Hyclone #SH30238.01), 1 mM Na pyruvate (Hyclone #SH30239.01), 1% Penicillin/Streptomycin (Hyclone SV30010.01)

**Complete Growth Medium:** Thaw Medium 1 (BPS Cat. #60187) plus 400 µg/ml of Geneticin (Life Technologies #11811031)

Cells should be grown at 37°C with 5% CO<sub>2</sub> using complete growth medium. It may be necessary to adjust the percentage of CO<sub>2</sub> in the incubator depending on the NaHCO<sub>3</sub> level in the basal medium.

**To thaw the cells,** it is recommended to quickly thaw the frozen cells from liquid nitrogen in a 37°C water-bath, transfer to a tube containing 10 ml of Thaw Medium 1 (**no Geneticin**), spin down cells, and resuspend cells in pre-warmed Thaw Medium 1 (**no Geneticin**). Transfer resuspended cells to a T25 flask and culture in a 37°C CO<sub>2</sub> incubator. At first passage, switch to complete growth medium (**contains Geneticin**). Cells should be split before they reach complete confluence.

**To passage the cells,** rinse cells with phosphate buffered saline (PBS), detach cells from culture vessel with 0.05% Trypsin/EDTA, and add Thaw Medium 1 containing Geneticin. Transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cell suspension into new culture vessels. Subcultivation ration: 1:10 to 1:20 weekly.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)



6042 Cornerstone Court W, Ste B  
San Diego, CA 92121  
**Tel:** 1.858.829.3082  
**Fax:** 1.858.481.8694  
**Email:** [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

### Vector and sequence

N-terminal FLAG tagged human TRPC7 was cloned into pIRES-neo expression vector.  
Polylinker: CMV-EcoRV-NheI-**TRPC7**-BamHI-NotI-BstXI-IRES-neomycin<sup>R</sup>

hTRPC7 sequence (accession number NM\_020389)

MLRNSTFKNMQRRHTTLREKGRRAIRGPAYMFNEKGTSLTPEEERFLDSA EYGNIPV  
VRKMLEESKTLNFCVDYMGQNALQLAVGNEHLEVTELLKKNLARVGDALLLAISKG  
YVRIVEAILNHPAFAAGQRLTLSPLEQELRDDDFYAYDEDGTRFSDITPIILAAHCQEYE  
IVHILLKLGARIERPHDYFCKCNECTEKQRKDSFSHSRMRMAYKGLASAA YLSLSSD  
PVLTALELSNELARLANIETEFKNDYRKL SMQCKDFVVGVLDCRDTEEVEAILNGDVNF  
QVWSDHHRPSLSRIKLAIKYEVKKFVAHPNCQQQLTMWYENLSGLRQQSIAVKFLAVF  
GVSIGLPFLAIAYWIAPCSKLGRTL RSPFMKFVAHAVSFTIFLGLLVNASDRFEGVKTLP  
NETFTDYPKQIFRVKTTQFSWTEMLIMKWV LGMWSECKEIWEEGPREYVHLWNLLD  
FGMLSIFVASFTARFMAFLKATEAQLYVDQH VQDDTLHNVS LPPPEVAYFTYARDKWWP  
SDPQIISEGLYAI AVVLSFSRIAYILPANESFGPLQISLGRTVKDIFKFMVIFIMVFVAFMIG  
MFNLYSYRGAKYNPAFTTVEESFKTLFWSIFGLSEVISVVLKYD HKFIENIGYVLYGVY  
NVTMVVLLNMLIAMINNSYQEIEEDADVEWKFARAKLWLSYFDEGRTL PAFPNLVPSP  
KSFYYLIMRIKMCLIKLCKSKAKSCENDLEMGM LNSKFKKTRYQAGMRNSENLTANNTL  
SKPTRYQKIMKRLIKRYVLKAQVDRENDEVNEGELKEIKQDISSRLRYELLEEK SQATGEL  
ADLIQQLSEKFGKLNKDH LRVNKGKDI

### References

1. Moran MM *et al* (2011) Transient receptor potential channels as therapeutic targets. *Nat Rev Drug Discov.* **10(8)**:601-620.
2. Vassort G. and Alvarez J. (2009) Transient receptor potential: a large family of new channels of which several are involved in cardiac arrhythmia. *Can. J. Physiol. Pharmacol.* **87**: 100-107.
3. Lievreumont J., Bird G.J., and Putney, Jr. J.W. (2004) Canonical transient receptor potential TRPC7 can function as both a receptor- and store-operated channel in HEK-293 cells. *Am J Physiol Cell Physiol* **287**: C1709-C1716.

### License Disclosure

Purchase of this cell line grants you with a 10-year license to use this cell line in your immediate laboratory, for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make the cell line available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit the use of this cell line in humans or for therapeutic or drug use. Inappropriate use or distribution of this cell line will result in revocation of the license and result in an immediate cease of sales and distribution of BPS products to your laboratory. Publications using this cell line should reference BPS Bioscience, Inc., San Diego.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082** Fax **1.858.481.8694**

Or you can Email us at: [info@bpsbioscience.com](mailto:info@bpsbioscience.com)

Please visit our website at: [www.bpsbioscience.com](http://www.bpsbioscience.com)