



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

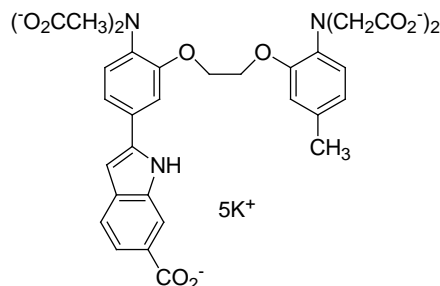
F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT AND SAFETY DATA SHEET

PRODUCT NAME: Indo-1, pentapotassium salt**CATALOG #:** 50041**MOLECULAR INFORMATION:** C₃₂H₂₆K₅N₃O₁₂
MWt: 840
[132319-21-9]**PROPERTIES:**

Color & Form	Light yellow solid
Purity	≥ 90% by HPLC
Solubility	Soluble in water (pH>6) or DMSO
Absorption/Emission	349 nm/475 nm (no Ca ²⁺); 330 nm/400 nm (high Ca ²⁺)
Extinction Coefficient	33,000 M ⁻¹ cm ⁻¹ (no Ca ²⁺); 33,000 M ⁻¹ cm ⁻¹ (high Ca ²⁺)

STORAGE AND HANDLING:

Store desiccated at 4°C. Protect from light, especially when in solution

APPLICATION:

Similar to fura-2, indo-1 is also a UV-excitable fluorescent Ca²⁺ sensor. However, different from fura-2, the fluorescent emission maximum undergoes a large blue shift from 482 nm to 398 nm upon Ca²⁺ binding. Thus, Ca²⁺ concentration can be determined by ratioing the fluorescence intensities at the two wavelengths. As with fura-2, this ratioing technique avoids problems associated with uneven dye distribution, cell or tissue thickness and photobleaching. Indo-1 has been widely used in flow cytometry studies.

Indo-1, pentapotassium salt is membrane-impermeant but can be loaded into cells via microinjection or scrape loading.

Ref.: Babcock, D.R., et al. *J. Biol. Chem.* **262**, 15041(1987).

TOXICITY: Unknown

FIRST AID:	Potentially harmful. Avoid prolonged or repeated exposure. Avoid getting in eyes, on skin, or on clothing. Wash thoroughly after handling. If eye or skin contact occurs, wash affected areas with plenty of water for 15 minutes and seek medical advice. In case of inhaling or swallowing, move individual to fresh air and seek medical advice immediately.
-------------------	-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

Disclaimer: Materials from Biotium are sold for research use only, and are not intended for food, drug, household, or cosmetic use. Biotium is not liable for any damage resulting from handling or contact with this product.
