



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

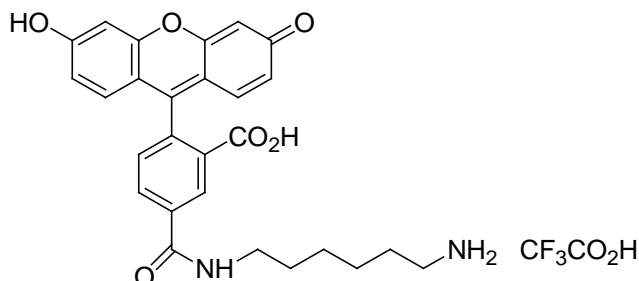
F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT AND SAFETY DATA SHEET

PRODUCT NAME: Fluorescein Hexylamine**CATALOG #:** 91058**MOLECULAR INFORMATION:** C₂₉H₂₇F₃N₂O₈
Mwt: 588.54**PROPERTIES:**

Color & Form Off-white solid
Purity ≥ 95% by TLC
Solubility Soluble in pH >6 or DMF or DMSO
λ_{ex}/λ_{em} 493/517nm
Extinction Coefficient 80,000

STORAGE AND HANDLING: Store desiccated at ≤ 4°C and protect from light**APPLICATION:** Fluorescein hexylamine is a fluorescein derivative with a primary amine group. The dye can be used to react with activated carboxylic acids, isothiocyanates, sulfonyl chlorides and other amine-reactive groups. This product could be used as a direct replacement to fluorescein cadaverine (**cat#: 92000**).**TOXICITY:** Unknown.**FIRST AID:** Potentially harmful. Avoid prolonged or repeated exposure. Avoid getting in eyes, on skin, or on clothing. Wash thoroughly after handling. If eye or skin contact occurs, wash affected areas with plenty of water for 15 minutes and seek medical advice. In case of inhaling or swallowing, move individual to fresh air and seek medical advice immediately.**Disclaimer:** Materials from Biotium are sold for research use only, and are not intended for food, drug, household, or cosmetic use. Biotium is not liable for any damage resulting from handling or contact with this product.