



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



## Mouse Anti-DYKDDDDK IgG conjugated to R-Phycoerythrin

Product Number D5-1718  
Amount 100 µg total protein  
Store at 4°C

### Form/ Storage

Supplied as a lyophilized powder. Upon receipt, store at 2-8°C in the dark. Phycobiliproteins are sensitive to freeze-thaw cycles: after reconstitution, store at 2-8°C in the dark – do not freeze.

### Handling

Avoid exposure to heat and light. Prior to use reconstitute to 1 ml with distilled deionized water, vortex and allow it to sit on ice for 20 minutes.

### Buffer


Upon reconstitution, the product is in 100 mM sodium phosphate (pH 7.4), 50 mM sucrose, 150 mM sodium chloride, 0.1% BSA as a stabilizer, and 2 mM sodium azide as a preservative. The concentration of the conjugate is 100 µg/ml.

### Stability

Lyophilized material is stable for one year. After product has been reconstituted, product should be stored at 2-8°C in the dark and be used within 3 months.

### Note

For research use only, not for diagnostic or therapeutic use.

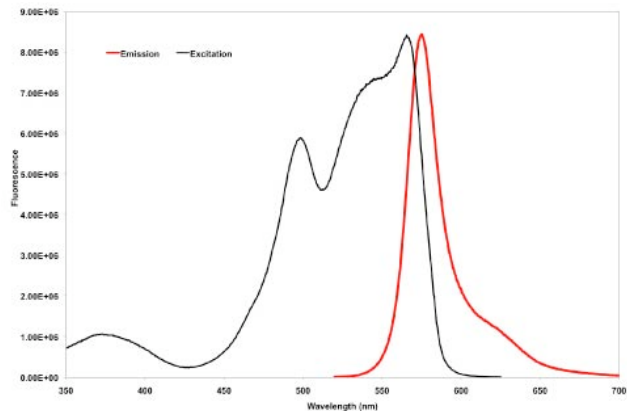

Fluorescence excitation and emission spectra of R-phycoerythrin in 100 mM sodium phosphate (pH 7.2) + 1 mM EDTA and 1 mM sodium azide. Emission scan was taken with excitation at 498 nm. Excitation scan was taken with emission at 575 nm. Scans were normalized to equalize peak heights.

### Spectral Characteristics

Visible absorption maxima 565>540>498  
Emission maximum 578 nm

### Concentration

After reconstitution to 1.0 ml

0.1 mg/mL

Fluor:Protein = ~1.5:1

