



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



tetranor-PGEM

Item No. 14840

CAS Registry No.: 24769-56-0
Formal Name: 9,15-dioxo-11 α -hydroxy-13,14-dihydro-2,3,4,5-tetranor-prostan-1,20-dioic acid
Synonym: tetranor-Prostaglandin E Metabolite
MF: C₁₆H₂₄O₇
FW: 328.4
Purity: \geq 98%
Stability: \geq 1 year at -80°C
Supplied as: A solution in methyl acetate

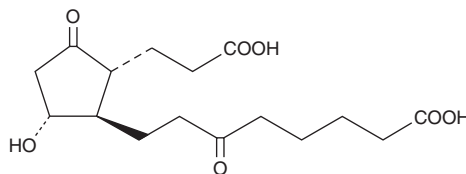

Laboratory Procedures

For long term storage, we suggest that tetranor-PGEM be stored as supplied at -80°C. It should be stable for at least one year.

tetranor-PGEM is supplied as a solution in methyl acetate. To change the solvent, simply evaporate the methyl acetate under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. The solubility of tetranor-PGEM in these solvents is approximately 50 mg/ml.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant since organic solvents may have physiological effects at low concentrations. If an organic solvent-free solution of tetranor-PGEM is needed, it can be prepared by evaporating the methyl acetate and directly dissolving the neat oil in aqueous buffers. The solubility of tetranor-PGEM in PBS (pH 7.2) is approximately 1 mg/ml. Store aqueous solutions of tetranor-PGEM on ice and use within 12 hours of preparation. We recommend making a fresh preparation each day.

Description

tetranor-PGEM is the major urinary metabolite of prostaglandin E₁ (PGE₁) and PGE₂, and is used as a urinary marker of PGE₂ biosynthesis.^{1,2} About 15% of an infused dose of PGE₂ appears as this metabolite in the urine of humans. Normal healthy males excrete 7-40 μ g of tetranor-PGEM over a 24-hour period.¹

References

1. Hamberg, M. Inhibition of prostaglandin synthesis in man. *Biochem. Biophys. Res. Commun.* **49**, 720-726 (1972).
2. Honda, H., Fukawa, K., and Sawabe, T. Influence of adjuvant arthritis on main urinary metabolites of prostaglandin F and E rats. *Prostaglandins* **19**, 259-269 (1980).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 03/03/2016

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM