



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



Calcein AM

Item No. 14948

CAS Registry No.: 148504-34-1

Formal Name: N,N'-[[3',6'-bis(acetyloxy)-3-oxospiro[isobenzofuran-1(3H),9'-[9H]xanthene]-2',7'-diyl]bis(methylene)]bis[N-[2-[(acetyloxy)methoxy]-2-oxoethyl]]-glycine, 1,1'-bis[(acetyloxy)methyl] ester

Synonyms:

Calcein Acetoxymethyl ester, NSC 689290

MF: $C_{46}H_{46}N_2O_{23}$

FW: 994.9

Purity: $\geq 90\%$


Ex./Em. Max: 494/517 nm

Supplied as: A crystalline solid

Storage: -20°C

Stability: ≥ 2 years

Special Conditions: Light sensitive

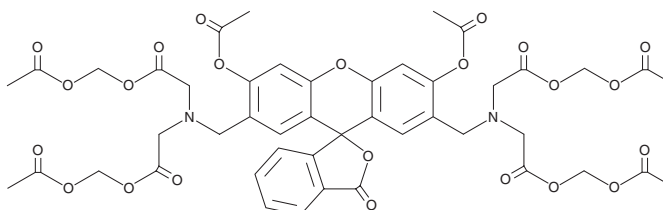

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Calcein AM is supplied as a crystalline solid. A stock solution may be made by dissolving the calcein AM in the solvent of choice, which should be purged with an inert gas. Calcein AM is soluble in organic solvents such as ethanol, methanol, DMSO, and dimethyl formamide. The solubility of calcein AM in these solvents is approximately 10 mg/ml.

Description

Calcein AM is a cell-permeable dye that, upon entering live cells, is cleaved by intracellular esterases, leaving membrane-impermeant calcein (Item No. 16221), a fluorescent indicator with absorption and emission maxima of 494 and 517 nm, respectively. It can be used as an indicator of cell viability, cell-cell communication, cytotoxicity, or changes in intracellular calcium, fluoride, iron, or mercury.^{1,2} As calcein is a substrate for multidrug resistance-associated protein 1 (MRP1) and multidrug resistance protein 3 (MDR3, P-glycoprotein 3), calcein AM is used to study these proteins and their regulation.³⁻⁵

References

1. Hynes, J., Floyd, S., Soini, A.E., *et al.* Fluorescence-based cell viability screening assays using water-soluble oxygen probes. *J. Biomol. Screen.* **8**, 264-272 (2003).
2. Sabnis, R.W. Handbook of biological dyes and stains: Synthesis and industrial applications. John Wiley & Sons, Inc., Hoboken, NJ, USA (2010).
3. Polli, J.W., Wring, S.A., Humphreys, J.E. *et al.* Rational use of *in vitro* P-glycoprotein assays in drug discovery. *J. Pharmacol. Exp. Ther.* **299**(2), 620-628 (2001).
4. Tiberghien, F. and Loor, F. Ranking of P-glycoprotein substrates and inhibitors by a calcein-AM fluorometry screening assay. *Anticancer Drugs* **7**(5), 568-578 (1996).
5. Szakics, G., Jakab, K., Antal, F., *et al.* Diagnostics of multidrug resistance in cancer. *Pathol. Oncol. Res.* **4**(4), 251-257 (2016).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 12/09/2020

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM