



**SZABO  
SCANDIC**

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](http://linkedin.com/company/szaboscandic)



# PRODUCT INFORMATION



## Enniatin B

Item No. 15382

CAS Registry No.: 917-13-5

Formal Name: cyclo[(2R)-2-hydroxy-3-methylbutanoyl-N-methyl-L-valyl-(2R)-2-hydroxy-3-methylbutanoyl-N-methyl-L-valyl-(2R)-2-hydroxy-3-methylbutanoyl-N-methyl-L-valyl]

Synonym: Antibiotic 86/88

MF: C<sub>33</sub>H<sub>57</sub>N<sub>3</sub>O<sub>9</sub>

FW: 639.8


Purity: ≥95%

Supplied as: A powder

Storage: -20°C

Stability: ≥2 years

Item Origin: Bacterium/Fusarium sp.

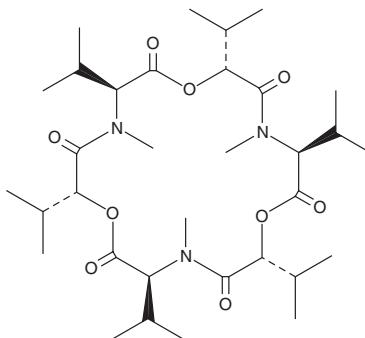

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

## Laboratory Procedures

Enniatin B is supplied as a powder. A stock solution may be made by dissolving the enniatin B in the solvent of choice. Enniatin B is soluble in organic solvents such as ethanol, methanol, DMSO, and dimethyl formamide, which should be purged with an inert gas.

## Description

Enniatins are cyclichexadepsipeptides commonly isolated from fungi. Many act as ionophores, forming pores in cellular membranes to allow selective ion transport.<sup>1</sup> Enniatin B is a relatively poor ionophore with some capacity to facilitate import of K<sup>+</sup> and Na<sup>+</sup> across membranes.<sup>2</sup> It inhibits the pleiotropic drug resistance protein 5 (Pdr5p) in yeast.<sup>3</sup> Through this mechanism enniatin B, at concentrations as low as 0.8 μM, augments the ability of cerulenin (Item No. 10005647) or cycloheximide (Item No. 14126) to impair cell proliferation in cells overexpressing Pdr5p, an effect that is not observed in cells lacking Pdr5p.<sup>3</sup> Like other enniatins, enniatin B inhibits acyl-CoA: cholesterol acyltransferase (IC<sub>50</sub> = 113 μM), blocking cholesterol ester formation.<sup>4</sup> Enniatin B (1 μM) also increases caspase activity and induces apoptosis in H4IE hepatoma cells and, when mixed with other enniatins, alters p53 signaling in human cancer cells.<sup>5,6</sup>

## References

1. Sy-Cordero, A.A., Pearce, C.J., and Oberlies, N.H. *J. Antibiot. (Tokyo)* **65**(11), 541-549 (2012).
2. Kamyar, M.R., Rawnduzi, P., Studenik, C.R., et al. *Arch. Biochem. Biophys.* **429**(2), 215-223 (2004).
3. Hiraga, K., Yamamoto, S., Fukuda, H., et al. *Biochem. Biophys. Res. Commun.* **328**(4), 1119-1125 (2005).
4. Tomoda, H., Huang, X.H., Cao, J., et al. *J. Antibiot. (Tokyo)* **45**(10), 1626-1632 (1992).
5. Wätjen, W., Debbab, A., Hohlfeld, A., et al. *Mol. Nutr. Food Res.* **53**(4), 431-440 (2009).
6. Dornetshuber, R., Heffeter, P., Kamyar, M.R., et al. *Chem. Res. Toxicol.* **20**(3), 465-473 (2007).

### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

## CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 - USA

PHONE: [800] 364-9897  
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM