



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

# PRODUCT INFORMATION



## Resolvin D3-d<sub>5</sub> Item No. 19512

**Formal Name:** 4S,11R,17S-trihydroxydocosa-5Z,7E,9E,13Z,15E,19Z-hexaenoic-21,21,22,22,22-d<sub>5</sub> acid

**Synonym:** RVD3-d<sub>5</sub>

**MF:** C<sub>22</sub>H<sub>27</sub>D<sub>5</sub>O<sub>5</sub>

**FW:** 381.5

**Chemical Purity:** ≥98% (Resolvin D3)


**Deuterium Incorporation:** ≥99% deuterated forms (d<sub>1</sub>-d<sub>5</sub>); ≤1% d<sub>0</sub>

**UV/Vis.:** λ<sub>max</sub>: 240, 262, 273, 283 nm

**Supplied as:** A solution in ethanol

**Storage:** -80°C

**Stability:** ≥1 year

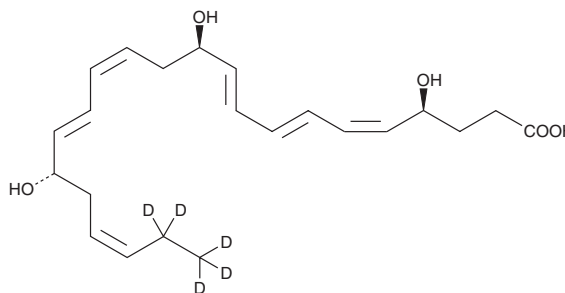

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

### Laboratory Procedures

Resolvin D3-d<sub>5</sub> (RVD3-d<sub>5</sub>) is intended for use as an internal standard for the quantification of RVD3 (Item No. 13834) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

RVD3-d<sub>5</sub> is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. It is recommended that this product be stored and handled in an ethanol solution. Resolvins can isomerize and degrade when put into freeze thaw conditions and/or in solvents such as dimethyl formamide or DMSO.

### Description

RVD3 is a DHA-derived product first identified in mouse inflammatory exudates.<sup>1</sup> It reduces neutrophil infiltration *in vivo* in both mouse peritonitis and dermal inflammation. In addition to suppressing leukocyte migration, RVD3 enhances macrophage phagocytosis and efferocytosis.<sup>2</sup> Unlike other resolvins, RVD3 appears late in resolution in mouse peritonitis.<sup>2,3</sup>

### References

1. Serhan, C.N., Hong, S., Gronert, K., *et al.* *J. Exp. Med.* **196(8)**,1025-1037 (2002).
2. Dalli, J., Winkler, J.W., Colas, R.A., *et al.* *Chem. Biol.* **20(2)**, 188-201 (2015).
3. Serhan, C.N. *Nature* **510(7503)**, 92-101 (2014).

#### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

#### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

#### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company,06/18/2019

#### CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA

**PHONE:** [800] 364-9897  
[734] 971-3335

**FAX:** [734] 971-3640

CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM