



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



α -Carotene

Item No. 19772

CAS Registry No.: 7488-99-5
Formal Name: (6'R)- β,ϵ -carotene
Synonym: all-*trans*- α -Carotene
MF: C₄₀H₅₆
FW: 536.9
Purity: $\geq 95\%$
UV/Vis.: λ_{\max} : 241, 281 nm
Supplied as: A solid
Storage: -20°C
Stability: ≥ 2 years

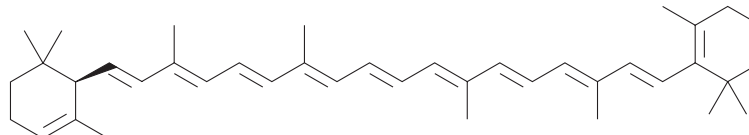

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

α -Carotene is supplied as a solid. A stock solution may be made by dissolving the α -carotene in the solvent of choice. α -Carotene is soluble in organic solvents such as DMSO and dimethyl formamide, which should be purged with an inert gas. The solubility of α -carotene in these solvents is approximately 1 and 0.1 mg/ml, respectively.

α -Carotene is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, α -carotene should first be dissolved in DMSO and then diluted with the aqueous buffer of choice. α -Carotene has a solubility of approximately 0.1 mg/ml in a 1:10 solution of DMSO:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

α -Carotene is a precursor of vitamin A (Item No. 20241) that has been found in various fruits and vegetables.¹ It inhibits proliferation of GOTO human neuroblastoma cells more potently than β -carotene (Item No. 16837) and halts the cell cycle at the G₀/G₁ phase concomitantly with a reduction in the mRNA expression of the protooncogene N-Myc.² It is also more potent than β -carotene in mouse models of skin and lung carcinogenesis and decreases the number of hepatomas in mice with spontaneous liver carcinogenesis when administered in drinking water at a concentration of 0.05%.³ α -Carotene levels are increased in patients with coronary heart disease and are inversely correlated with the risk of estrogen receptor-negative breast cancer.^{4,5}

References

1. Bushway, R.J. and Wilson, A.M. Determination of α - and β -carotene in fruit and vegetables by high performance liquid chromatography. *Can. Inst. Food Sci. Technol. J.* **15(3)**, 165-169 (1982).
2. Murakoshi, M., Takayusa, J., Kimura, O., *et al.* Inhibitory effects of α -carotene on proliferation of the human neuroblastoma cell line GOTO. *J. Natl. Cancer Inst.* **81(21)**, 1649-1652 (1989).
3. Murakoshi, M., Nishino, H., Satomi, Y., *et al.* Potent preventive action of α -carotene against carcinogenesis: spontaneous liver carcinogenesis and promoting stage of lung and skin carcinogenesis in mice are suppressed more effectively by α -carotene than by β -carotene. *Cancer Res.* **52(23)**, 6583-6587 (1992).
4. Kontush, A., Spranger, T., Reich, A., *et al.* Lipophilic antioxidants in blood plasma as markers of atherosclerosis: The role of α -carotene and γ -tocopherol. *Atherosclerosis* **144(1)**, 117-122 (1999).
5. Zhang, X., Spiegelman, D., Baglietto, L., *et al.* Carotenoid intakes and risk of breast cancer defined by estrogen receptor and progesterone receptor status: A pooled analysis of 18 prospective cohort studies. *Am. J. Clin. Nutr.* **95(3)**, 713-725 (2012).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 07/05/2018

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD

ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM