



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



Leukotriene B₄ Pathway LC-MS Mixture Item No. 22640


Supplied as: A solution in ethanol (1 µg/ml of each compound)
Storage: -20°C
Stability: ≥5 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Description and Contents

The leukotriene B₄ pathway LC-MS mixture contains 5-lipoxygenase pathway metabolites derived from arachidonic acid (Item No. 90010). The mixture is supplied in an amber ampule in which the headspace has been purged with argon to prevent lipid oxidation. This product has been designed for direct use in LC-MS applications. The solution may be used as a system suitability standard or tuning standard. After opening, we recommend that the mixture be transferred immediately to a 1 ml glass screw cap vial, to prevent solvent evaporation, and stored at -20°C. The mixture should be discarded after multiple freeze/thaw cycles.

This mixture contains 18-carboxy dinor leukotriene B₄ (Item No. 20170), 20-carboxy leukotriene B₄ (Item No. 20180), 20-hydroxy leukotriene B₄ (Item No. 20190), leukotriene B₄ (Item No. 20110), and 5(S)-HETE (Item No. 34230).

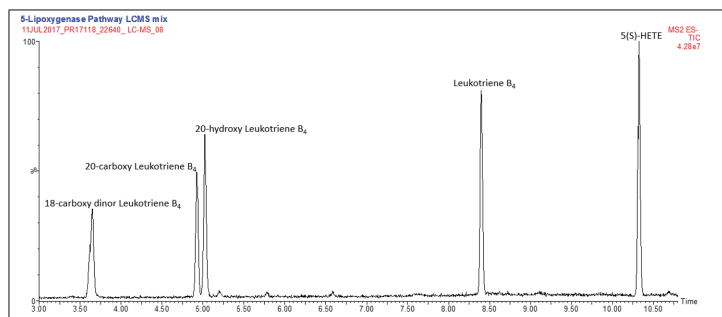

Item Number: 22640		Leukotriene B ₄ Pathway LC-MS Mixture			
Item Number	Item Name	Formula:	Mass (Da)	[MH] ⁺ (m/z)	RT (min)
20170	18-carboxy dinor Leukotriene B ₄	C ₁₈ H ₂₆ O ₆	338.4	337.3	3.65
20180	20-carboxy Leukotriene B ₄	C ₂₀ H ₃₀ O ₆	366.5	365.4	4.92
20190	20-hydroxy Leukotriene B ₄	C ₂₀ H ₃₂ O ₅	352.5	351.4	5.02
20110	Leukotriene B ₄	C ₂₀ H ₃₂ O ₄	336.5	335.4	8.40
34230	5(S)-HETE	C ₂₀ H ₃₂ O ₃	320.5	319.3	10.32

LC-MS Instrument: Waters Acquity UPLC-Xevo TQD
Mobile Phase A: Water + 0.1% Formic Acid
Mobile Phase B: Acetonitrile + 0.1% Formic Acid
Column: Waters BEH C8, 2.1 x 100 mm, 1.7 µm Flow Rate: 400 µl/min; Injection Volume: 5 µl
Sample Preparation: 200 µl of the mixture was dried under nitrogen and reconstituted in 200 µl 1:1 ACN:Water
Negative Electrospray Ionization Full Mass Scan: 317-369 m/z

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA
This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY
Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 06/16/2020

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA
PHONE: [800] 364-9897
[734] 971-3335
FAX: [734] 971-3640
CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM