



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



C6 Ceramide (d18:1/6:0) Alkyne

Item No. 24519

Formal Name: N-((2S,3R,E)-1,3-dihydroxyoctadec-4-en-17-yn-2-yl)hexanamide

Synonyms: Cer(d18:1/6:0) Alkyne, C6 Ceramide Alkyne, N-hexanoyl-D-erythro-Sphingosine Alkyne

MF: C₂₄H₄₃NO₃


FW: 393.6

Purity: ≥95%

Supplied as: A solution in ethanol

Storage: -20°C

Stability: ≥1 year

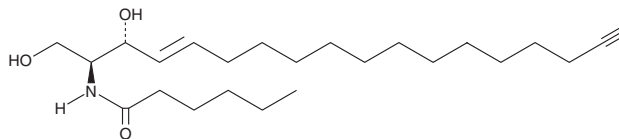

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

C6 Ceramide (d18:1/6:0) alkyne is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as dimethyl formamide (DMF) purged with an inert gas can be used. The solubility of C6 ceramide (d18:1/6:0) alkyne in DMF is approximately 0.15 mg/ml.

Description

C6 Ceramide alkyne is an ω-alkyne derivative of C6 ceramide (Item No. 62525). The terminal alkyne group can be used in click chemistry reactions to tag C6 ceramide with fluorescent or biotinylated labels for analysis of its metabolism and biological activity.

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 05/21/2019

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM