

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



CALXGLUK

Item No. 24784

CAS Registry No.: 1820596-85-7

Formal Name: (2R,2'R,2"R,3S,3'S,3"S,4S,4'S,4"S,5R,5'R,

5"R,6R,6"R,6"R)-6,6',6"-(((((7²-(heptyloxy)-

1²,3²,5²-trimethoxy-1,3,5,7(1,3)tetrabenzenacyclooctaphane-1⁵,3⁵,5⁵triyl)tris(ethane-2,1-diyl))tris(1H-1,2,3triazole-1,4-diyl))tris(methylene))tris(oxy)) tris(2-(hydroxymethyl)tetrahydro-2H-


pyran-3,4,5-triol) $C_{71}H_{95}N_9O_{22}$

MF: FW: 1,426.6 **Purity:** ≥95%

Supplied as: A crystalline solid

-20°C Storage: Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

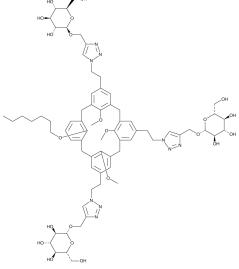


Laboratory Procedures

CALXGLUK is supplied as a crystalline solid. A stock solution may be made by dissolving the CALXGLUK in the solvent of choice, which should be purged with an inert gas. CALXGLUK is soluble in organic solvents such as methanol and DMSO. CALXGLUK is also soluble in water at a concentration of approximately 2.5 mM. We do not recommend storing the aqueous solution for more than one day.


Description

CALXGLUK is a calixarene detergent that can be used to solubilize membrane proteins.¹ It has a critical micelle concentration (CMC) of 0.025 mM. CALXGLUK has been used to solubilize a membrane GPCR (MP-A) and a membrane transporter (MP-B) from yeast alone and in combination with Triton X-100 (Item No. 600217) and digitonin (Item No. 14952).

Images

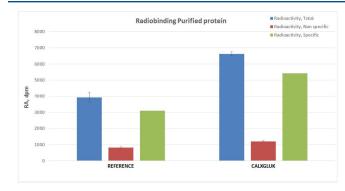

Binding of radioligand on target protein, purified in reference detergent with or without CALXGLUK as an additive. Purified protein was incubated with radioligand in absence (total, blue bars) or presence (Non Specific signal, red bars) of an excess of cold ligand. After filtration on GF/C membranes and washing, scintillation agent was added and radioactivity was detected using a Microbeta2. Specific radioactivity (green bars) corresponds to (total signal) – (non-specific signal).

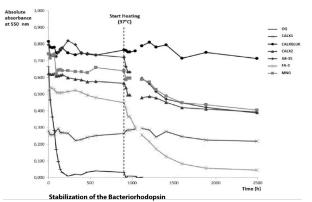

After solubilization and purification in OG, reagents we added as stabilizers. Bacteriorhodopsin's activity was monitored over time at first 25°C and then 37°C by measuring absorbance at 560nm.

1. Desuzinges Mandon, E., Agez, M., Pellegrin, R., et al. Anal. Biochem. 517, 40-49 (2017).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 10/08/2019

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA **PHONE:** [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM