



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



13-methyl Myristic Acid methyl ester

Item No. 24812

CAS Registry No.: 5129-59-9
Formal Name: 13-methyl-tetradecanoic acid, methyl ester
Synonyms: Methyl Isopentadecanoate,
13-methyl Tetradecanoate,
13-methyl Tetradecanoic Acid methyl ester

MF: C₁₆H₃₂O₂


FW: 256.4

Purity: ≥98%

Supplied as: A solution in ethanol

Storage: -20°C

Stability: ≥1 year

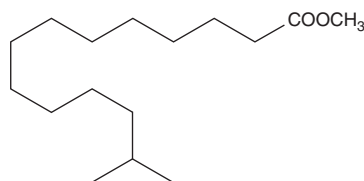

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

13-methyl Myristic acid methyl ester is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. The solubility of 13-methyl myristic acid methyl ester in these solvents is approximately 20 and 30 mg/ml, respectively.

Description

13-methyl Myristic acid methyl ester is a methylated bacterial fatty acid methyl ester.^{1,2} It has been used as a marker of bacteria in sewage sludges and of bacterial contaminants in pharmaceutical manufacturing by GC-MS.^{3,4}

References

1. Murugan, K., Prabu, R.V., Sangeetha, S., *et al.* Antiviral activity of *Cardiospermum Halicacabum* L. extract against coinfecting agents HIV and HBV. *J. Herbs Spices Med. Plants* **17(4)**, 403-418 (2011).
2. de Blas, E., Almendros, G., and Sanz, J. Molecular characterization of lipid fractions from extremely water-repellent pine and eucalyptus forest soils. *Geoderma* **206(2013)**, 75-84 (2013).
3. Jardé, E., Mansuy, L., and Faure, P. Organic markers in the lipidic fraction of sewage sludges. *Water Res.* **39(7)**, 1215-1232 (2005).
4. Olson, W.P., Groves, M.J., and Klegerman, M.E. Identifying bacterial contaminants in a pharmaceutical manufacturing facility by gas chromatographic fatty acid analysis. *Pharm. Tech.* **14(2)**, 32-36 (1990).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 04/02/2019

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM