

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



LgtC (N. meningitidis) Glycosyltransferase

Item No. 25002

Overview and Properties

Synonyms: Lipopolysaccharyl α-galactosyltransferase, Lipopolysaccharyl-1,4-galactosyltransferase

Source: Active recombinant C-terminal His-tagged LgtC expressed in E. coli

Amino Acids: 1-286 **Uniprot No.:** Q93EK7 Molecular Weight: 33.4 kDa

Storage: -80°C (as supplied)

Stability:

Purity: ≥90% estimated by SDS-PAGE

50 mM Tris HCl, pH 7.5, with 3 mM DTT Supplied in:

Protein

Concentration: batch specific mg/ml batch specific U/ml Activity: Specific Activity: batch specific U/mg


Unit Definition: One unit is defined as the amount of LgtC enzyme require to produce 1 µmol of

globotriose (Gb3) from lactose per minute at pH 7.5 at 30°C in the presence of 10 mM

of UDP-galactose.

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Images

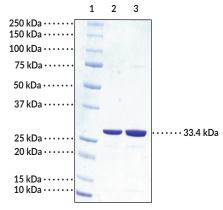

Lane 1: MW Markers Lane 2: LgtC (2 µg) Lane 3: LgtC (4 µg)

SDS-PAGE Analysis of LgtC.

Representative gel image shown; actual purity may vary between each batch.



Lane 2: Galactose (control)

Lane 3: Negative control (no enzyme control)

Lane 4: Reaction mixture with enzyme

TLC Analysis of Lactose Conversion by Purified LgtC. Assay was performed in 500 μl volume containing 25 µg of enzyme with lactose as substrate (2.0 mM) in the assay buffer (100 mM Tris HCI buffer, pH 7.5, containing 5 mM DTT, 1 mM MnCl₂). The reaction was incubated in a thermoshaker (300 rpm) at 30°C for a total of 15 minutes. Analysis was done in a TLC plate (Analtech, 2.5 x 10 cm) using the solvents ethyl acetate:methanol:water:acetic acid at a ratio of 5:2:1.5:0.5 v/v.

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA
This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 02/29/2024

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.CAYMANCHEM.COM

PRODUCT INFORMATION



Description

Lipopolysaccharide α -1,4-galactosyltransferase (LgtC) is a retaining glycosyltransferase and member of glycosyltransferase family 8 (GT8) that is found in *N. meningitidis*, an opportunistic bacterium that can cause meningitis. It is composed of an N-terminal mixed α/β domain, which contains the active site, and a C-terminal helical domain, which is responsible for attachment to cell membranes. LgtC transfers a galactosyl moiety from UDP-galactose (UDP-Gal) to the terminal lactose moiety on lipooligosaccharides (LOSs) while retaining the configuration of the donor sugar glycosidic bond. Isolated LgtC has been used in the analysis of newborn dried blood spots as a potential screening assay for galactosemia, an inborn error of metabolism characterized by a deficiency in galactose-1-phosphate uridylyltransferase (GALT), the enzyme that converts UDP-glucose to UDP-galactose. Cayman's LgtC (*N. meningitidis*) protein can be used for biocatalysis and enzyme activity assay applications.

References

- 1. Ly, H.D., Lougheed, B., Wakarchuk, W.W., et al. Mechanistic studies of a retaining α-galactosyltransferase from *Neisseria meningitidis*. *Biochemistry* **41(16)**, 5075-5085 (2002).
- 2. Persson, K., Ly, H.D., Dieckelmann, M., et al. Crystal structure of the retaining galactosyltransferase LgtC from 'Neisseria meningitidis' in complex with donor and acceptor sugar analogs. *Nat. Struct. Biol.* **8(2)**, 166-175 (2001).
- 3. Chan, P.H.W., Weissbach, S., Okon, M., *et al.* Nuclear magnetic resonance spectral assignments of α-1,4-galactosyltransferase LgtC from *Neisseria meningitidis*: Substrate binding and multiple conformational states. *Biochemistry* **51(41)**, 8278-8292 (2012).
- 4. Hong, X., Kumar, A.B., Scott, C.R., *et al.* Multiplex tandem mass spectrometry assay for newborn screening of X-linked adrenoleukodystrophy, biotinidase deficiency, and galactosemia with flexibility to assay other enzyme assays and biomarkers. *Mol. Genet. Metab.* **124(2)**, 101-108 (2018).

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897