



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



Acridine Orange 10-nonyl (bromide)

Item No. 25545

CAS Registry No.: 75168-11-5
Formal Name: 3,6-bis(dimethylamino)-10-nonyl-acridinium, monobromide
Synonym: 10-nonyl Acridine Orange
MF: C₂₆H₃₈N₃• Br
FW: 472.5
Purity: ≥90%
Ex./Em. Max: 490/510 nm
Supplied as: A solid
Storage: -20°C
Stability: ≥4 years

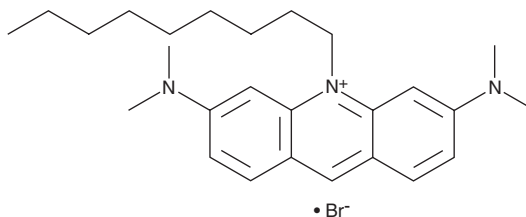

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Acridine orange 10-nonyl (bromide) is supplied as a solid. A stock solution may be made by dissolving the acridine orange 10-nonyl (bromide) in the solvent of choice. Acridine orange 10-nonyl (bromide) is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide, which should be purged with an inert gas.

Description

Acridine orange 10-nonyl is a mitochondrial dye and derivative of acridine orange (Item No. 14338) that incorporates into the cardiolipin-rich inner mitochondrial membrane.¹ It selectively binds to cardiolipin-containing liposomes over phosphatidylserine- and phosphatidylinositol-containing liposomes *in vitro* (K_s = 0.002, 7,000, and 7,000 mM, respectively). Acridine orange 10-nonyl binding to mitochondria is abolished by the mitochondrial uncoupler FCCP (Item No. 15218) in rat cortical astrocytes, neonatal cardiomyocytes, and isolated brain mitochondria, indicating dependence on mitochondrial membrane potential.² Acridine orange 10-nonyl displays excitation/emission maxima of 490/510 nm, respectively, and has been used to quantify cardiolipin content and mitochondrial mass in isolated mitochondria and cellular systems.

References

1. Petit, J.M., Maftah, A., Ratinaud, M.H., *et al.* 10N-nonyl acridine orange interacts with cardiolipin and allows the quantification of this phospholipid in isolated mitochondria. *Eur. J. Biochem.* **209(1)**, 267-273 (1992).
2. Jacobson, J., Duchon, M.R., and Heales, S.J. Intracellular distribution of the fluorescent dye nonyl acridine orange responds to the mitochondrial membrane potential: Implications for assays of cardiolipin and mitochondrial mass. *J. Neurochem.* **82(2)**, 224-233 (2002).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/28/2022

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM